



Cadenza CD-C18

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Vitamin B12

ビタミンB12のLC-MS/MS分析

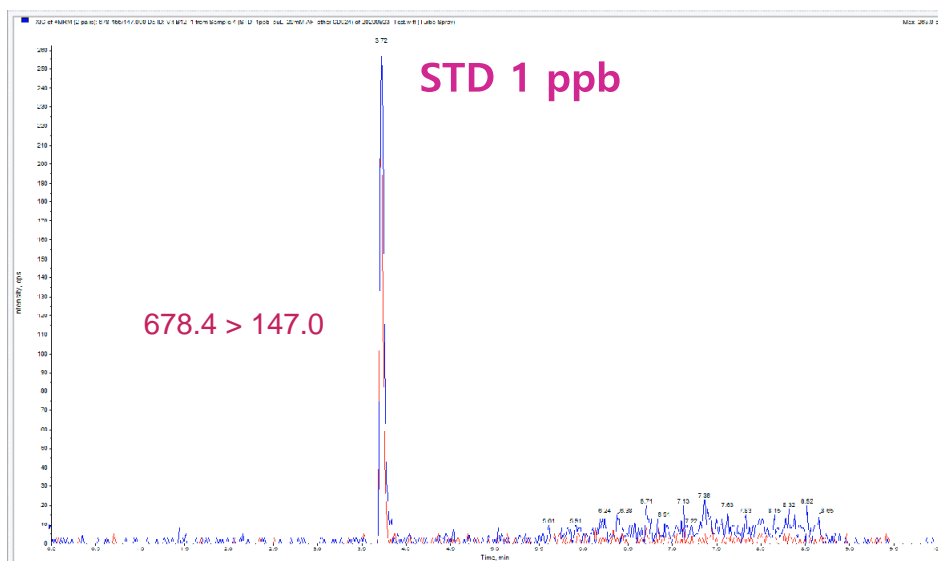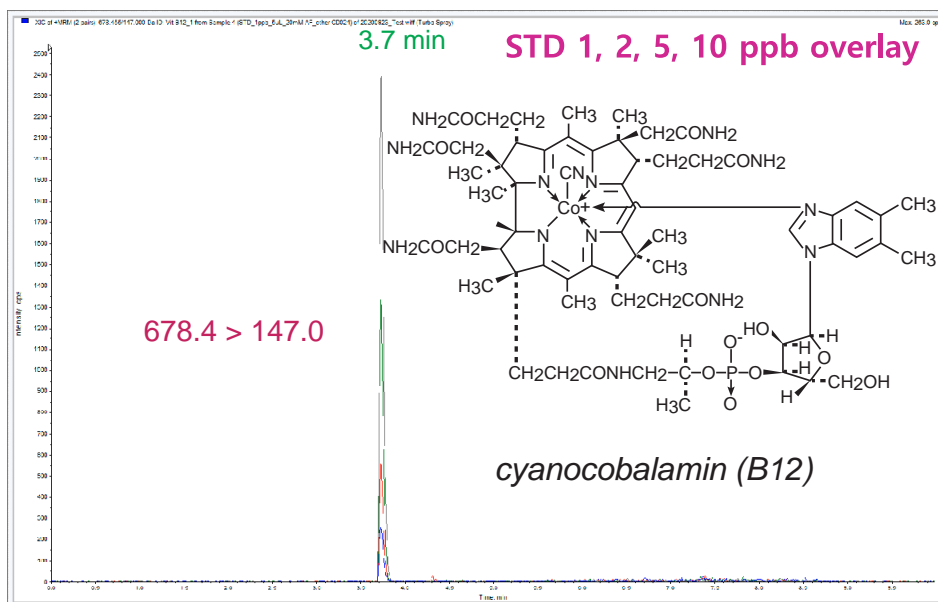

Cadenza CD-C18, 100 x 2 mm

A: 20mM ammonium formate

B: acetonitrile

5%B (0-1min), 5-80%B (1-4min), 80%B (4-7min),

80-5%B (7-7.1min), 5%B (7.1-10min)

0.3mL/min, 35°C, Triple Quad. MS (ESI, Positive)

5µL (1, 2, 5, 10 ppb)

Courtesy of PEAKMAN SP, KOREA