



Unison UK-C18

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Stevia sweeteners

ステビア甘味料の LC-MS/MS 分析

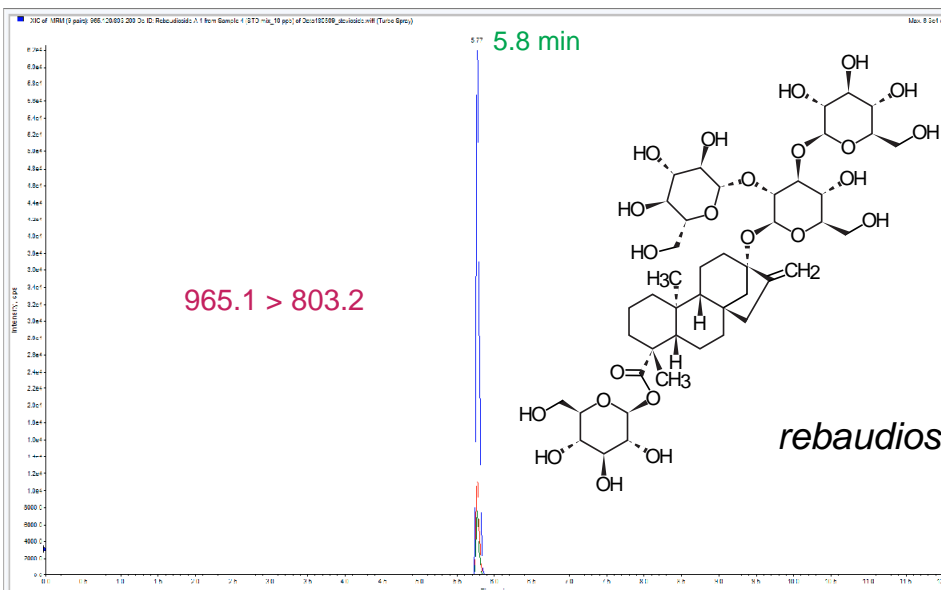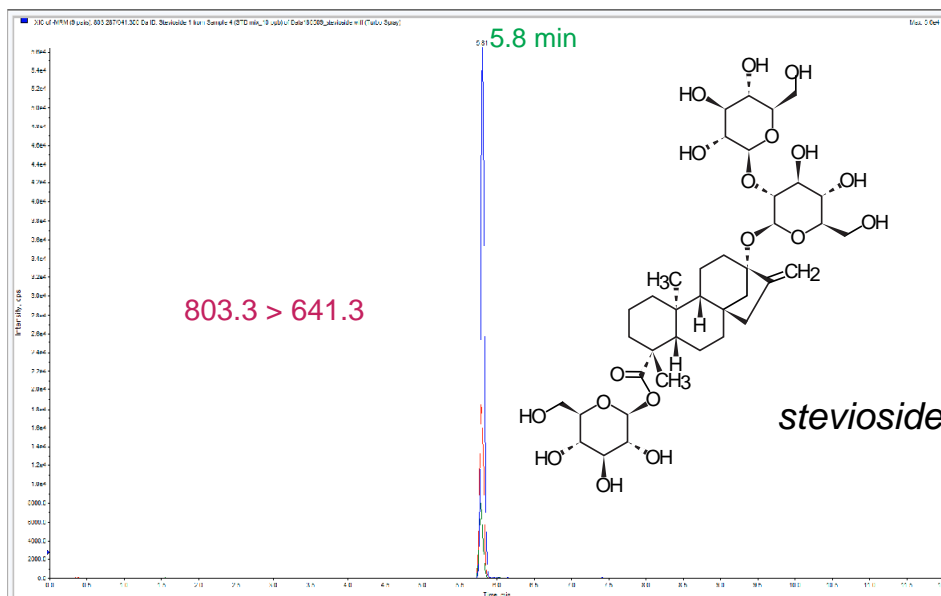

Unison UK-C18, 100 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid in acetonitrile

10%B (0-2min), 10-90%B (2-5min), 90%B (5-8min),

90-10%B (8-8.1min), 10%B (8.1-12min)

0.2mL/min, 40°C, Triple Quad. MS (ESI, Negative)

5µL (10 ppb)

Courtesy of PEAKMAN SP, KOREA