


Cadenza CD-C18

50 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Pesticides (1)

殺虫剤の LC-MS/MS 分析 (1)

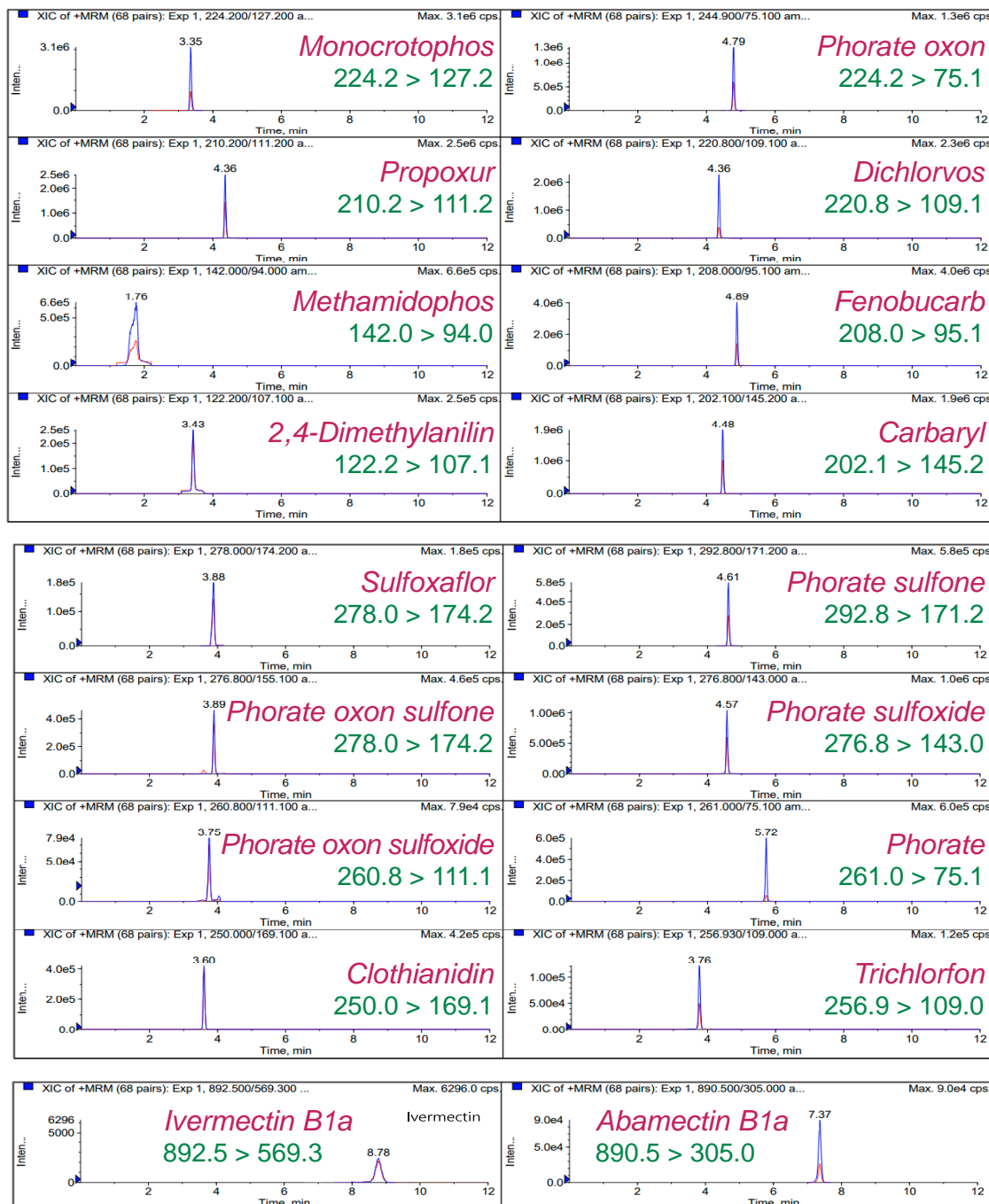

Cadenza CD-C18, 50 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid + 5mM ammonium acetate in methanol

5-90%B (0-3min), 90%B (3-7min), 90-5%B (7-7.1min), 5%B (7.1-12min)

0.2mL/min, 40°C, Triple Quad. MS (ESI, Positive), 2µL

Courtesy of PEAKMAN SP, KOREA