


Cadenza CD-C18

50 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Imidacloprid compounds

イミダクロプリド類のLC-MS/MS分析

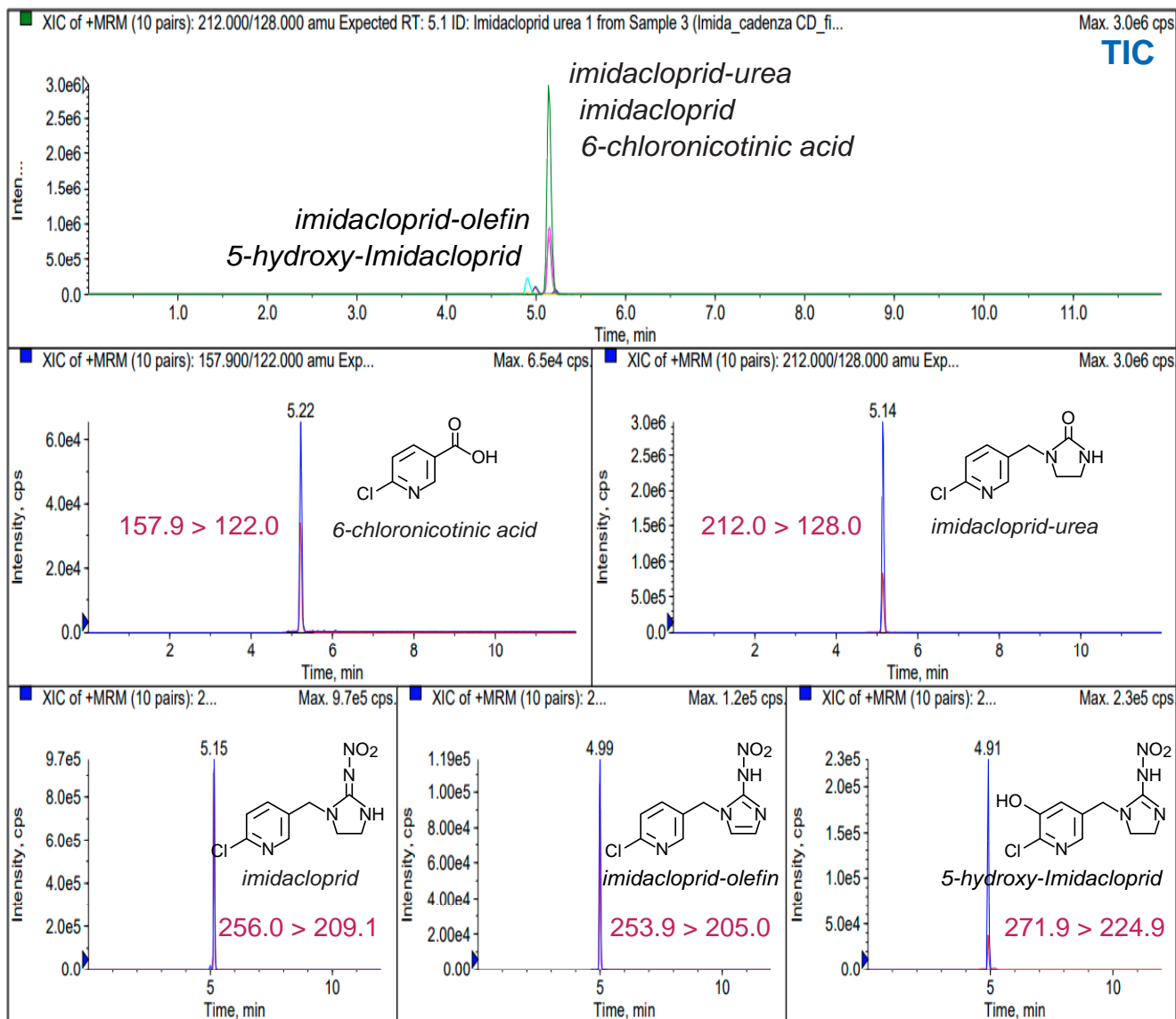

Cadenza CD-C18, 50 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid + 5mM ammonium acetate in methanol

5%B (0-2min), 5-90%B (2-4min), 90-70%B (4-4.5min), 70%B (4.5-6min)

70-5%B (6-7min), 5%B (7-12min)

0.2mL/min, 40°C, Triple Quad. MS (ESI, Positive)

2µL

Courtesy of PEAKMAN SP, KOREA