




Unison UK-C18

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Flonicamid Metabolites

フロニカミド代謝物の LC-MS/MS 分析

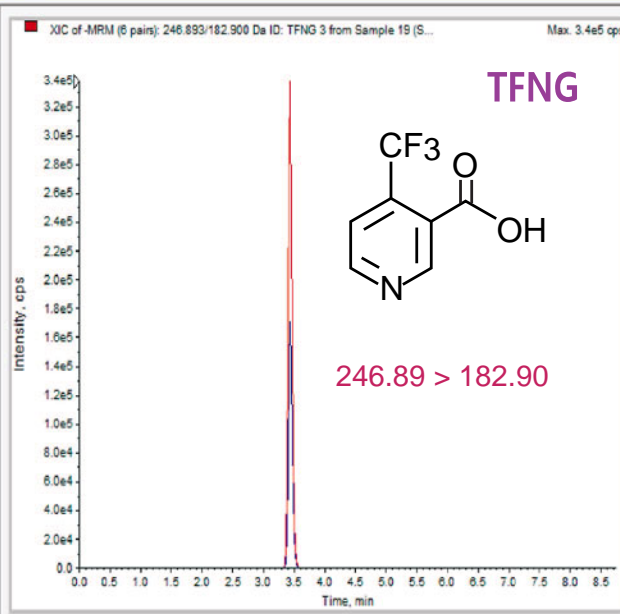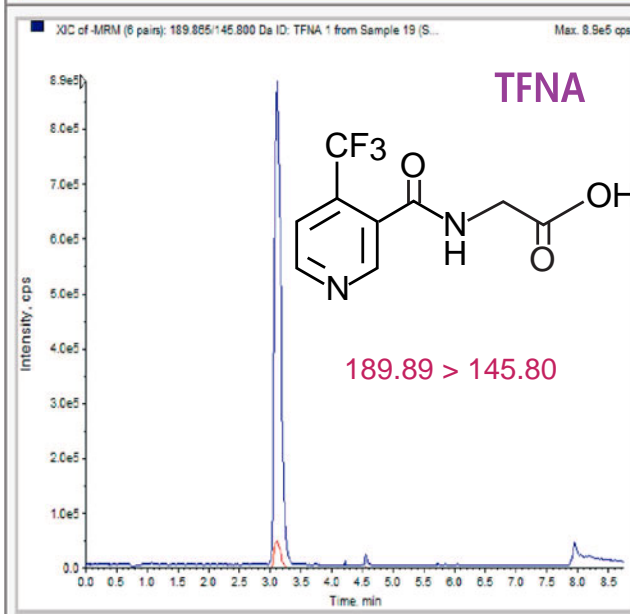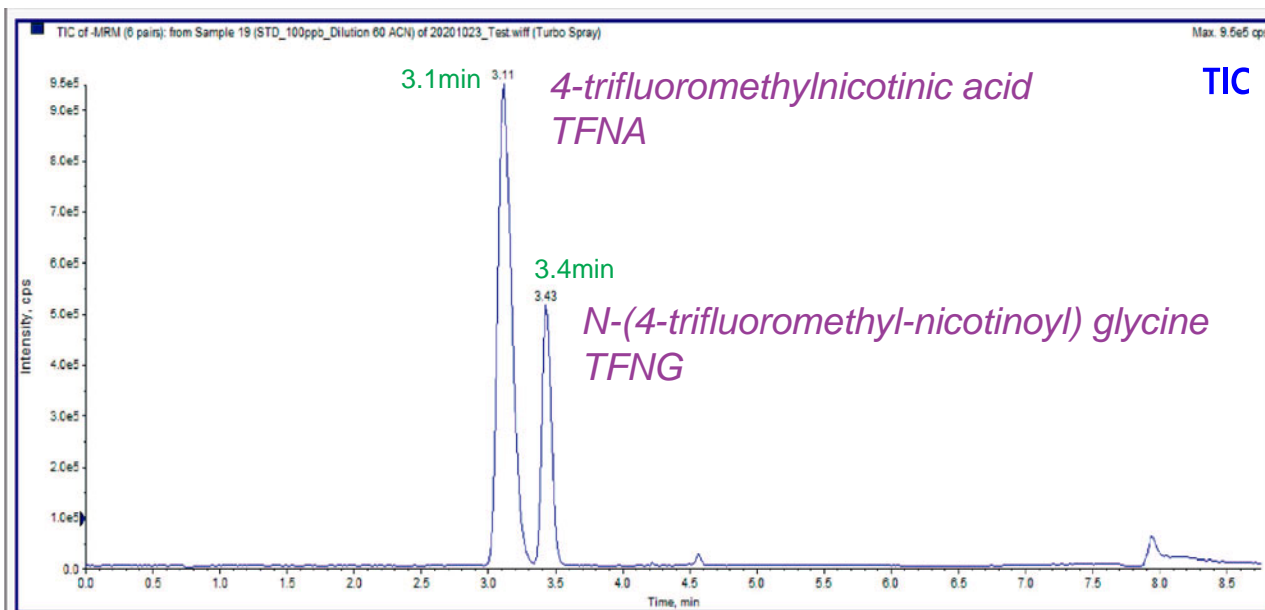

Unison UK-C18, 100 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid in acetonitrile

10%B (0-2min), 10-40%B (2-4min), 40-80%B (4-5min), 80%B (5-7min),  
80-10%B (7-7.1min), 10%B (7.1-10min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Negative)

5µL (100 ppb)

Courtesy of PEAKMAN SP, KOREA