


Cadenza CD-C18 HT

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Pesticides (3)

殺虫剤の LC-MS/MS 分析 (3)

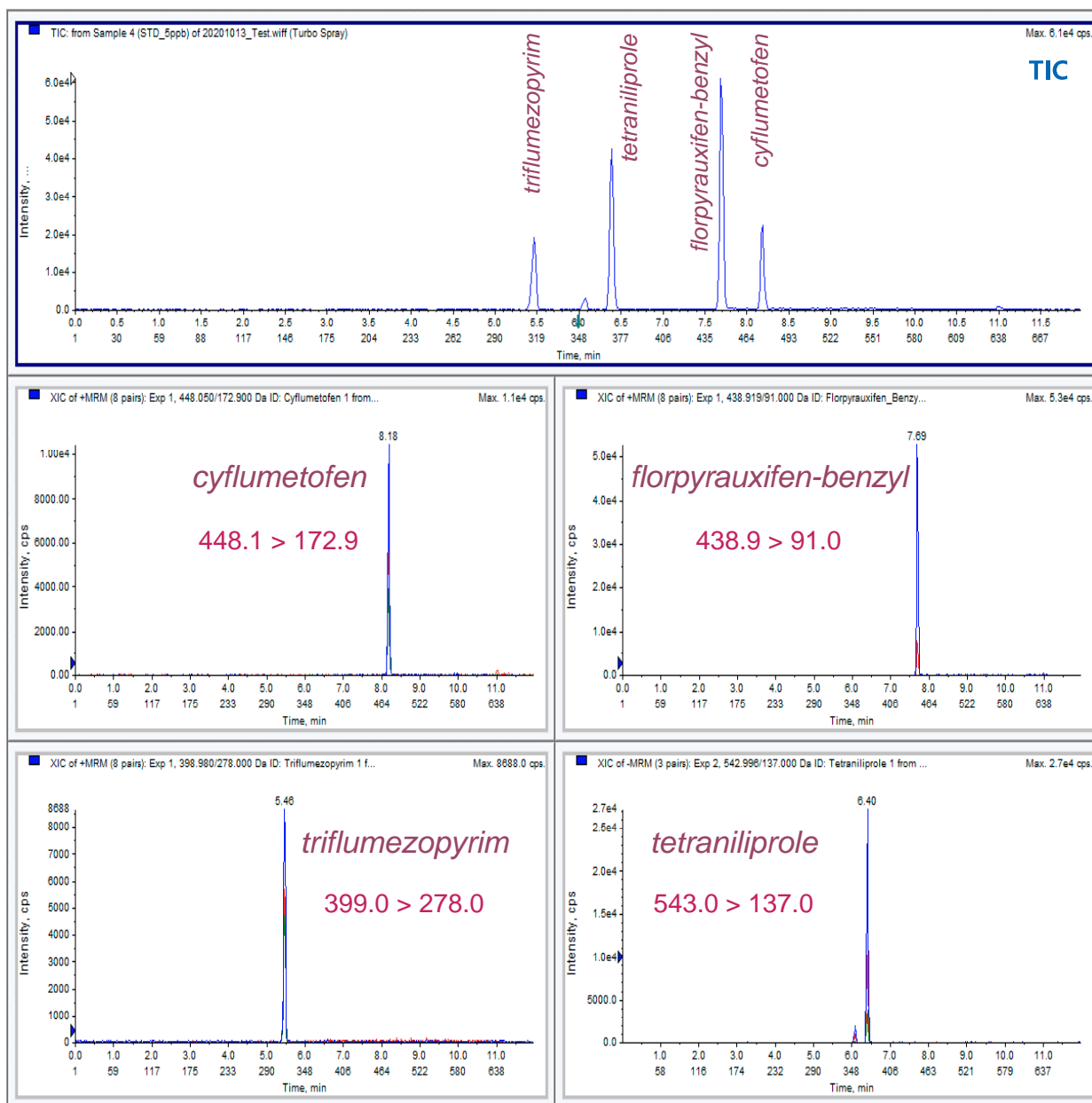

Cadenza CD-C18 HT, 150 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid in acetonitrile

10%B (0-1min), 10-95%B (1-6min), 95%B (6-9min),

95-10%B (9-9.1min), 10%B (9.1-12min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive/Negative)

5µL (5 ppb)

Courtesy of PEAKMAN SP, KOREA