

Unison UK-C18


50 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Pesticides (4)

殺虫剤の LC-MS/MS 分析 (4)

STD 5 mix 5, 10, 20, 40, 80 ppb Overlay

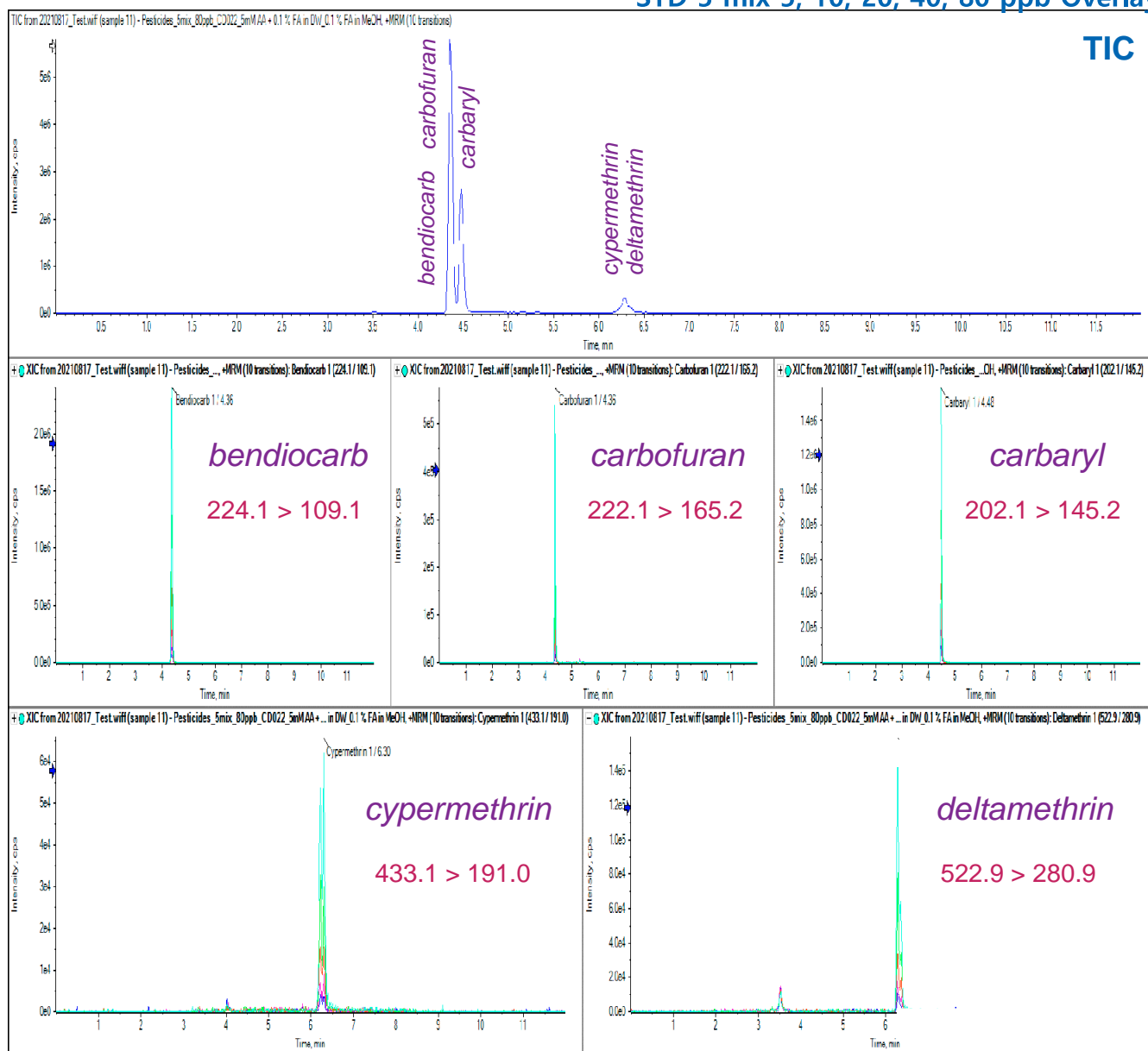

Unison UK-C18, 50 x 2 mm

A: 0.1% formic acid + 5mM ammonium acetate in water

B: 0.1% formic acid in methanol

5%B (0-0.5min), 5-90%B (0.5-4min), 90%B (4-7min),

90-5%B (7-7.1min), 5%B (7.1-12min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)

2µL (5, 10, 20, 40, 80 ppb)

Courtesy of PEAKMAN SP, KOREA