
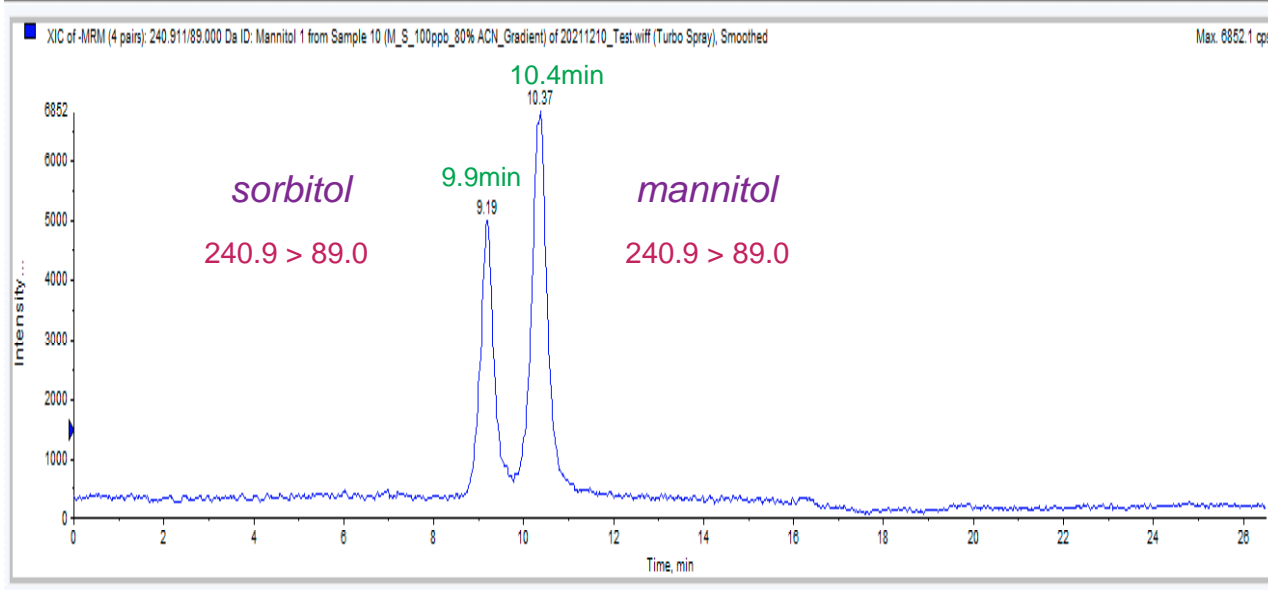


Unison UK-Amino HT 150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Sugar alcohols

糖アルコールの LC-MS/MS 分析

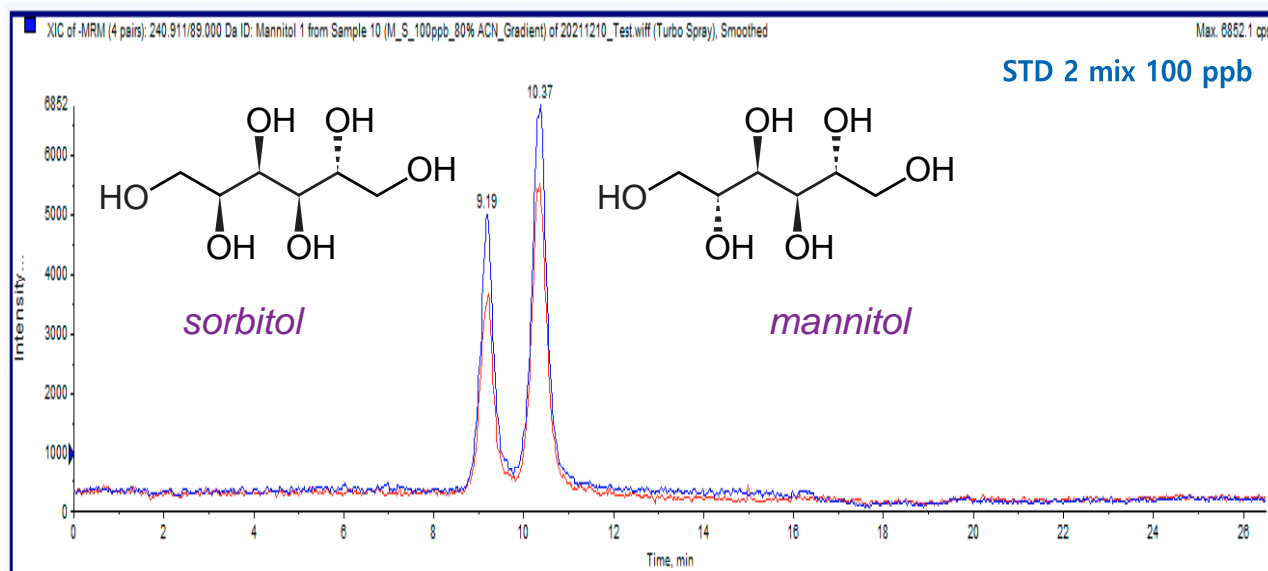

Unison UK-Amino HT, 150 x 2 mm

A: acetonitrile

B: 10mM ammonium acetate

10%B (0-7min), 10-12%B (7-12min), 12-80%B (12-14min),

80%B (14-17min), 80-10%B (17-17.1min), 10%B (17.1-26min)

0.2mL/min, 40°C, Triple Quad. MS (ESI, Negative)

5µL (100 ppb)

Courtesy of PEAKMAN SP, KOREA