




Cadenza CD-C18

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Cannabinoid

カンナビノイドの LC-MS/MS 分析

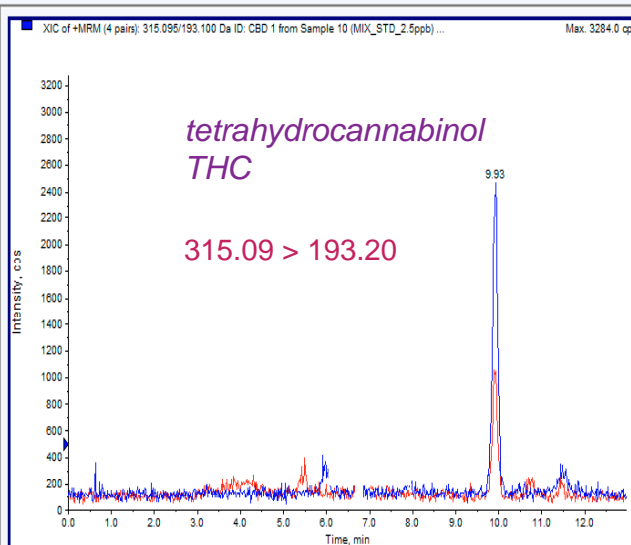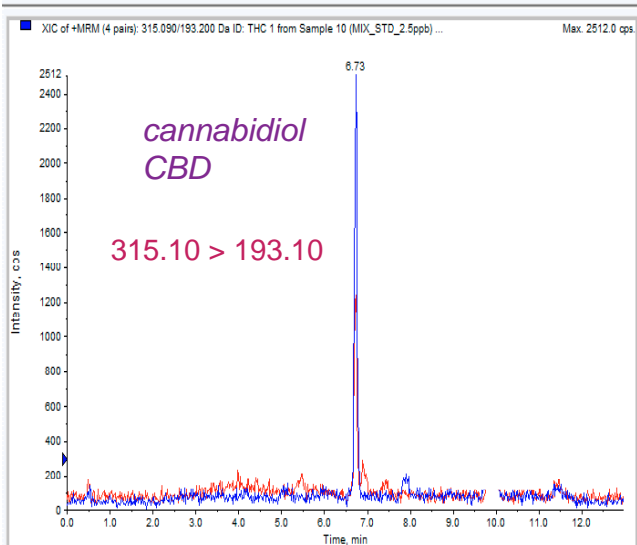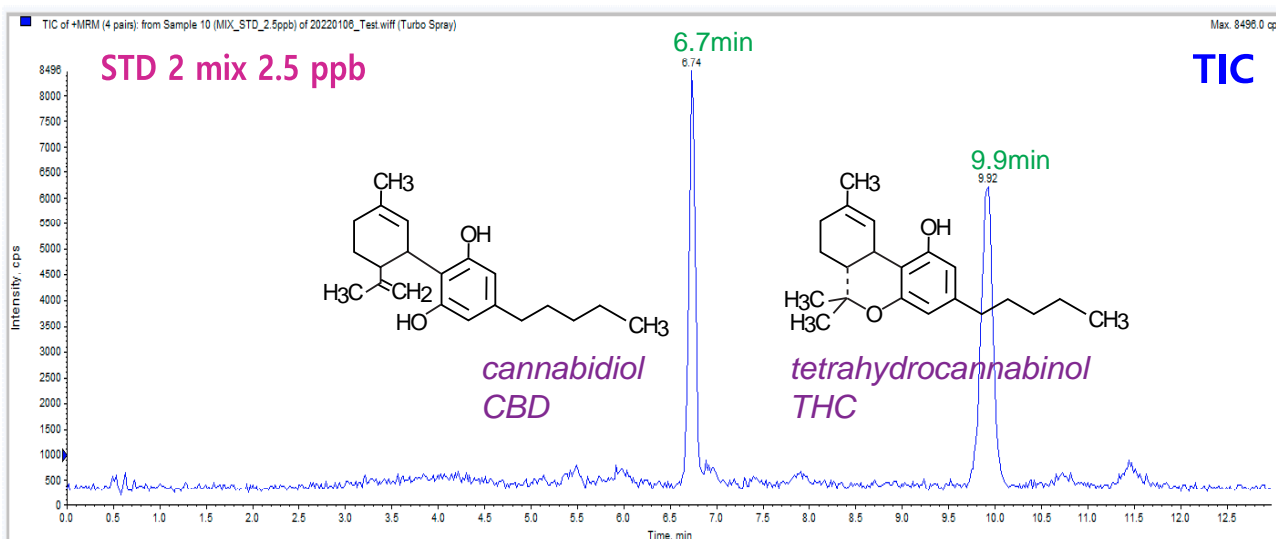

Cadenza CD-C18, 100 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid in acetonitrile

30%B (0-0.5min), 30-70%B (0.5-2.5min), 70%B (2.5-10min),

70-30%B (10-10.1min), 30%B (10.1-13.1min)

0.4mL/min, 40°C, Triple Quad. MS (ESI, Positive)

5µL (2.5 ppb)

Courtesy of PEAKMAN SP, KOREA