


Unison UK-C8

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Preservatives

保存料の LC-MS/MS 分析

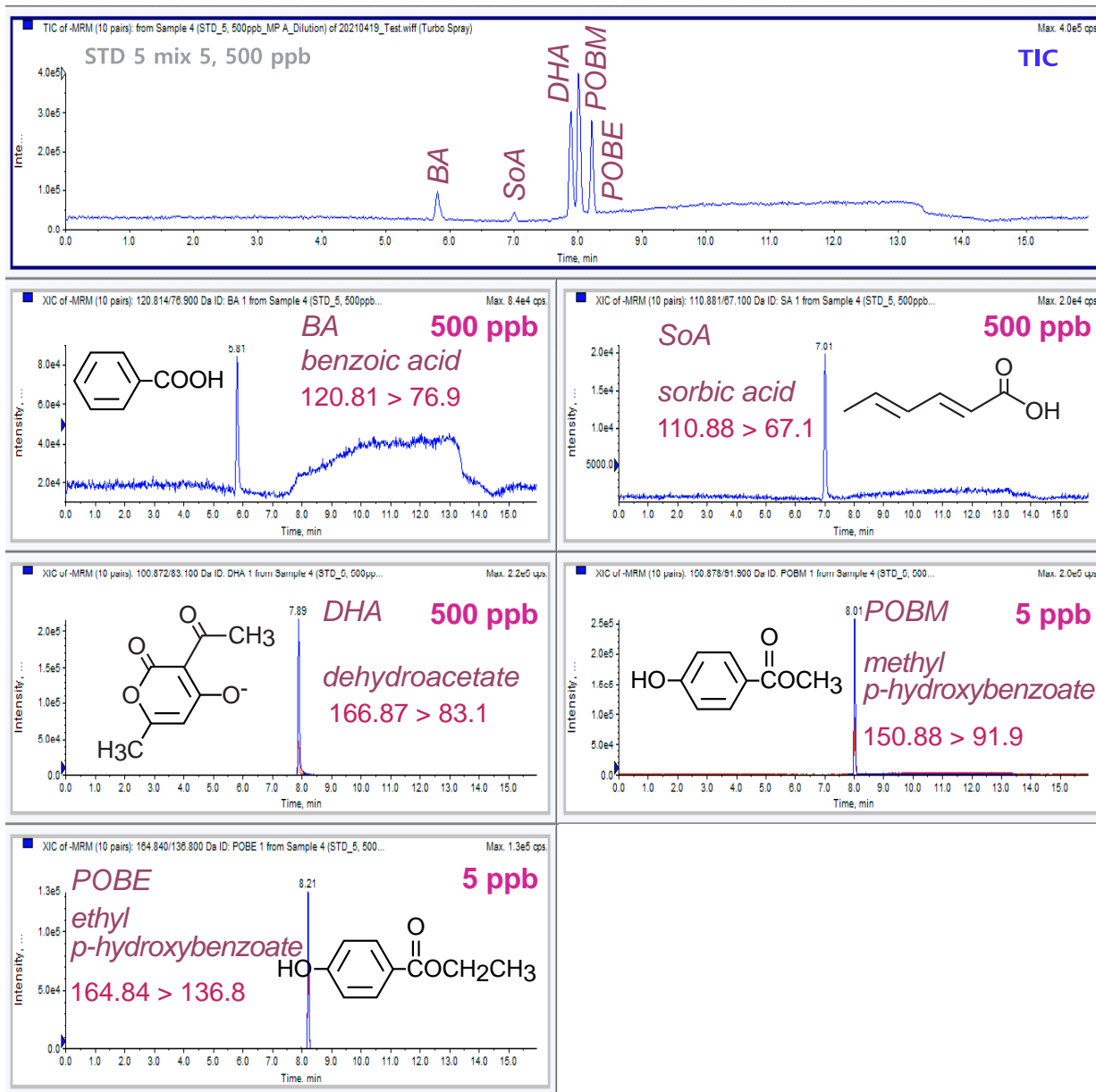

Unison UK-C8, 100 x 2 mm

A: 0.1% acetic acid + 20mM ammonium acetate

B: acetonitrile

5%B (0-1min), 5-25%B (1-5min), 25-90%B (5-6min), 90%B (6-11min),

90-5%B (11-11.1min), 5%B (11.1-16min)

0.2mL/min, 40°C, Triple Quad. MS (ESI, Negative)

2µL (5 or 500 ppb)

Courtesy of PEAKMAN SP, KOREA