


Cadenza CD-C18

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Butylated hydroxytoluene (BHT)

ブチル化ヒドロキシトルエンの LC-MS/MS 分析

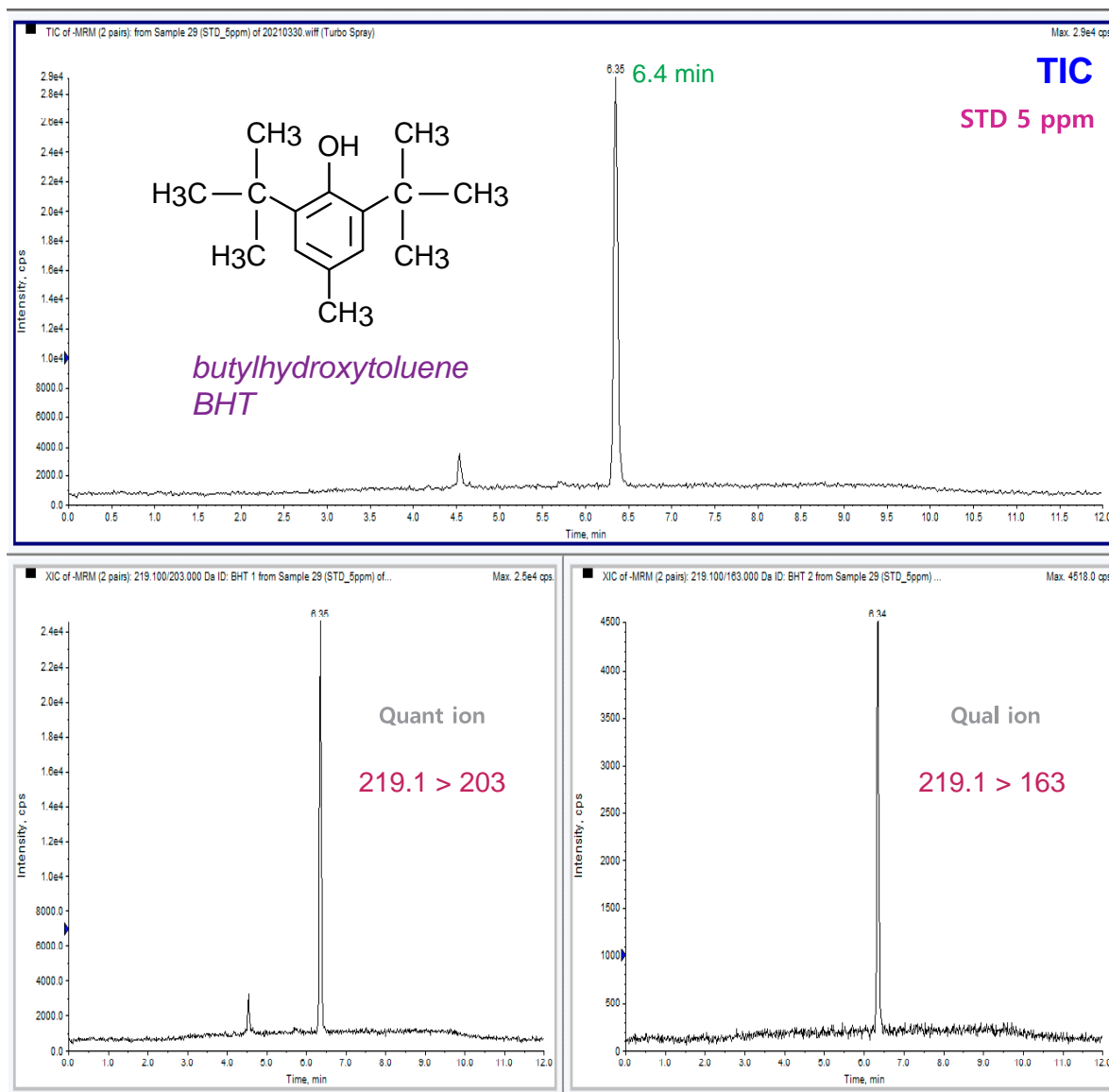

Cadenza CD-C18, 100 x 2 mm

A: 10mM ammonium acetate / acetonitrile = 90 / 10

B: 10mM ammonium acetate / acetonitrile = 10 / 90

10%B (0-0.5min), 10-100%B (0.5-3min), 100%B (3-8min),

100-10%B (8-8.1min), 10%B (8.1-12min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Negative)

5µL (5 ppb)

Courtesy of PEAKMAN SP, KOREA