


**Cadenza CD-C18 HT** 150 x 2 mm **HPLC COLUMN APPLICATION**

**LC-MS/MS: Antibiotics (1)**  
**抗生物質の LC-MS/MS 分析 (1)**

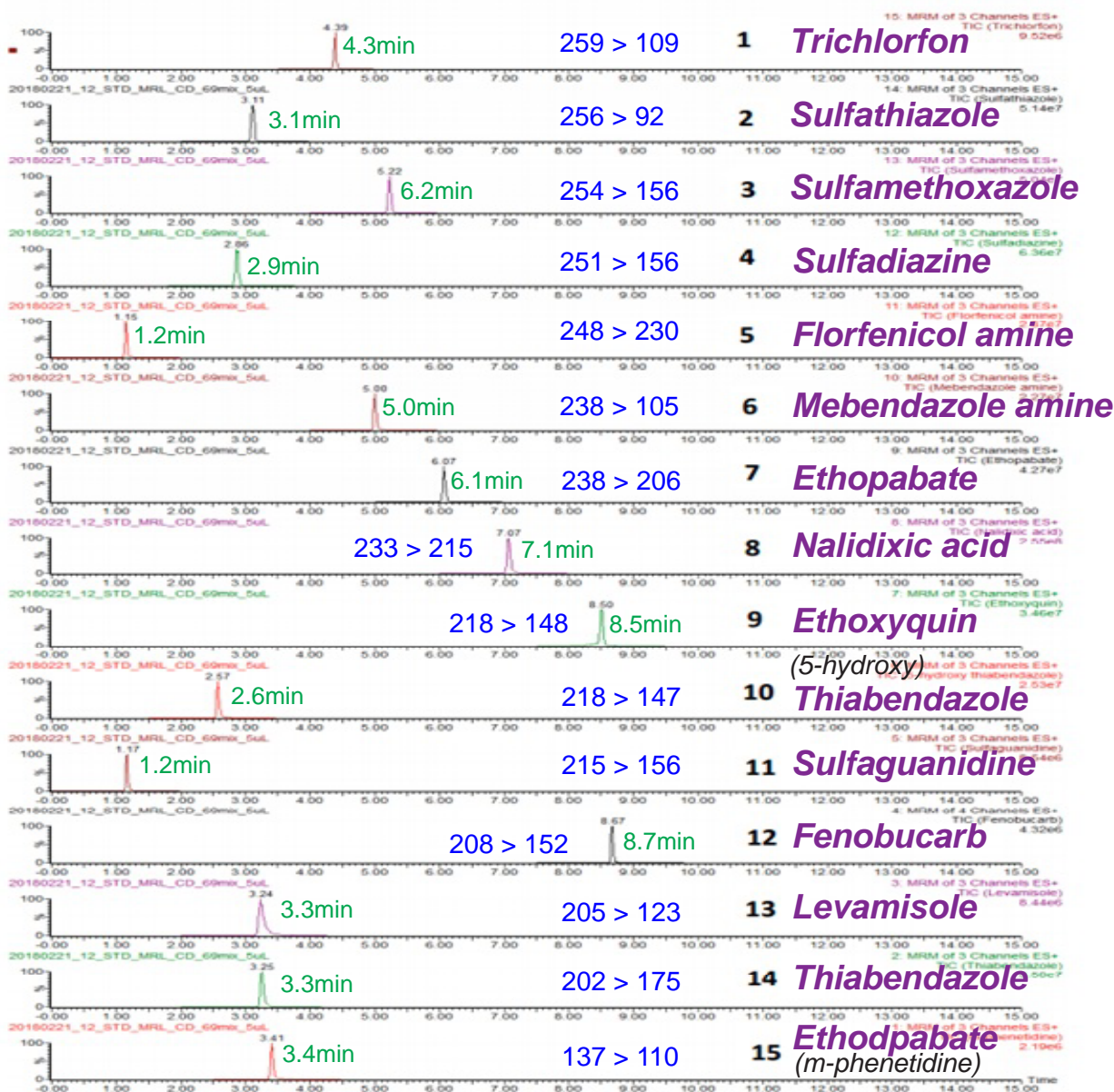

Cadenza CD-C18 HT, 150 x 2 mm  
 A: 0.1% formic acid in water  
 B: 0.1% formic acid + 2mM ammonium formate in acetonitrile  
 10%B (0-1min), 10-40%B (1-6min), 40-95%B (6-10min),  
 95%B (10-15min), 95-10%B (15-15.1min), 10%B (15.1-20min)  
 0.4mL/min, 40°C, Triple Quad. MS (ESI, Positive/Negative)  
 5µL

Courtesy of PEAKMAN SP, KOREA