




Intrada Amino Acid

50 x 3 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Ornithine

オルニチンの LC-MS/MS 分析

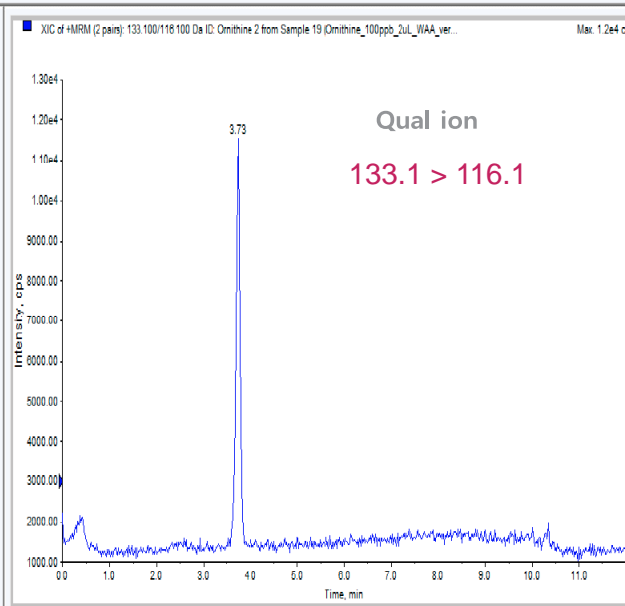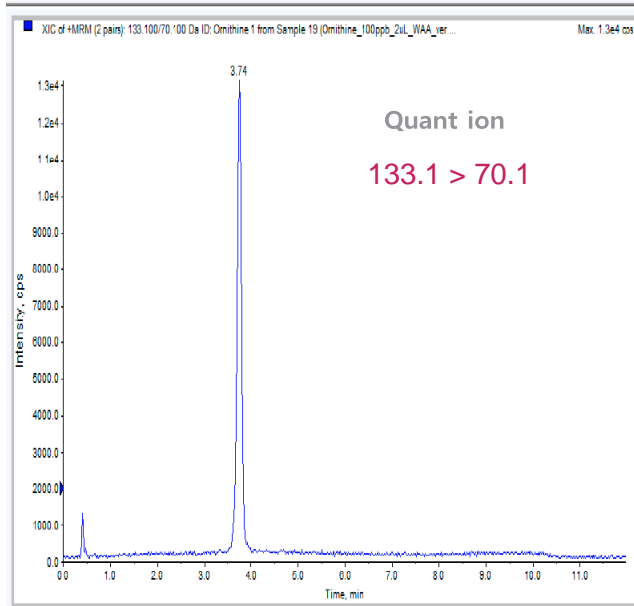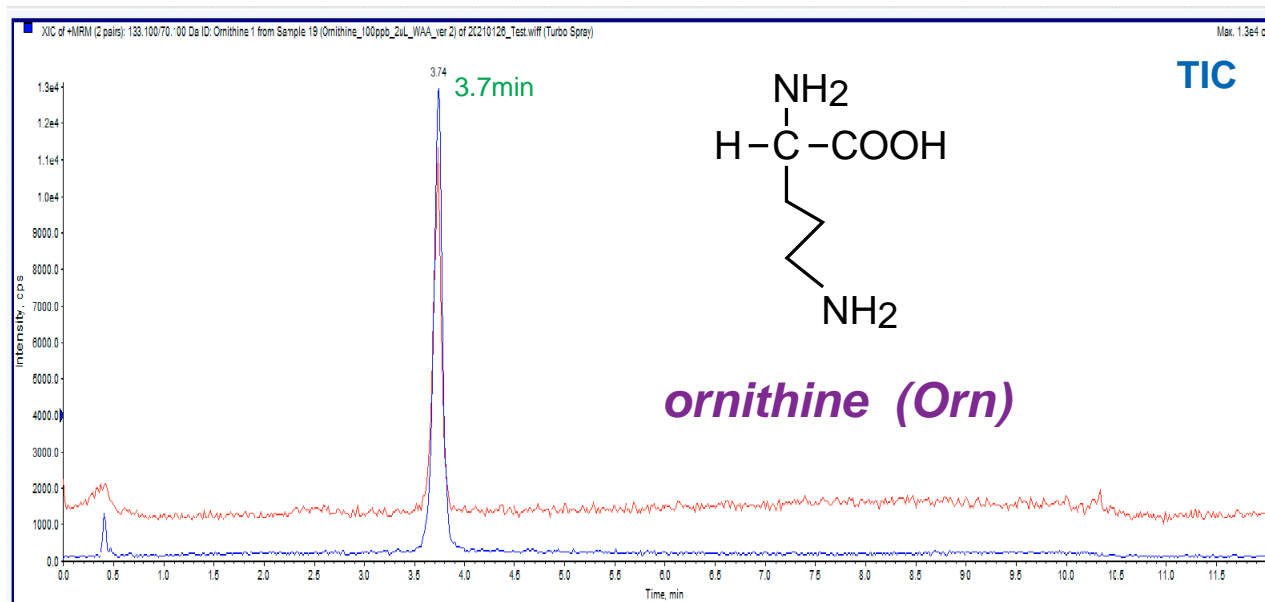

Intrada Amino Acid, 50 x 3 mm
A: 0.1% formic acid in acetonitrile
B: 100mM ammonium formate
60%B (0-1min), 60-100%B (1-7min), 100%B (7-9min),
100-60%B (9-9.1min), 60%B (9.1-12min)
0.5mL/min, 40°C, Triple Quad. MS (ESI, Positive)
2µL (5 ppb)

Courtesy of PEAKMAN SP, KOREA