





Unison UK-C18

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Acrylamide

アクリルアミドの LC-MS/MS 分析

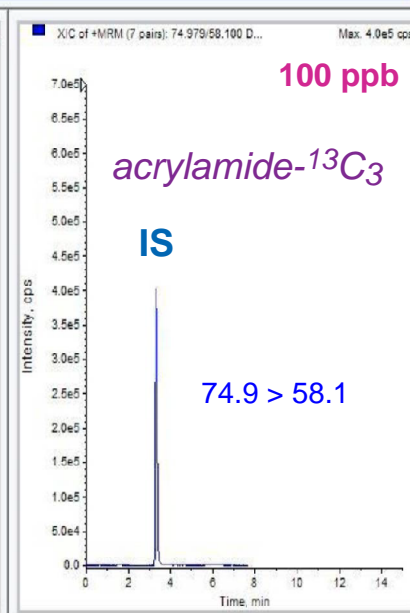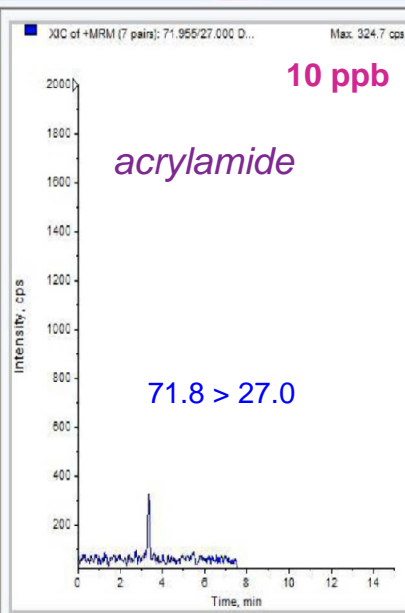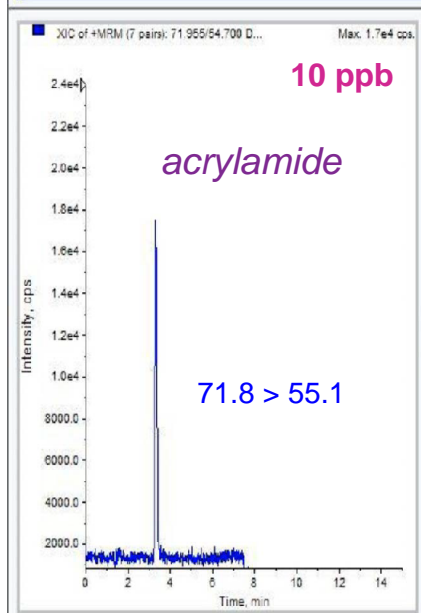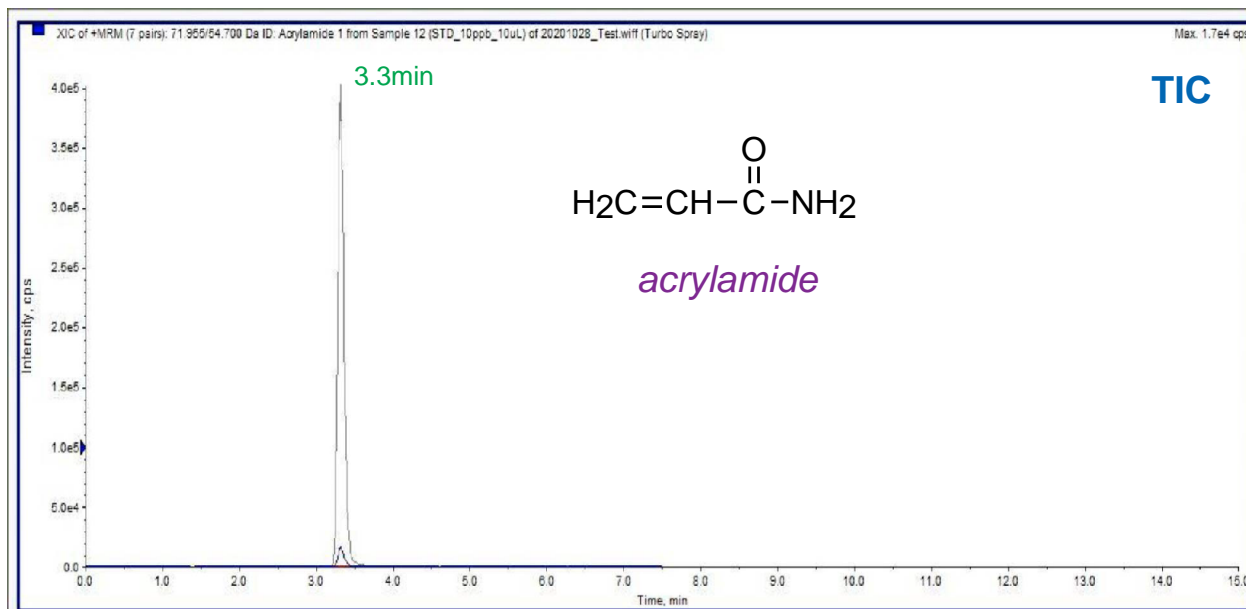

Unison UK-C18, 150 x 2 mm
 A: 0.1% formic acid in water
 B: methanol
 0%B (0-2min), 0-90%B (2-6min), 90%B (6-8min),
 90-0%B (8-8.1min), 0%B (8.1-15min)
 0.2mL/min, 40°C, Triple Quad. MS (ESI, Positive)
 10µL (5 ppb)

Courtesy of PEAKMAN SP, KOREA