




Unison UK-C18 HT

50 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Carotenoids

カロテノイドの LC-MS/MS 分析

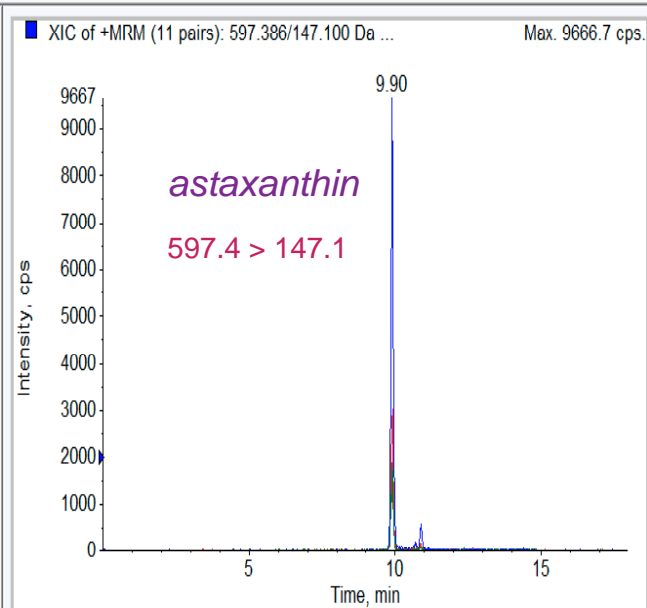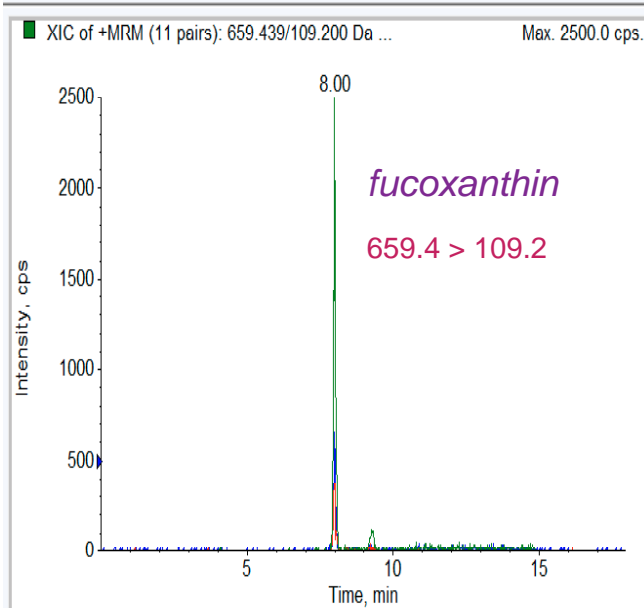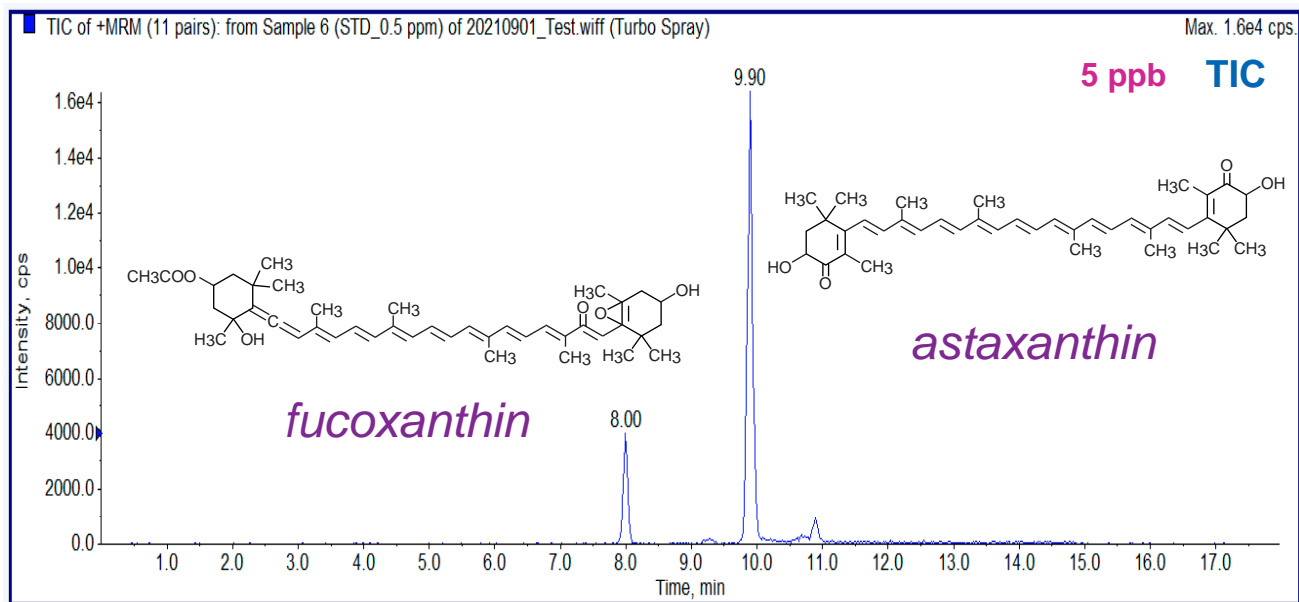

Unison UK-C18 HT, 50 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid in acetonitrile

10%B (0-0.5min), 10-70%B (0.5-3min), 70-95%B (3-11min),

95%B (11-13min), 95-10%B (13-13.1min), 10%B (13.1-18min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)

5µL (5 ppb)

Courtesy of PEAKMAN SP, KOREA