




**Cadenza CD-C18 HT**

150 x 2 mm

**HPLC COLUMN APPLICATION**

**LC-MS/MS: Ethoxyquin**

エトキシキンの LC-MS/MS 分析

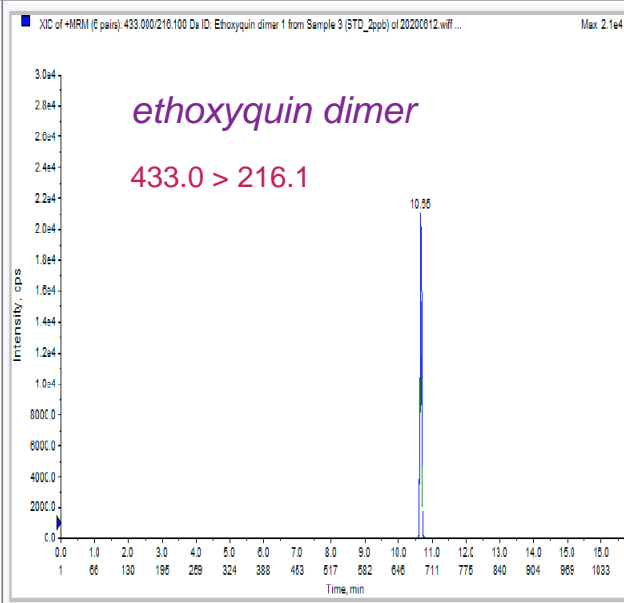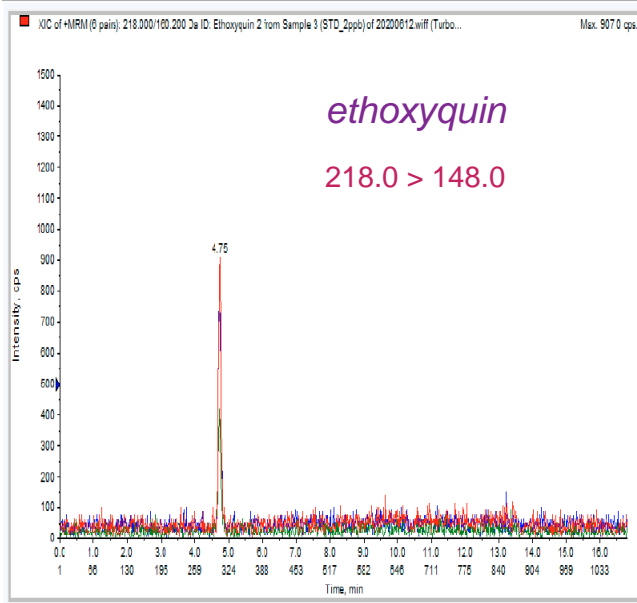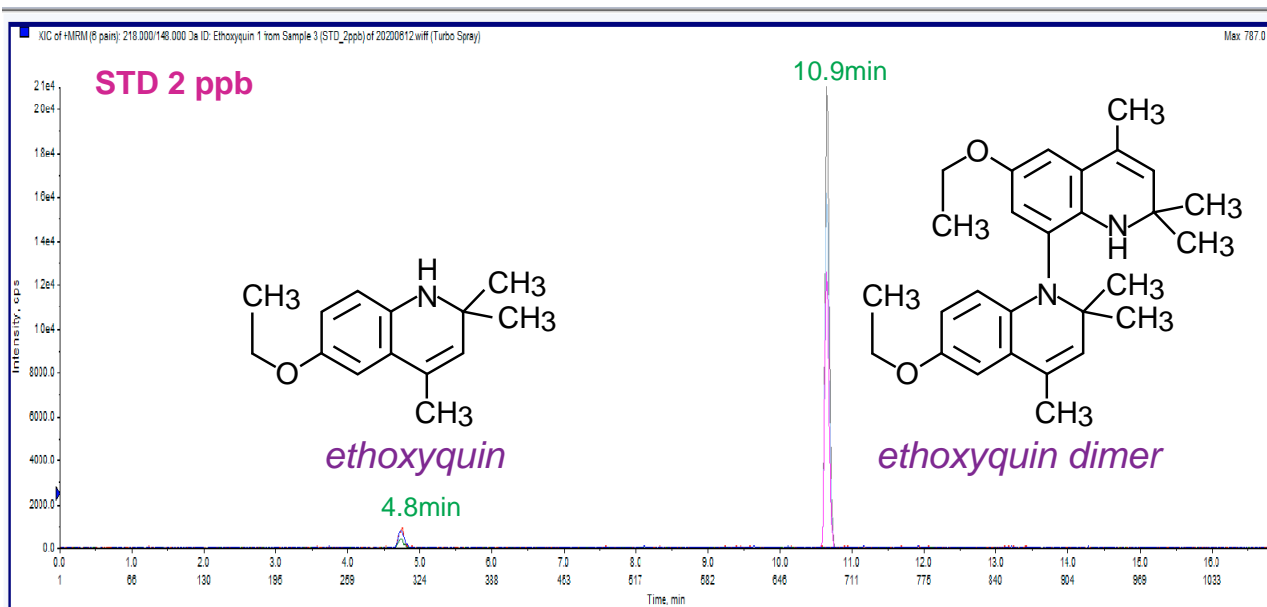

Cadenza CD-C18 HT, 150 x 2 mm

A: 10mM ammonium formate

B: acetonitrile

60%B (0-1min), 60-100%B (1-8min), 100%B (8-11min),

100-60%B (11-11.1min), 60%B (11.1-16min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)

5µL (2 ppb)

Courtesy of PEAKMAN SP, KOREA