

Unison UK-C18 HT


100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Antibiotics (9)

抗生物質の LC-MS/MS 分析 (9)

Antibiotics 107mix 10min method

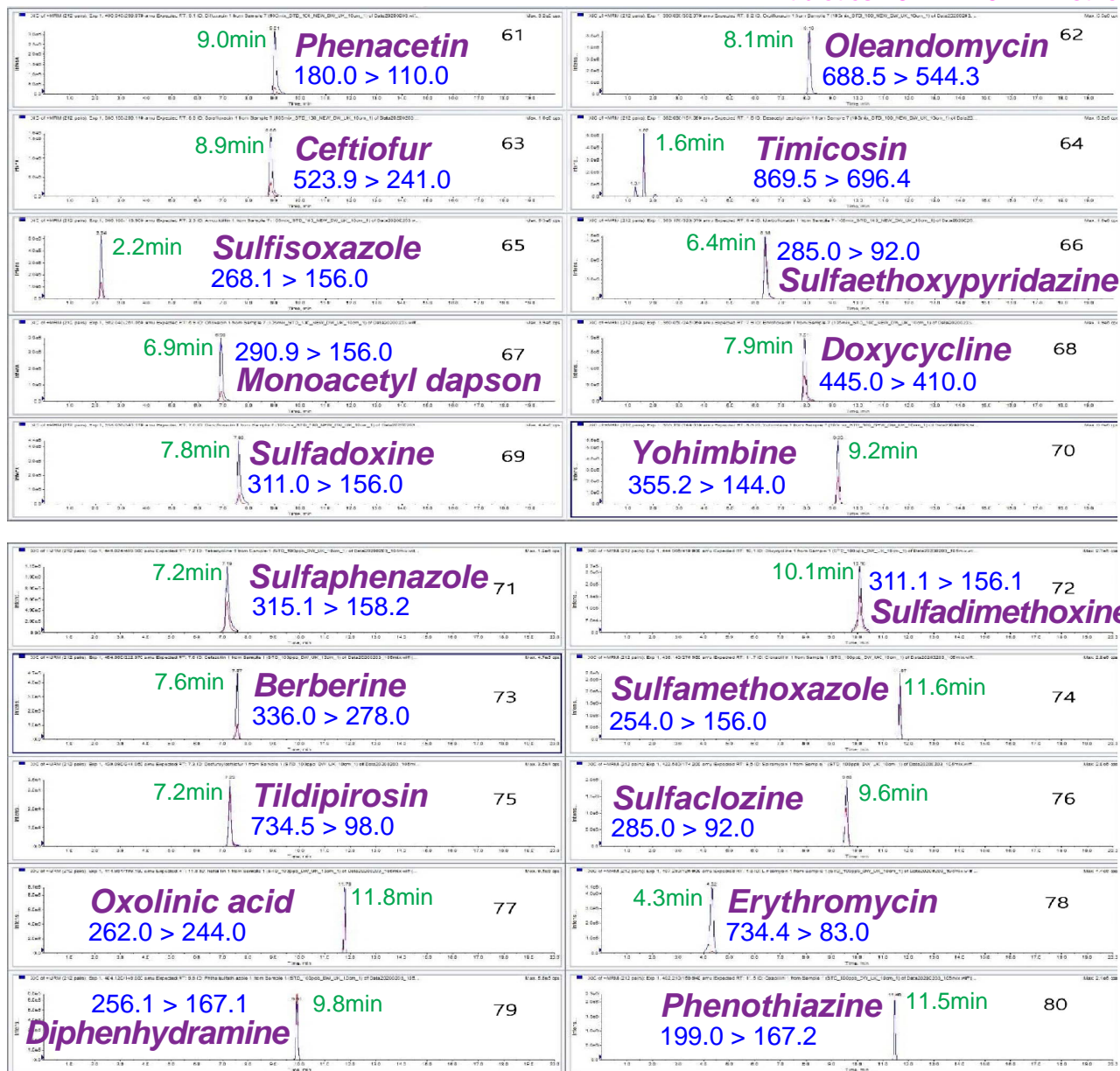

Unison UK-C18 HT, 100 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid in acetonitrile

10%B (0-0.5min), 10-15%B (0.5-3min), 15-95%B (3-6min),

95%B (6-8min), 95-10%B (8-8.1min), 10%B (8.1-10min)

0.4mL/min, 40°C, Triple Quad. MS (ESI, Positive/Negative)

5µL

Courtesy of PEAKMAN SP, KOREA