







Unison UK-C18 HT

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Pesticides (5)

殺虫剤の LC-MS/MS 分析 (5)

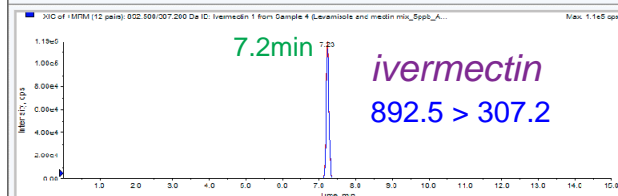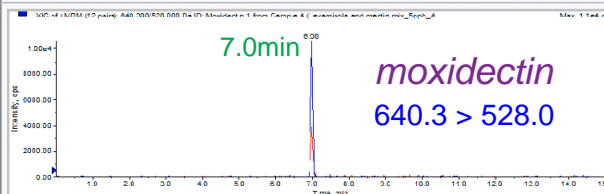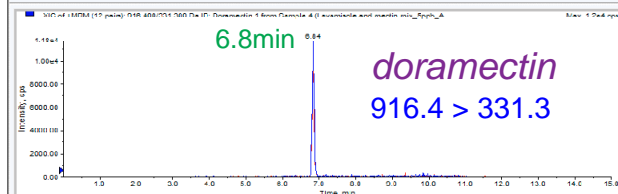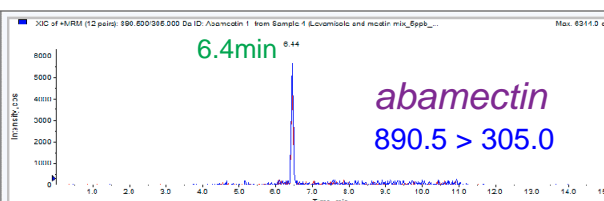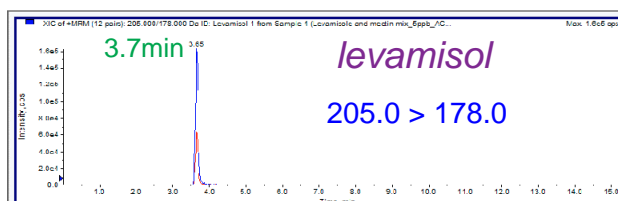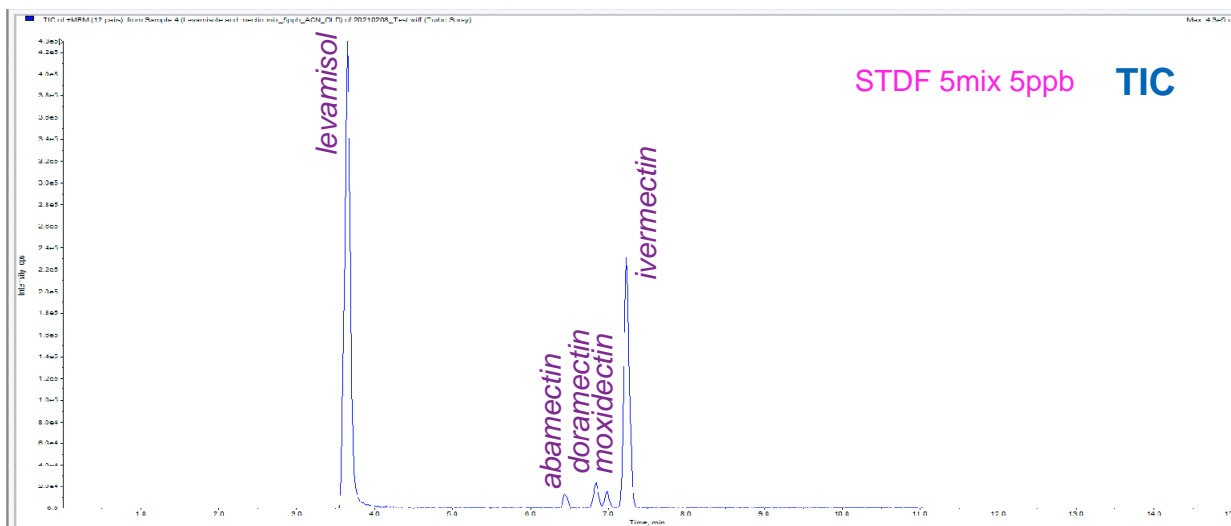

Unison UK-C18 HT, 100 x 2 mm

A: 0.1% formic acid + 20mM ammonium formate in water

B: 0.2% formic acid in methanol

10%B (0-1min), 10-95%B (1-2min), 95%B (2-9min),

95-10%B (9-9.1min), 10%B (9.1-15min)

0.25mL/min, 40°C, Triple Quad. MS (ESI, Positive)

2µL (5 ppb)

Courtesy of PEAKMAN SP, KOREA