


Cadenza CD-C18 HT

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Antiprotozoal agents

抗原虫薬の LC-MS/MS 分析

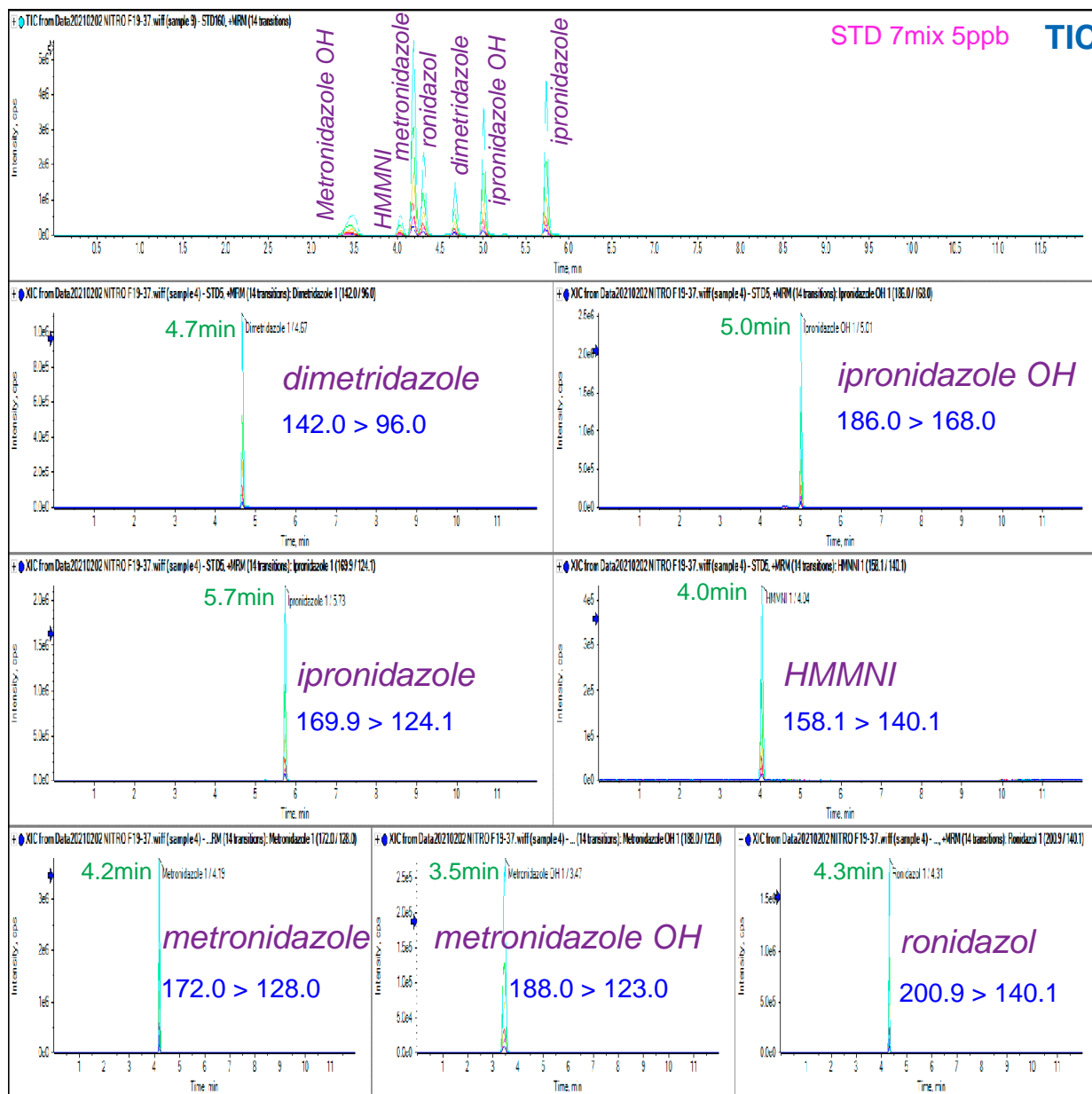

Cadenza CD-C18 HT, 150 x 2 mm
 A: 5mM ammonium acetate (adjust pH 5.0)
 B: acetonitrile
 5%B (0-1min), 5-95%B (1-6min), 95%B (6-8min),
 95-5%B (8-8.1min), 5%B (8.1-12min)
 0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)
 2µL (5 ppb)

Courtesy of PEAKMAN SP, KOREA