



**Cadenza CD-C18**

150 x 2 mm

**HPLC COLUMN APPLICATION**

**LC-MS/MS: Sulfonamides**

サルファ剤の LC-MS/MS 分析

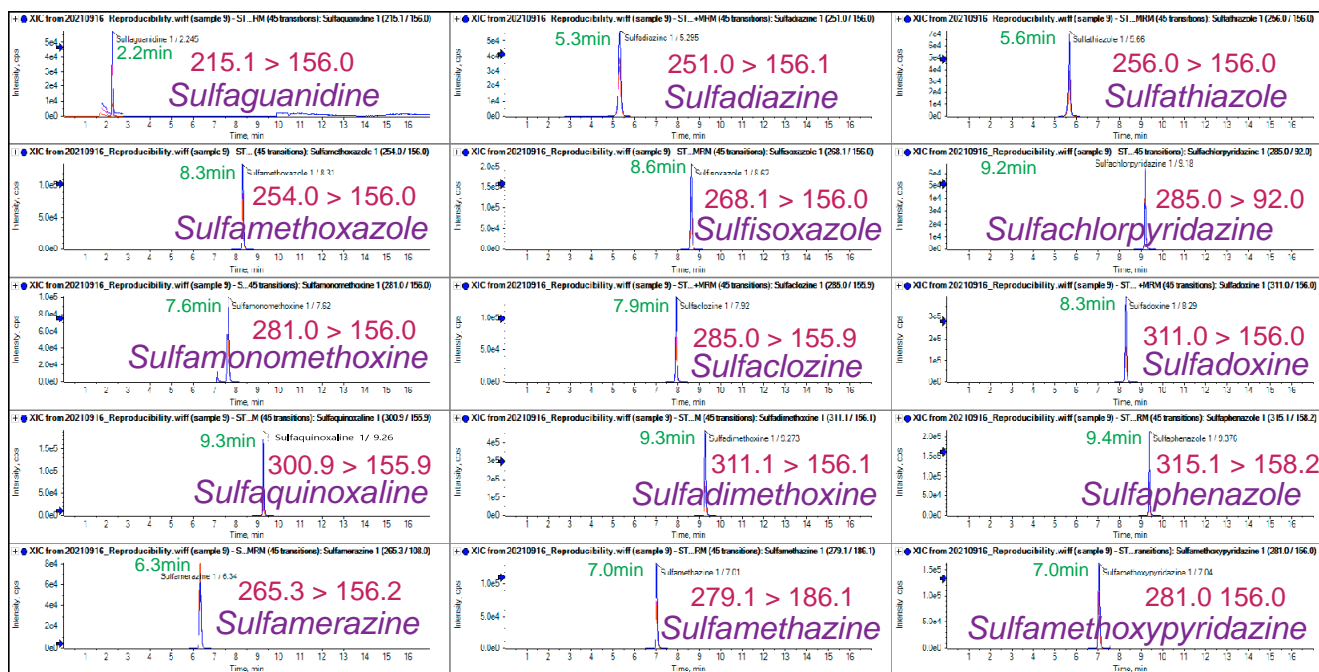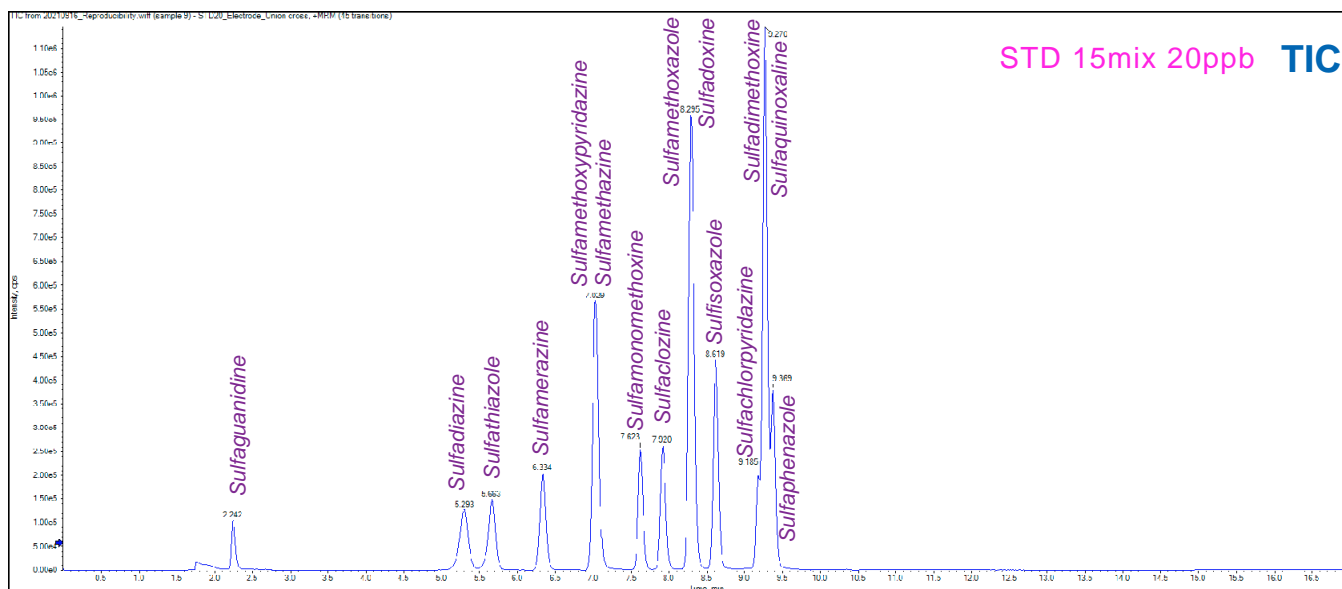

Cadenza CD-C18, 150 x 2 mm  
 A: 0.1% formic acid in water  
 B: 0.1% formic acid in acetonitrile  
 10%B (0-1min), 10-40%B (1-6min), 40-95%B (6-10min),  
 95%B (10-12min), 95-10%B (12-12.1min), 10%B (12.1-17min)  
 0.2mL/min, 40°C, Triple Quad. MS (ESI, Positive)  
 2µL (20 ppb)

Courtesy of PEAKMAN SP, KOREA