

Cadenza CD-C18 HT


150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Bacitracin

バシトラシンの LC-MS/MS 分析

STD 5, 10, 20, 40, 80 ppb Overlay

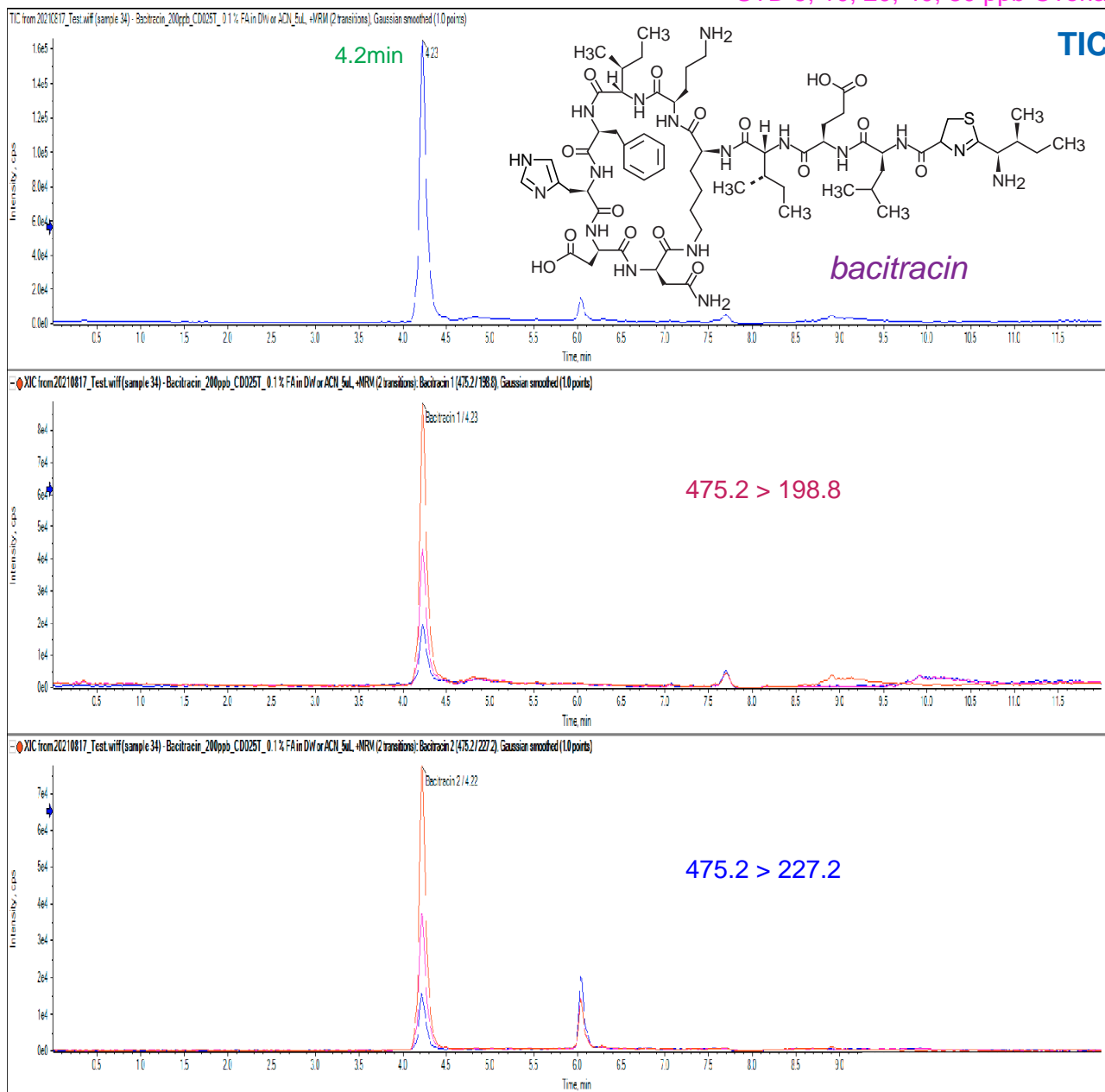

Cadenza CD-C18 HT, 150 x 2 mm  
 A: 0.1% formic acid in water  
 B: 0.1% formic acid in acetonitrile  
 10%B (0-1min), 10-90%B (1-5min), 90%B (5-7min),  
 90-10%B (7-7.1min), 10%B (7.1-12min)  
 0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)  
 5µL (5, 10, 20, 40, 80 ppb)

Courtesy of PEAKMAN SP, KOREA