


Cadenza CD-C18 HT

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Cephem antibiotics

セフェム系抗生物質の LC-MS/MS 分析

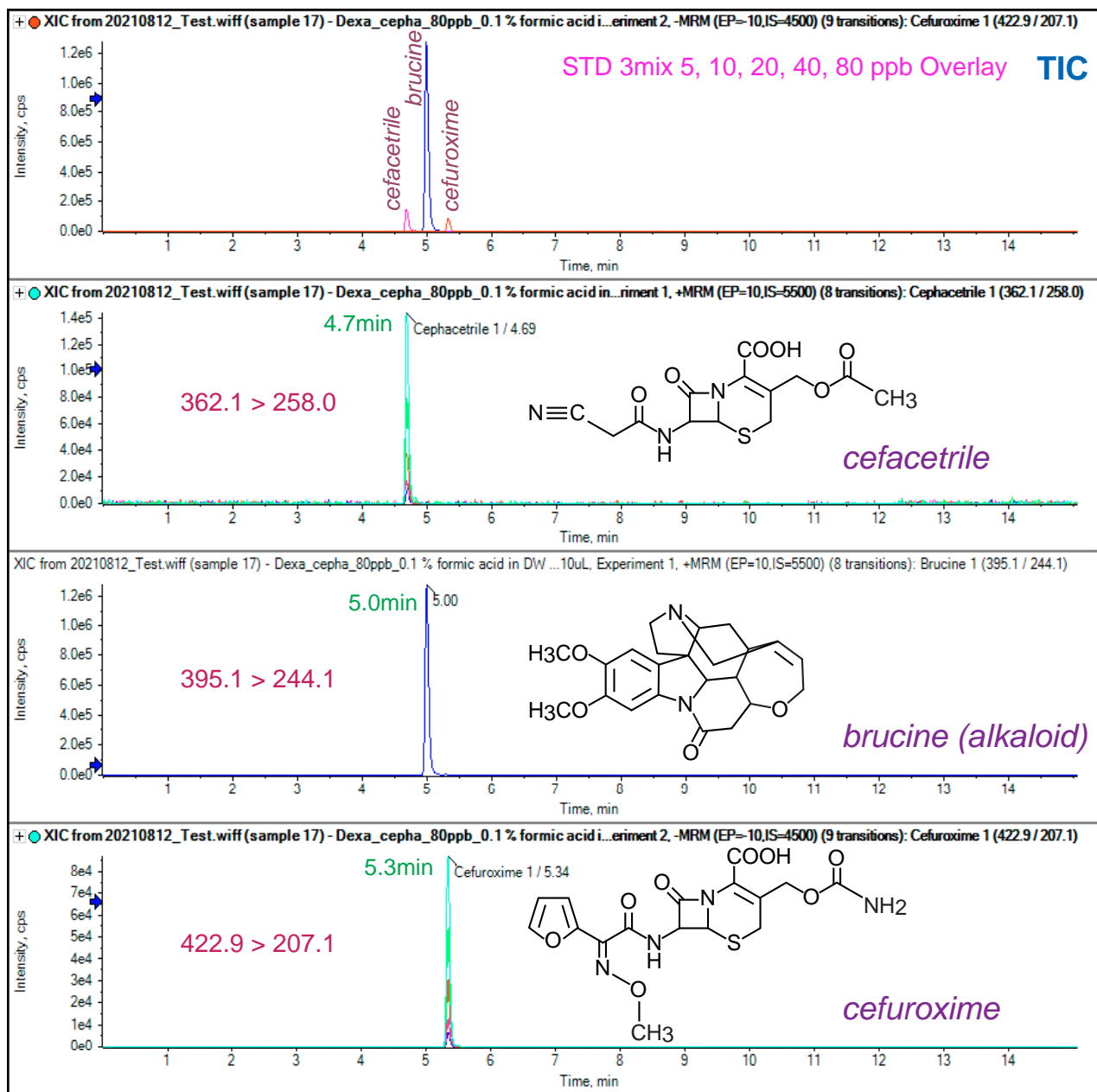

Cadenza CD-C18 HT, 150 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid in acetonitrile

10%B (0-1min), 10-90%B (1-7min), 90%B (7-10min),

90-10%B (10-10.1min), 10%B (10.1-15min)

0.25mL/min, 40°C, Triple Quad. MS (ESI, Negative)

10µL (5, 10, 20, 40, 80 ppb)

Courtesy of PEAKMAN SP, KOREA