


Cadenza CD-C18

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Penicillin antibiotics

ペニシリン系抗菌薬のLC-MS/MS分析

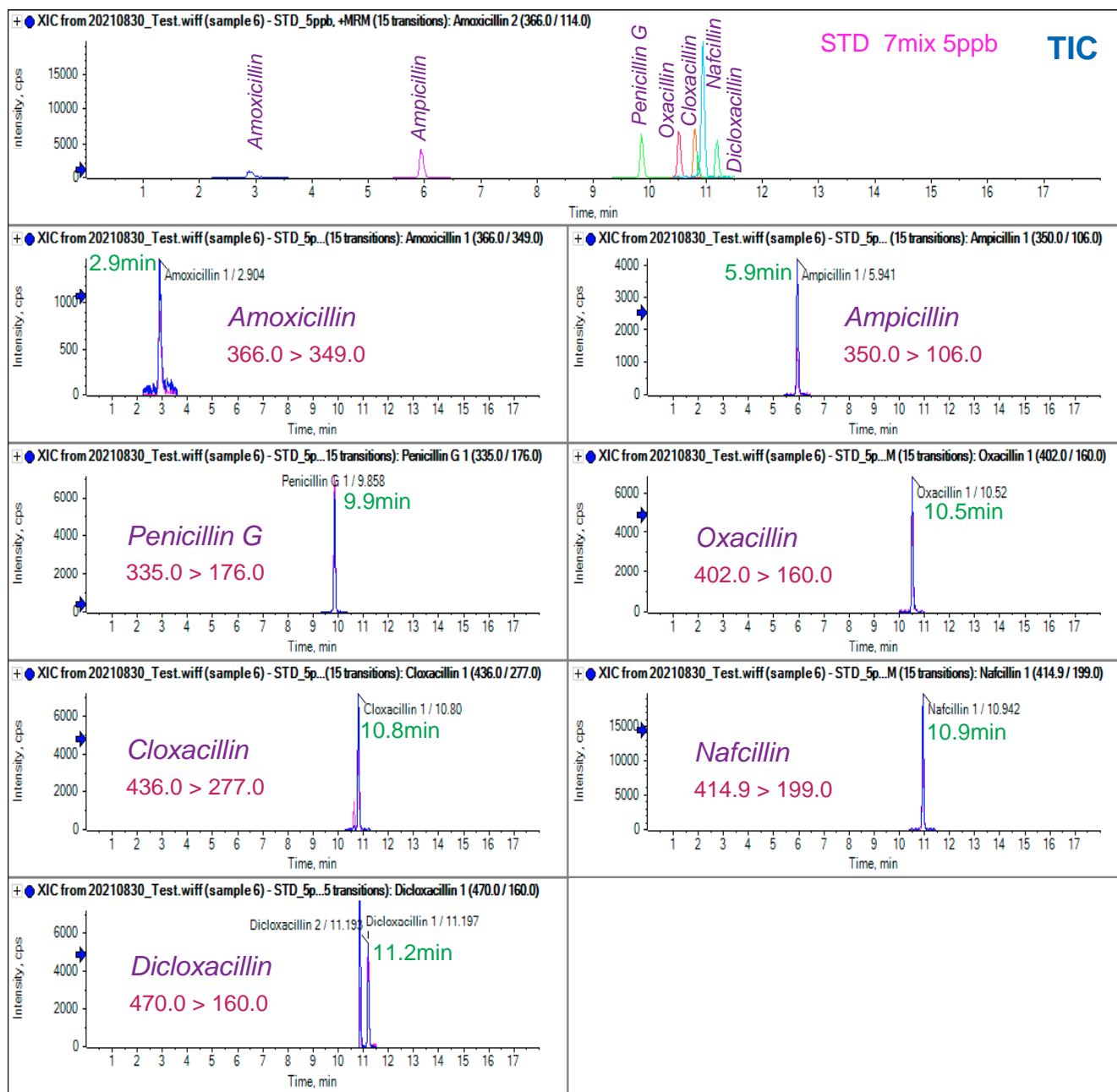

Cadenza CD-C18, 150 x 2 mm

A: 0.1% formic acid in water

B: 0.1% formic acid in acetonitrile

10%B (0-1min), 10-40%B (1-6min), 40-95%B (6-10min),

95%B (10-12min), 95-10%B (12-12.1min), 10%B (12.1-18min)

0.2mL/min, 40°C, Triple Quad. MS (ESI, Positive)

5µL (5 ppb)

Courtesy of PEAKMAN SP, KOREA