







Unison UK-C8HT

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Polyether antibiotics

ポリエーテル系抗菌薬の LC-MS/MS 分析

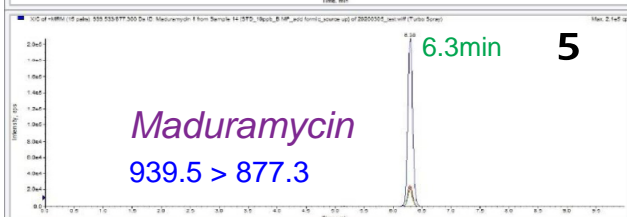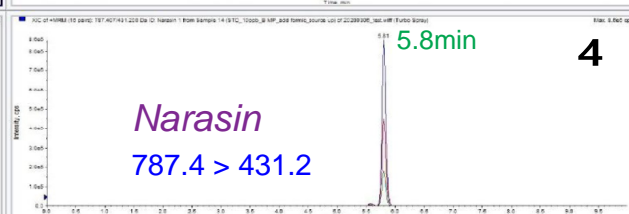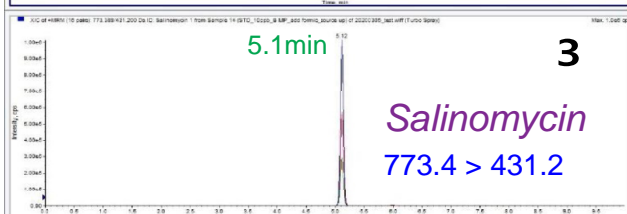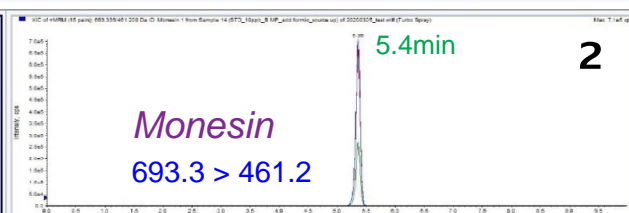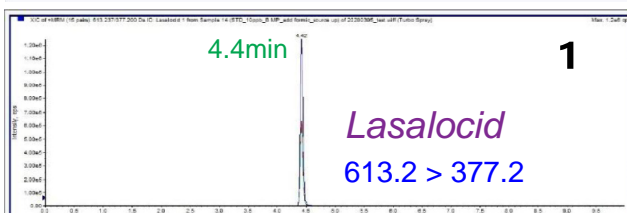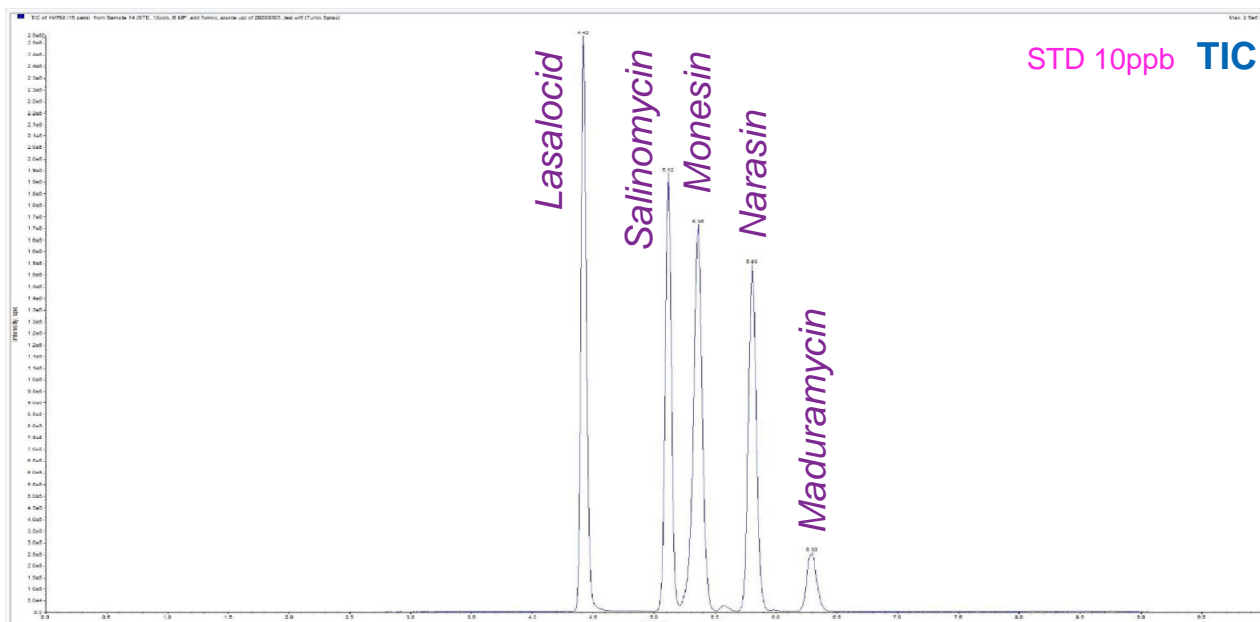

Unison UK-C8 HT, 100 x 2 mm
 A: 0.1% formic acid in water
 B: 0.1% formic acid in acetonitrile
 60%B (0-0.5min), 60-95%B (0.5-2.5min), 95%B (2.5-7.5min),
 95-60%B (7.5-7.6min), 60%B (7.6-10min)
 0.4mL/min, 35°C, Triple Quad. MS (ESI, Positive)
 2µL (10 ppb)

Courtesy of PEAKMAN SP, KOREA