


Unison UK-C18

50 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Estrogenic substances

エストロゲン様物質の LC-MS/MS 分析

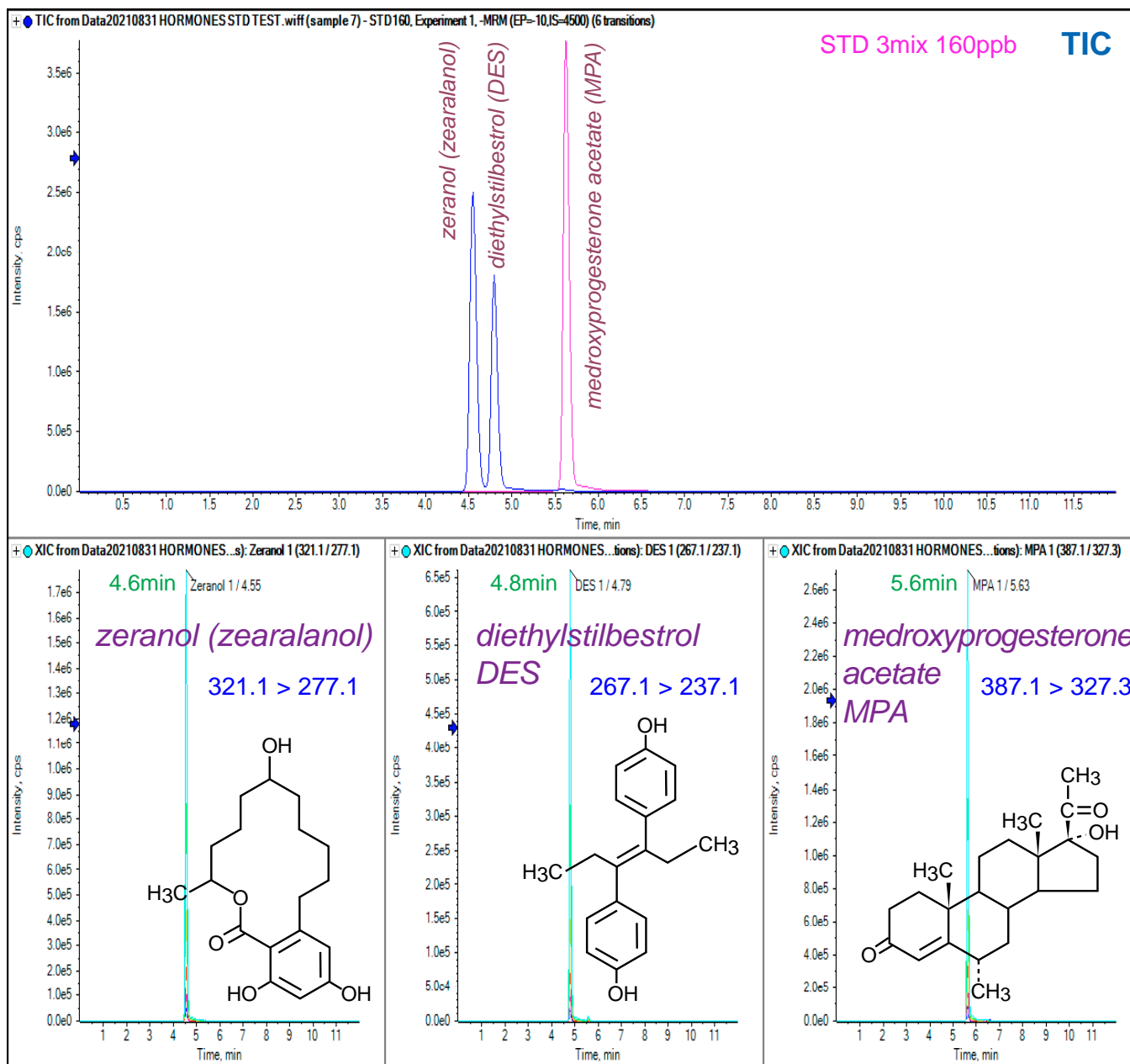

Unison UK-C18, 50 x 2 mm

A: 10mM ammonium acetate

B: methanol

50%B (0-0.5min), 50-90%B (0.5-4min), 90%B (4-7.5min),

90-50%B (7.5-7.6min), 50%B (7.6-12min)

0.2mL/min, 40°C, Triple Quad. MS (ESI, Negative)

2µL (160 ppb)

Courtesy of PEAKMAN SP, KOREA