




Cadenza CD-C18 HT

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Herbicides

除草剤の LC-MS/MS 分析

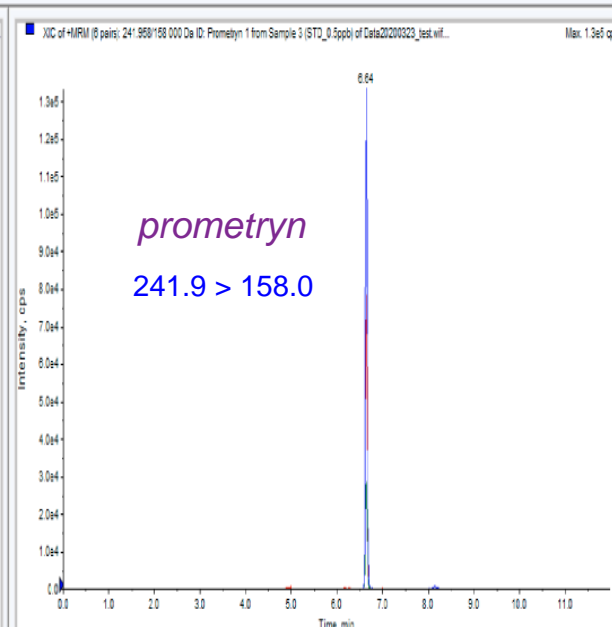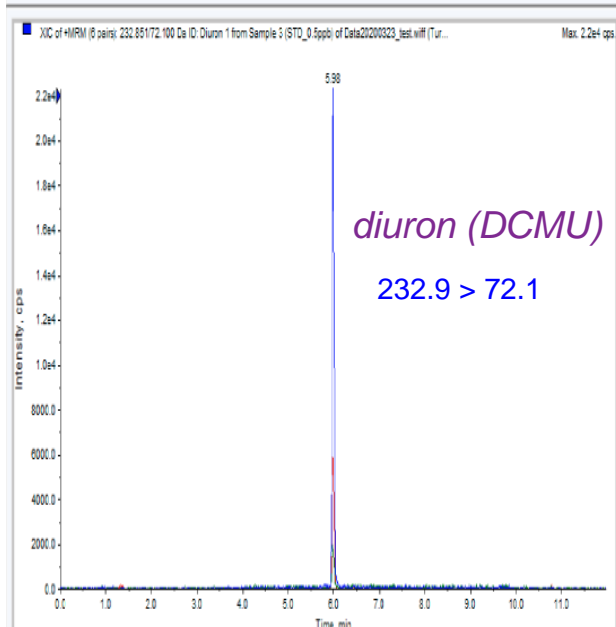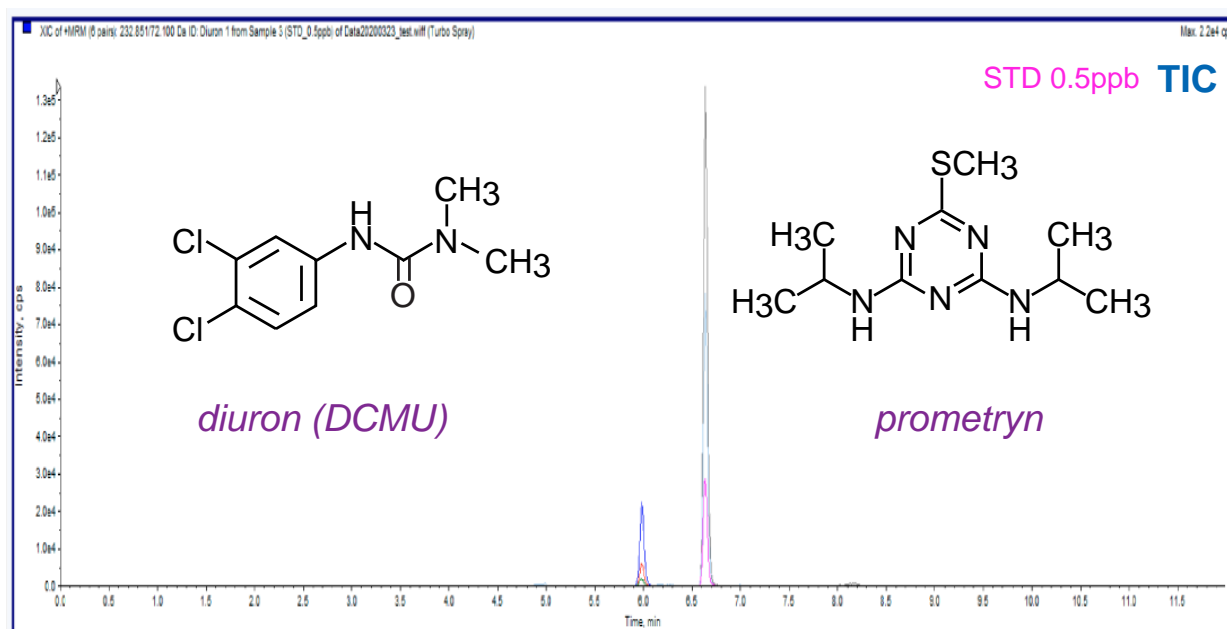

Cadenza CD-C18 HT, 150 x 2 mm

A: 5mM ammonium formate

B: 5mM ammonium formate in methanol

15%B (0-1min), 15-40%B (1-1.5min), 40-95%B (1.5-5min),

95%B (5-8min), 95-15%B (8-8.1min), 15%B (8.1-12min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)

5µL (0.5 ppb)

Courtesy of PEAKMAN SP, KOREA