



Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Mirtazapine (Antidepressant)

ミルタザピン (抗うつ薬)

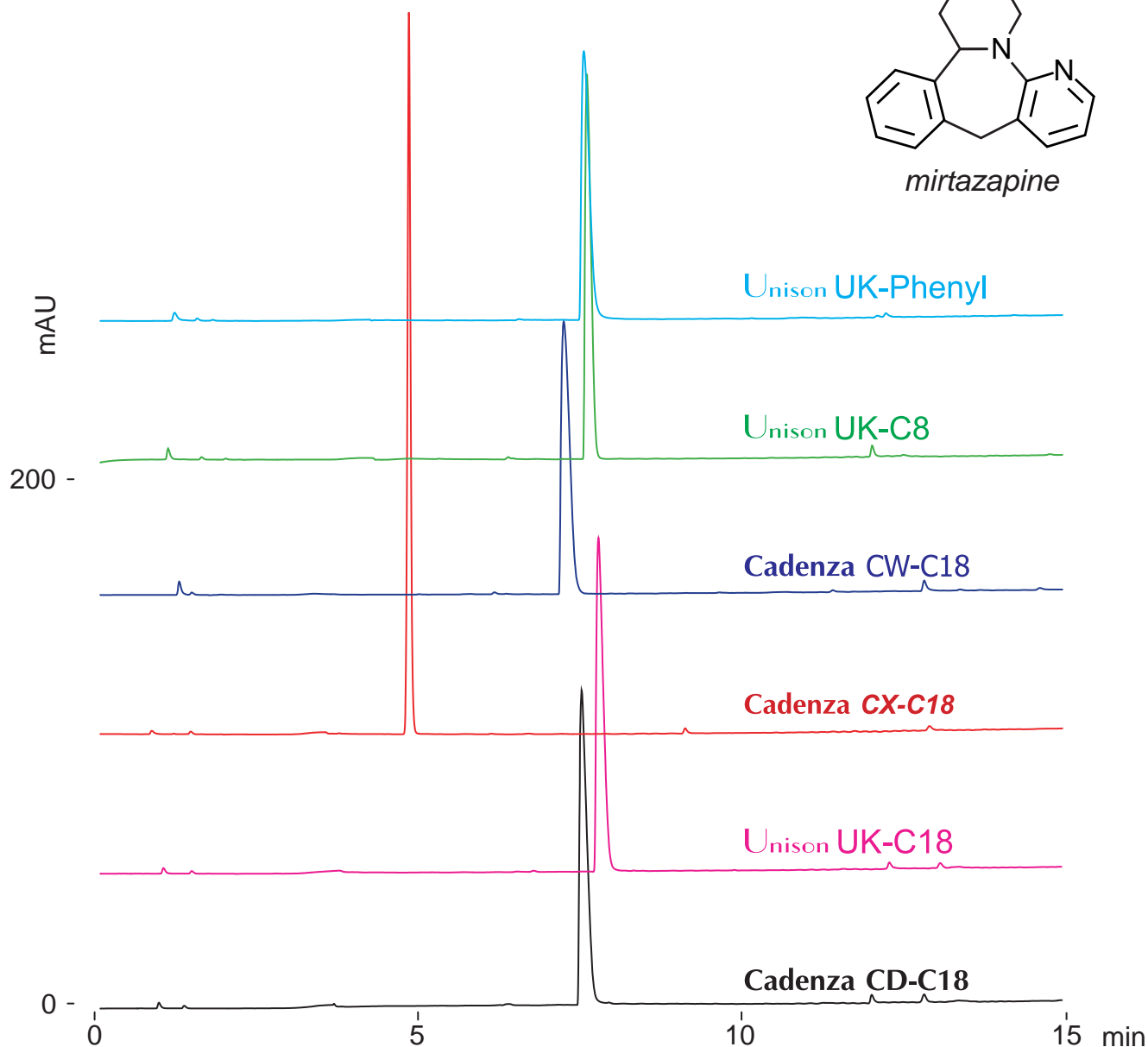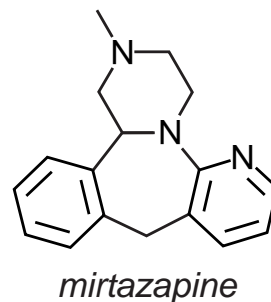

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 50 / 50 / 0.5

0%B (0-1min), 0-90%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.25mg/mL, water / formic acid = 100 / 1)