


Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Mirtazapine (Antidepressant)

ミルタザピン (抗うつ薬)

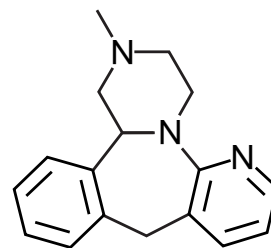

mirtazapine

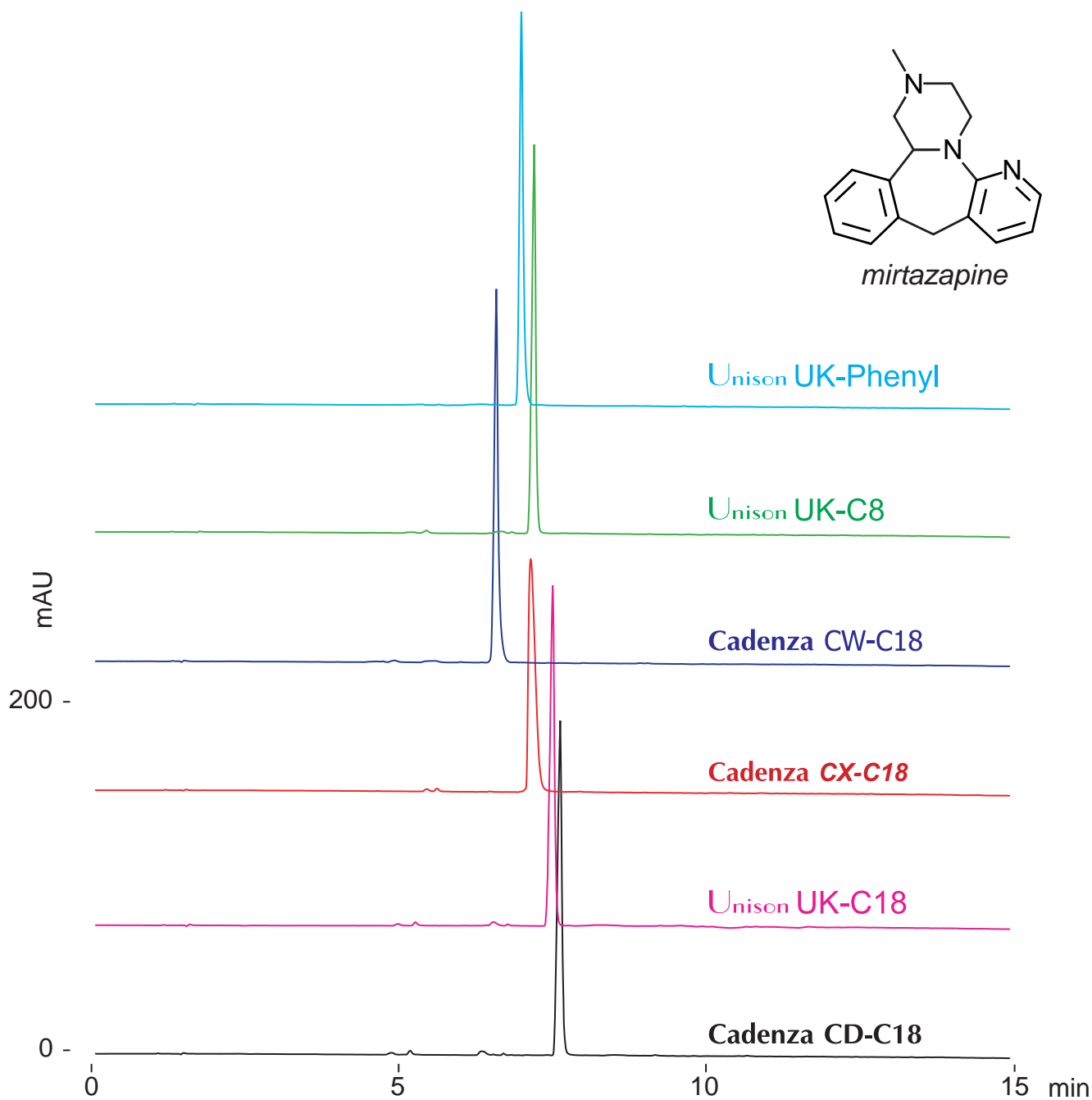

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.2mg/mL, 100mM ammonium formate / acetonitrile = 80 / 20)