



Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Venlafaxine (Antidepressant)

ベンラファキシン (抗うつ薬)

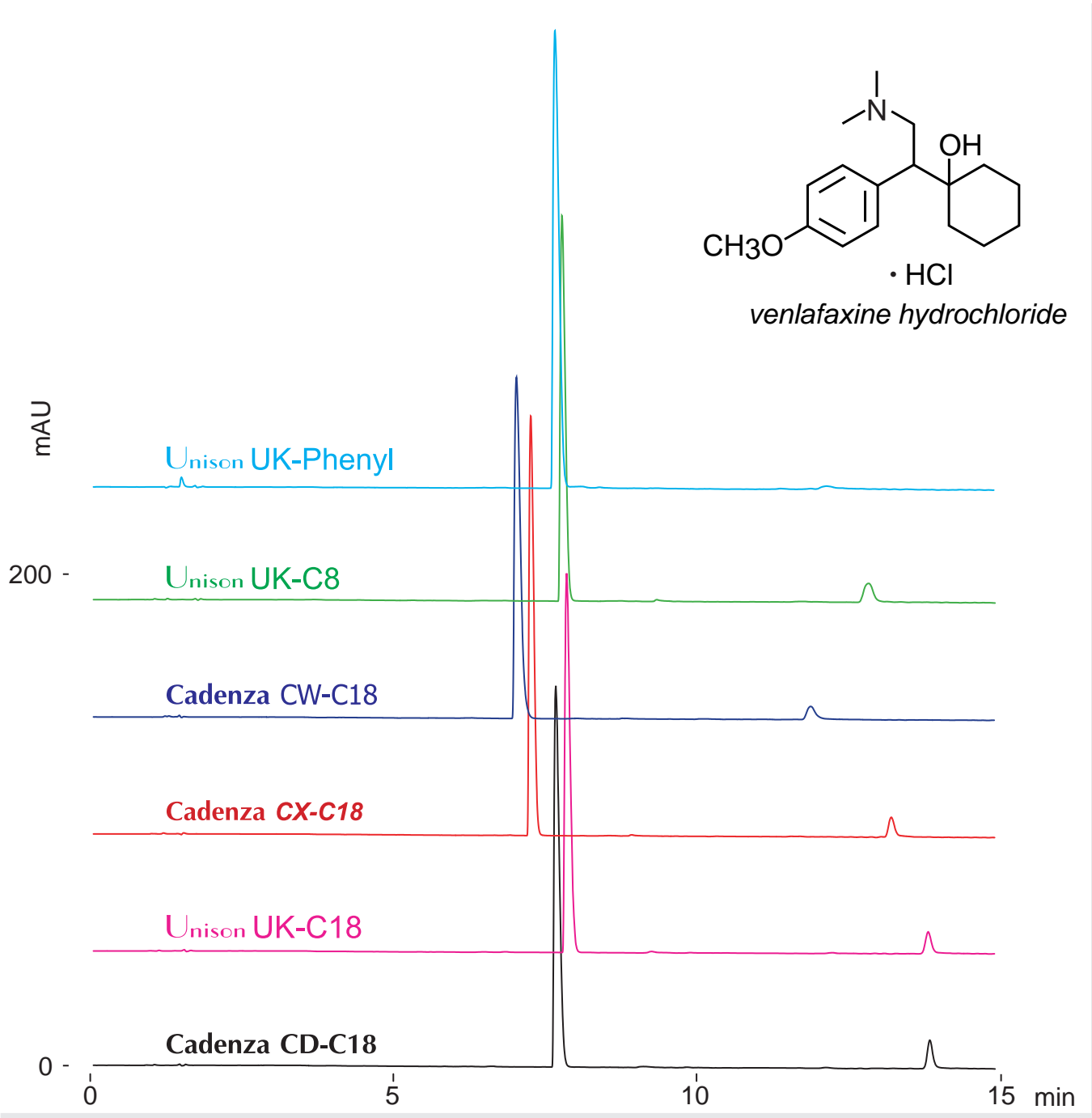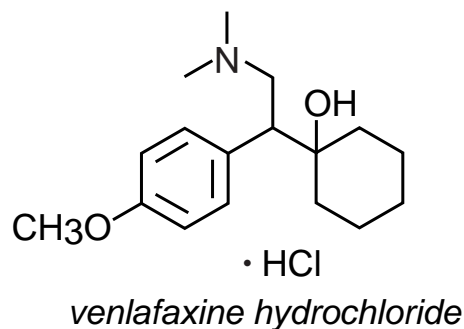

100 x 4.6 mm
 A: 100mM ammonium formate / acetonitrile = 90 / 10
 B: 100mM ammonium formate / acetonitrile = 10 / 90
 0%B (0-1min), 0-75%B (1-15min)
 1.0mL/min (9MPa), 37°C, 270nm
 10µL (1mg/mL, 100mM ammonium formate)