



Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Rasagiline (Antiparkinsonian drug)

ラサジリン (パーキンソン病治療薬)

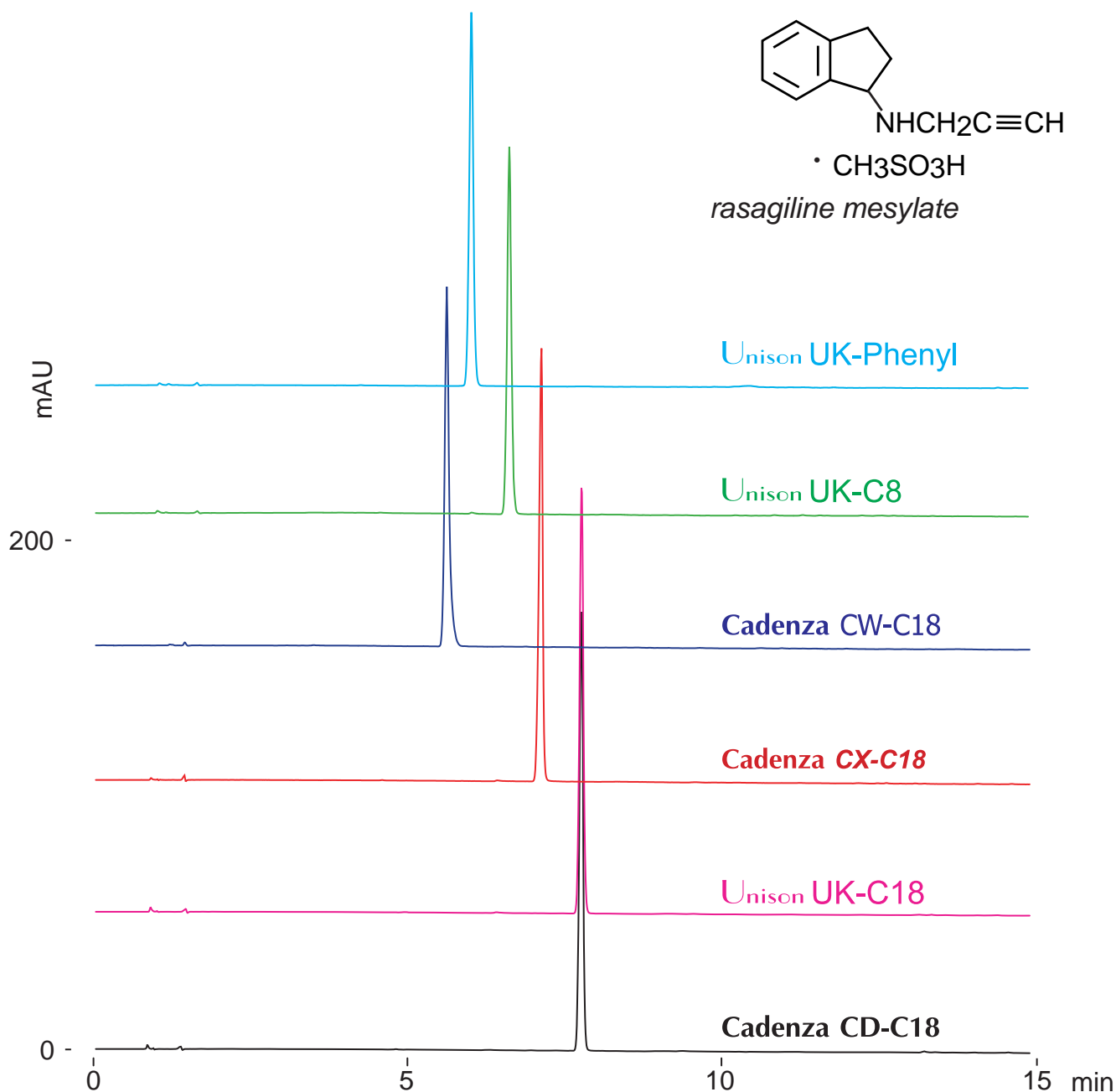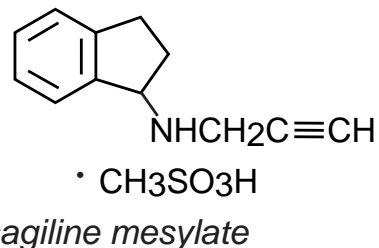

100 x 4.6 mm

A: 100mM ammonium acetate / acetonitrile = 90 / 10

B: 100mM ammonium acetate / acetonitrile = 10 / 90

20-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (1mg/mL, 100mM ammonium acetate)