



**Cadenza CX-C18**

150 x 3 mm

**HPLC COLUMN APPLICATION**

**Rasagiline (Antiparkinsonian drug)**

ラサジリン (パーキンソン病治療薬)

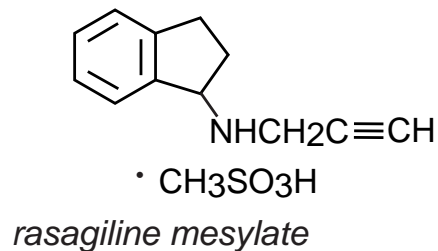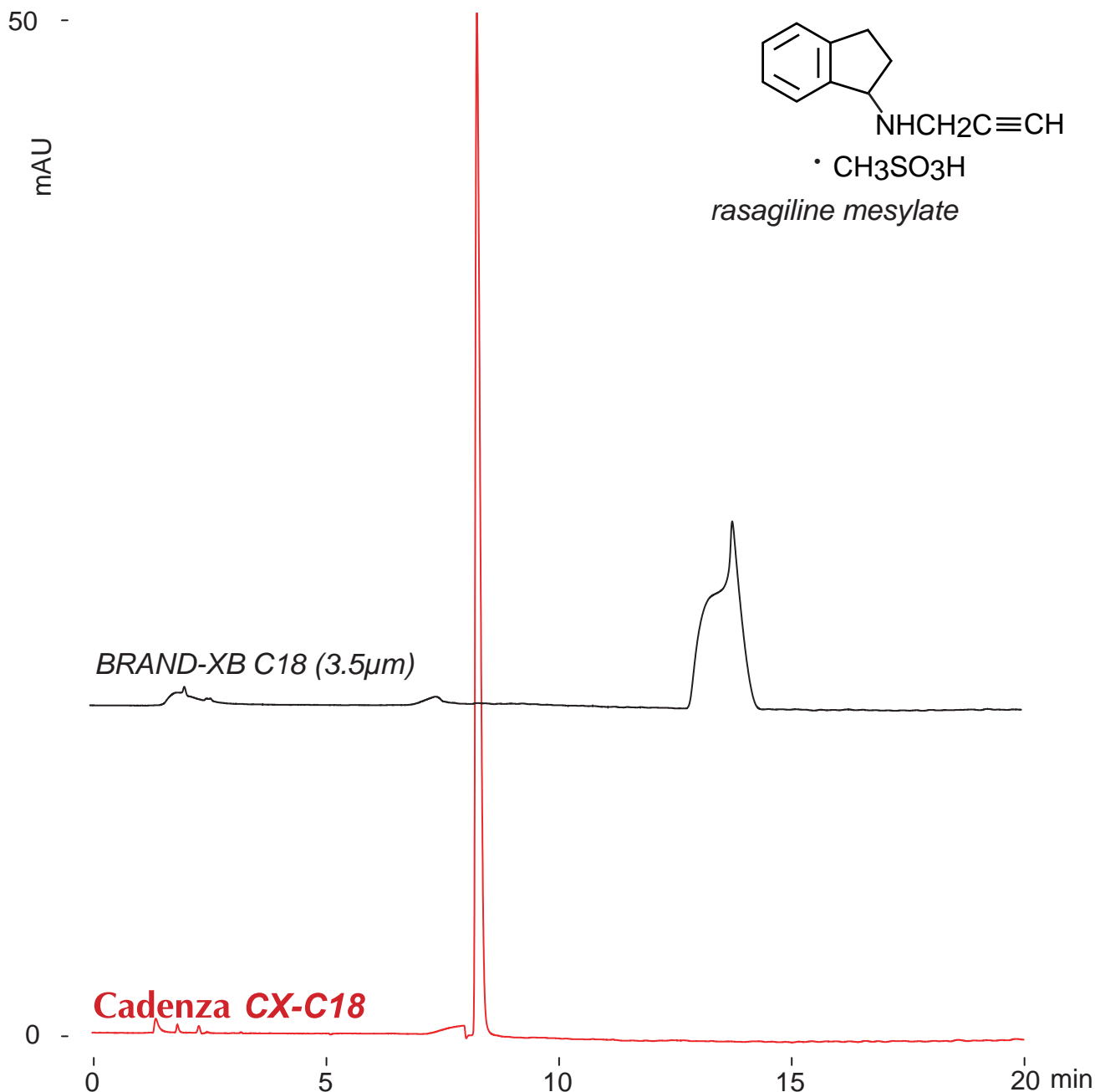

150 x 3 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

0%B (0-1min), 0-25%B (1-20min)

0.4mL/min (9MPa), 37°C, 260nm

2µL (1mg/mL, water / formic acid = 100 / 1)