


DACAPO DX-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION


Rasagiline (Antiparkinsonian drug)

ラサジリン (パーキンソン病治療薬)

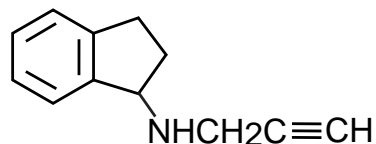

• CH₃SO₃H

rasagiline mesylate

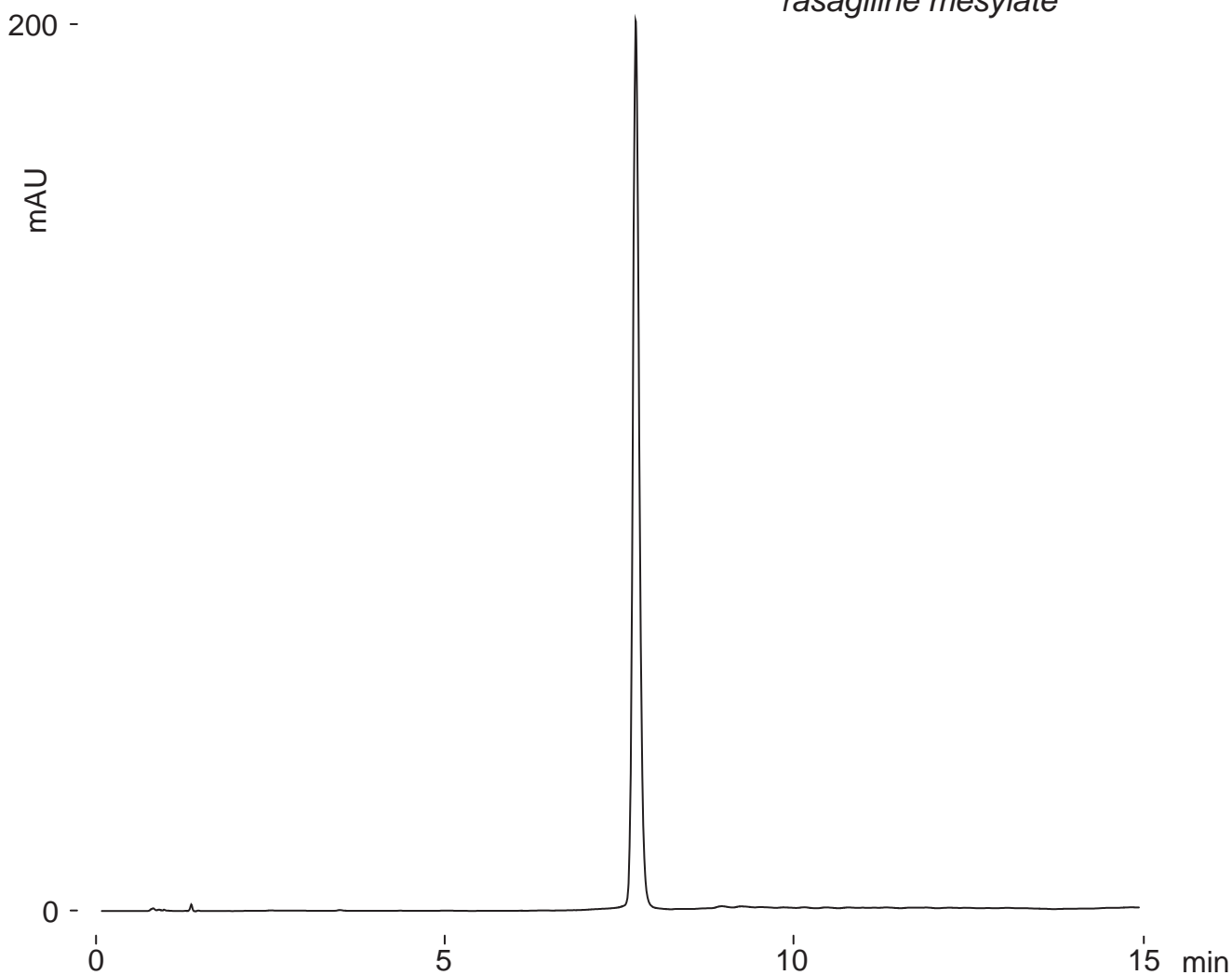

Dacapo DX-C18, 100 x 4.6 mm
A: acetonitrile / 10mM NH₄OH = 10 / 90
B: acetonitrile / 10mM NH₄OH = 50 / 50
0-100%B (0-15min)
1mL/min (12MPa), 37°C, 260nm
10μL (1mg/mL, 100mM NH₄OH)