



Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Aztreonam (Monobactam antibiotic)

アズトレオナム (モノバクタム系抗菌薬)

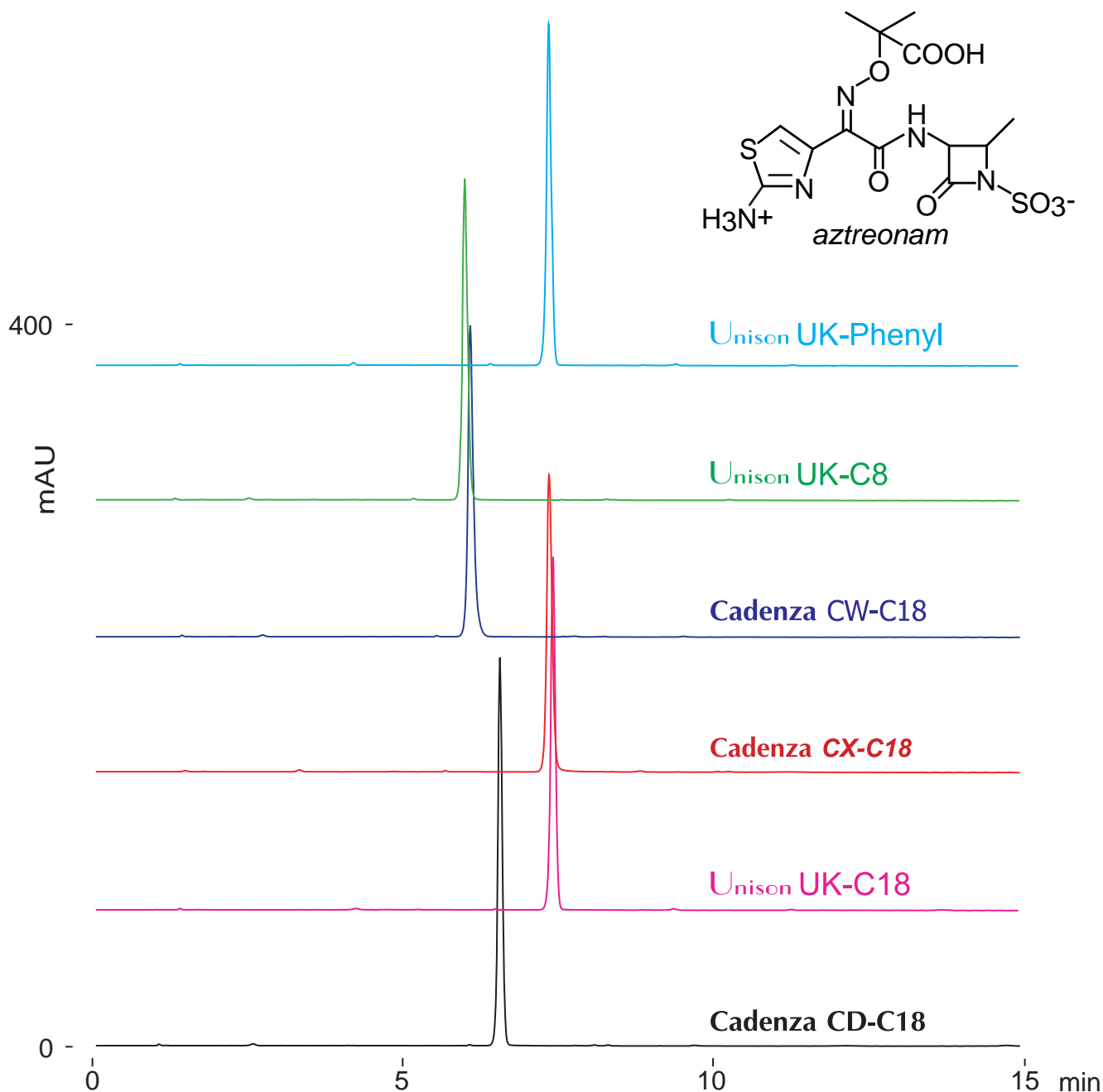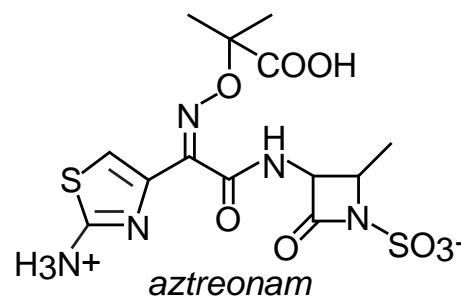

100 x 4.6 mm

A: 100mM ammonium acetate

B: 100mM ammonium acetate / acetonitrile = 70 / 30

0%B (0-1min), 0-60%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10μL (0.15mg/mL, 100mM ammonium acetate)