
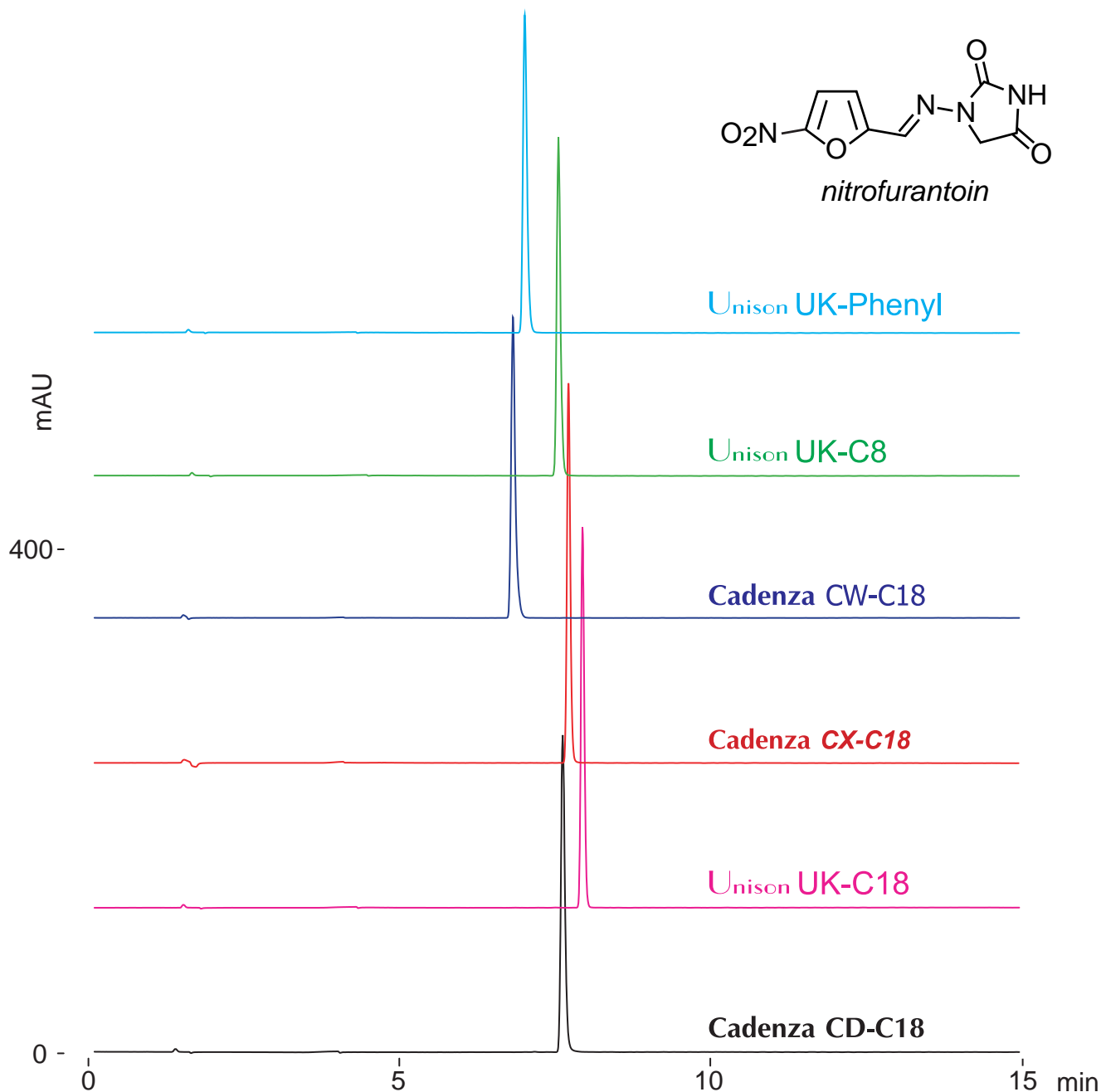


Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Nitrofurantoin (5-Nitrofuran antibacterial)
 ニトロフラントイン (5-ニトロフラン系抗菌薬)

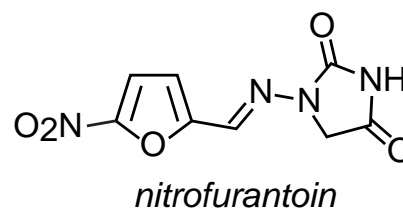

100 x 4.6 mm
 A: water / formic acid = 100 / 0.1
 B: water / acetonitrile / formic acid = 10 / 90 / 0.1
 0%B (0-1min), 0-75%B (1-15min)
 1.0mL/min (9MPa), 37°C, 360nm
 10µL (0.05mg/mL, water / acetonitrile / formic acid = 50 / 50 / 1)