
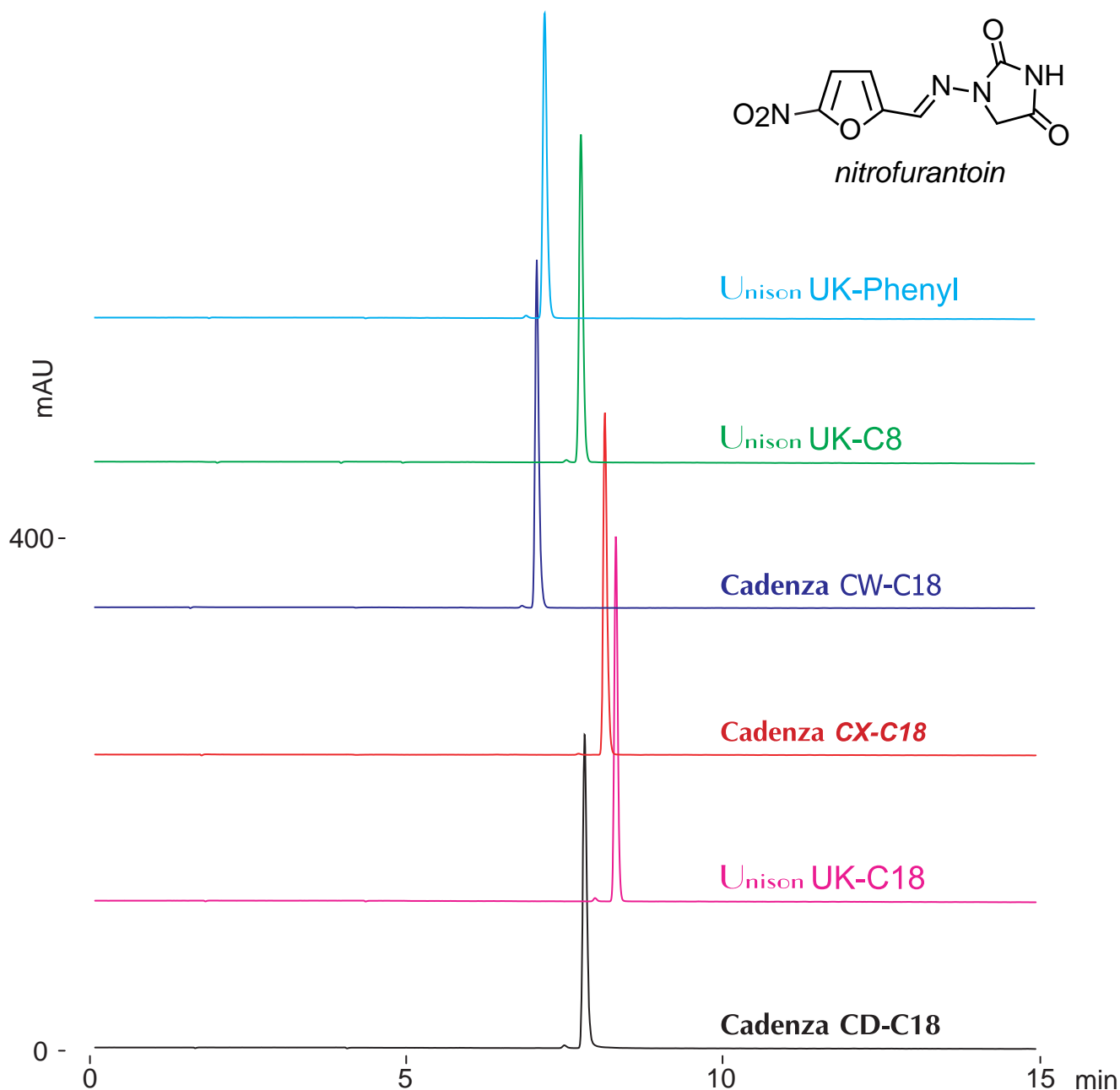


Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Nitrofurantoin (5-Nitrofuran antibacterial)
 ニトロフラントイン (5-ニトロフラン系抗菌薬)

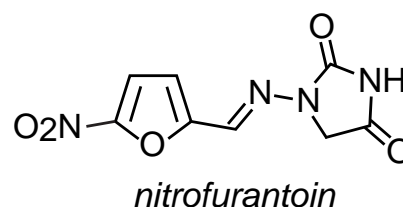

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0%B (0-1min), 0-65%B (1-15min)

1.0mL/min (9MPa), 37°C, 360nm

10µL (0.05mg/mL, 100mM ammonium formate / acetonitrile = 40/ 60)