



Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Phlorizin (Dihydrochalcone glycoside)

フロリジン (ジヒドロカルコン配糖体)

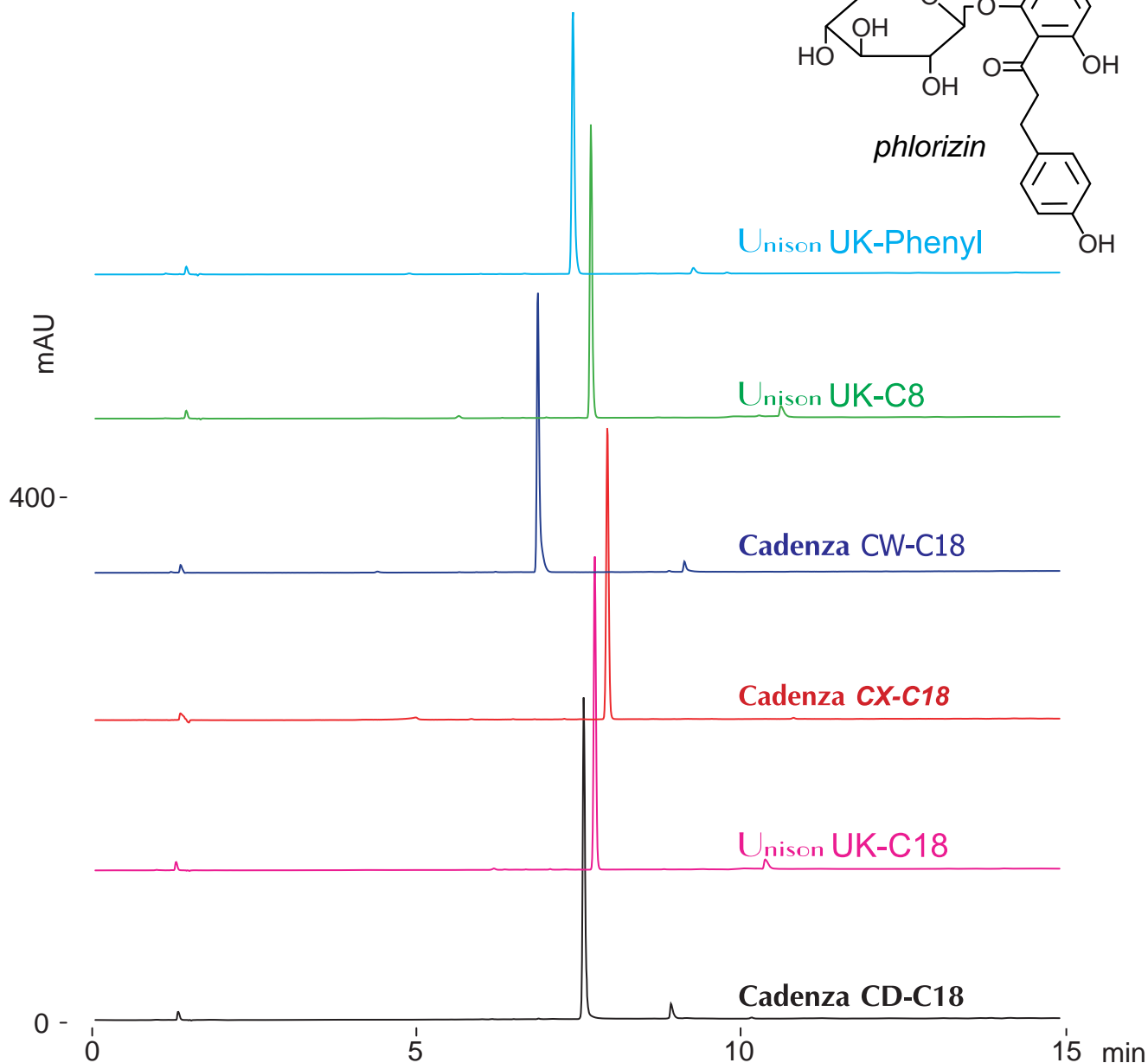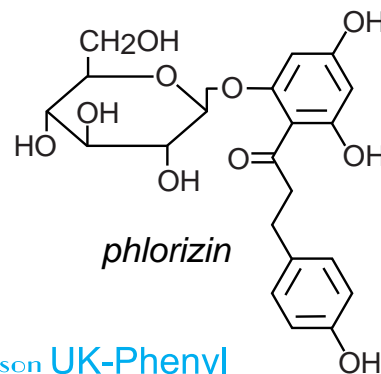

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0%B (0-1min), 0-65%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.15mg/mL, water / acetonitrile / formic acid = 80 / 20 / 1)