
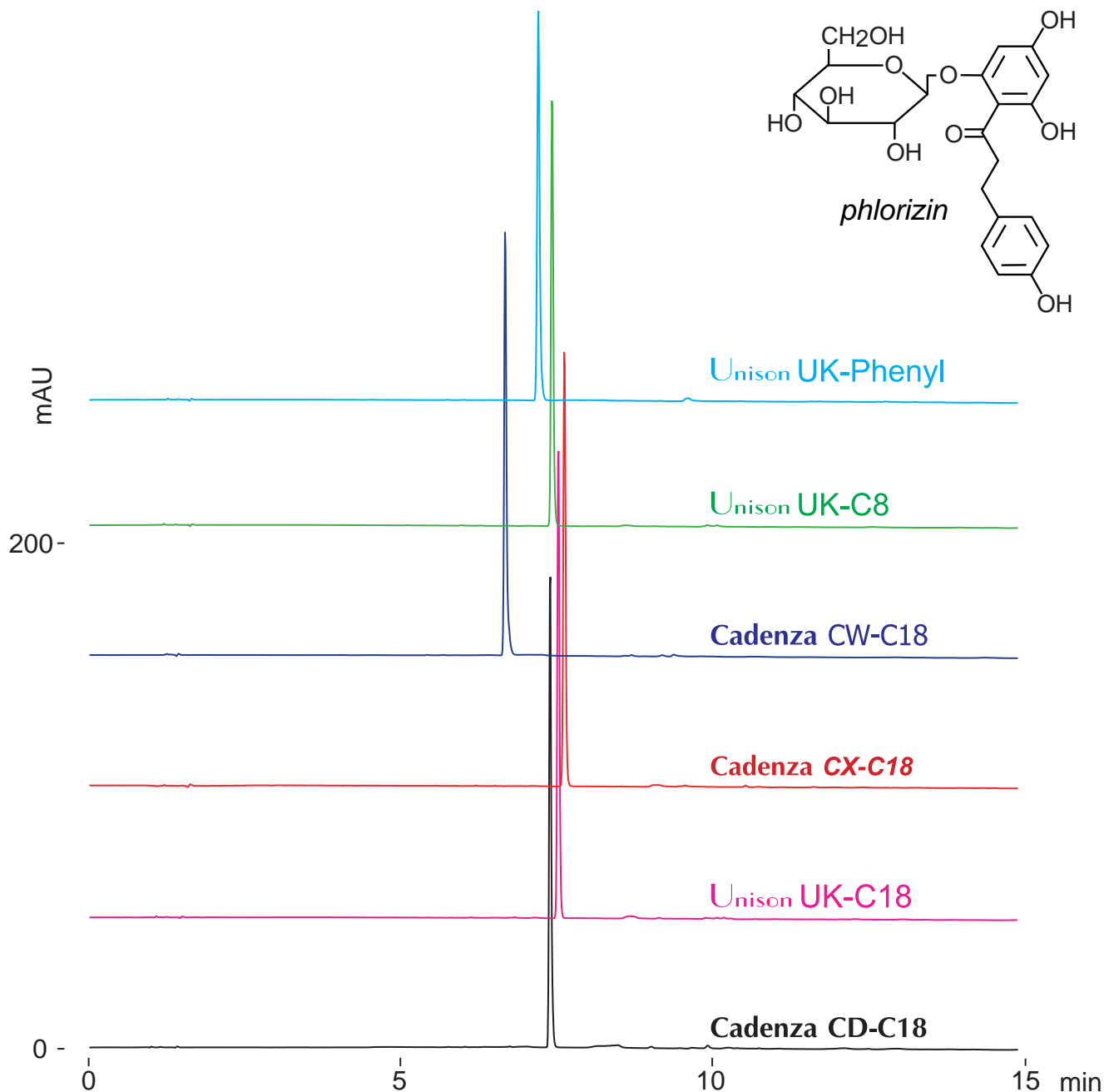


Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

**Phlorizin (Dihydrochalcone glycoside)**  
 フロリジン (ジヒドロカルコン配糖体)

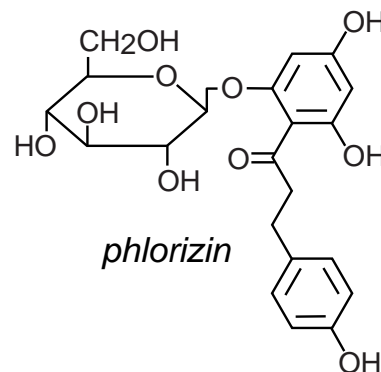

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0%B (0-1min), 0-65%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.15mg/mL, 100mM ammonium formate / acetonitrile = 80 / 20)