



Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Flonicamid (Pyridine carboxamide pesticides)

フロニカミド (ピリジンカルボキサミド系殺虫剤)

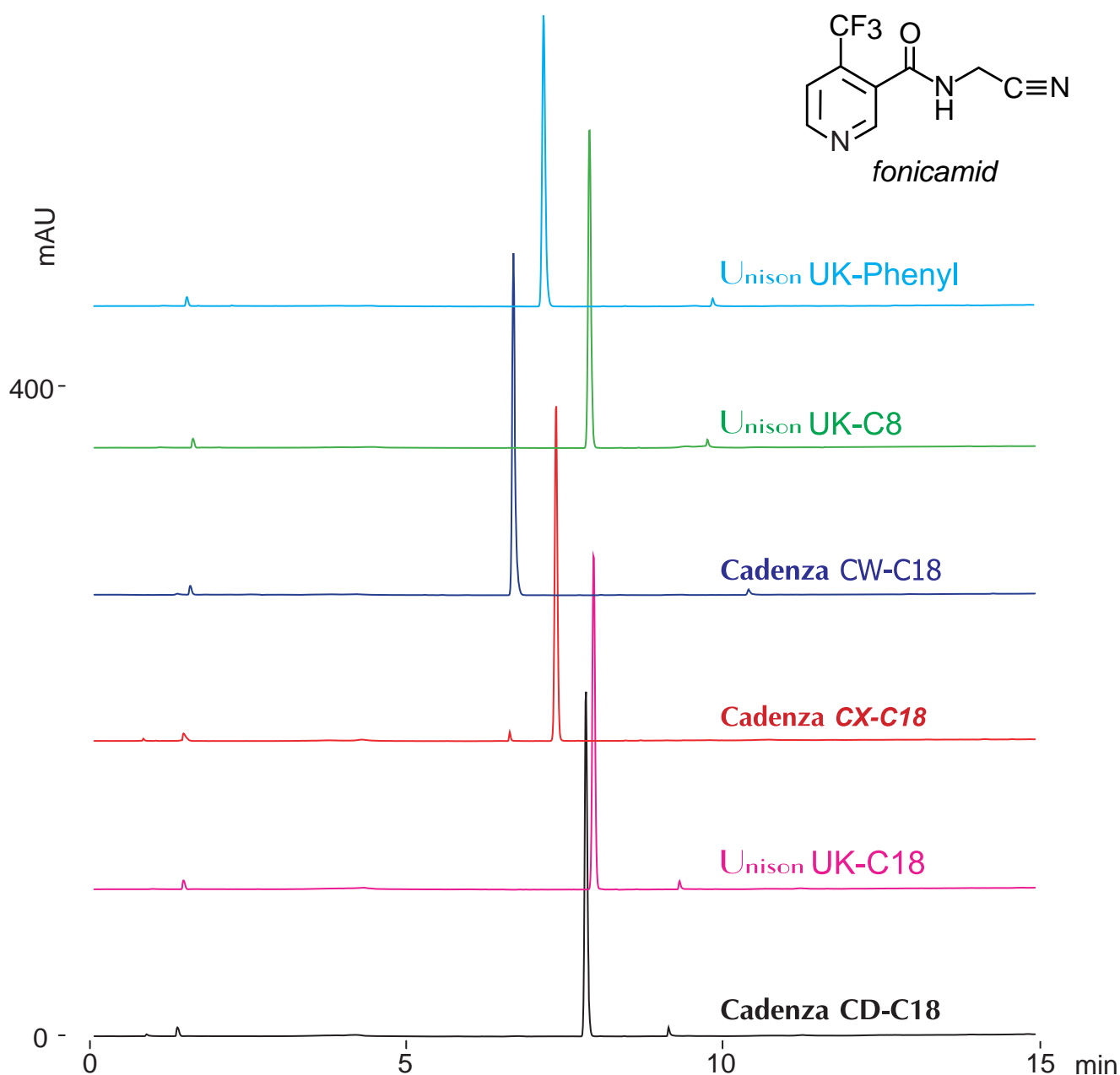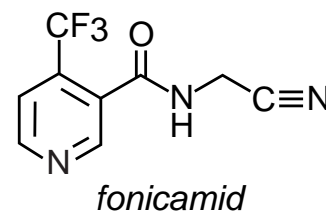

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0%B (0-1min), 0-70%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.15mg/mL, water / formic acid = 100 / 1)