



Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Abietic acid (Diterpenoids)

アビエチン酸 (ジテルペノイド)

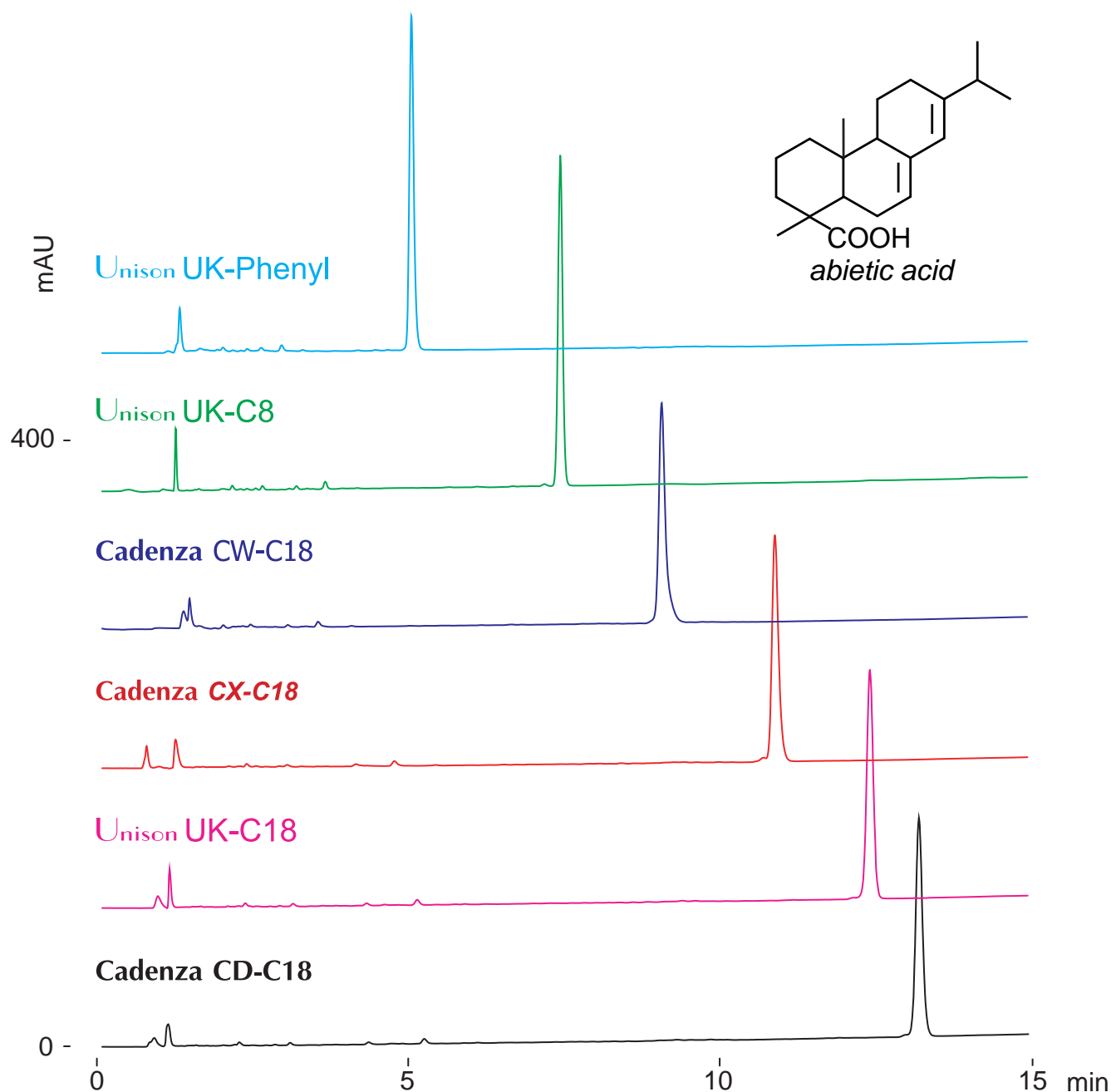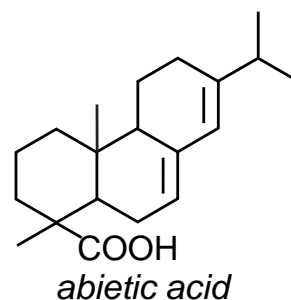

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

10-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 250nm

10µL (0.1mg/mL, water / acetonitrile / formic acid = 25 / 75 / 1)