
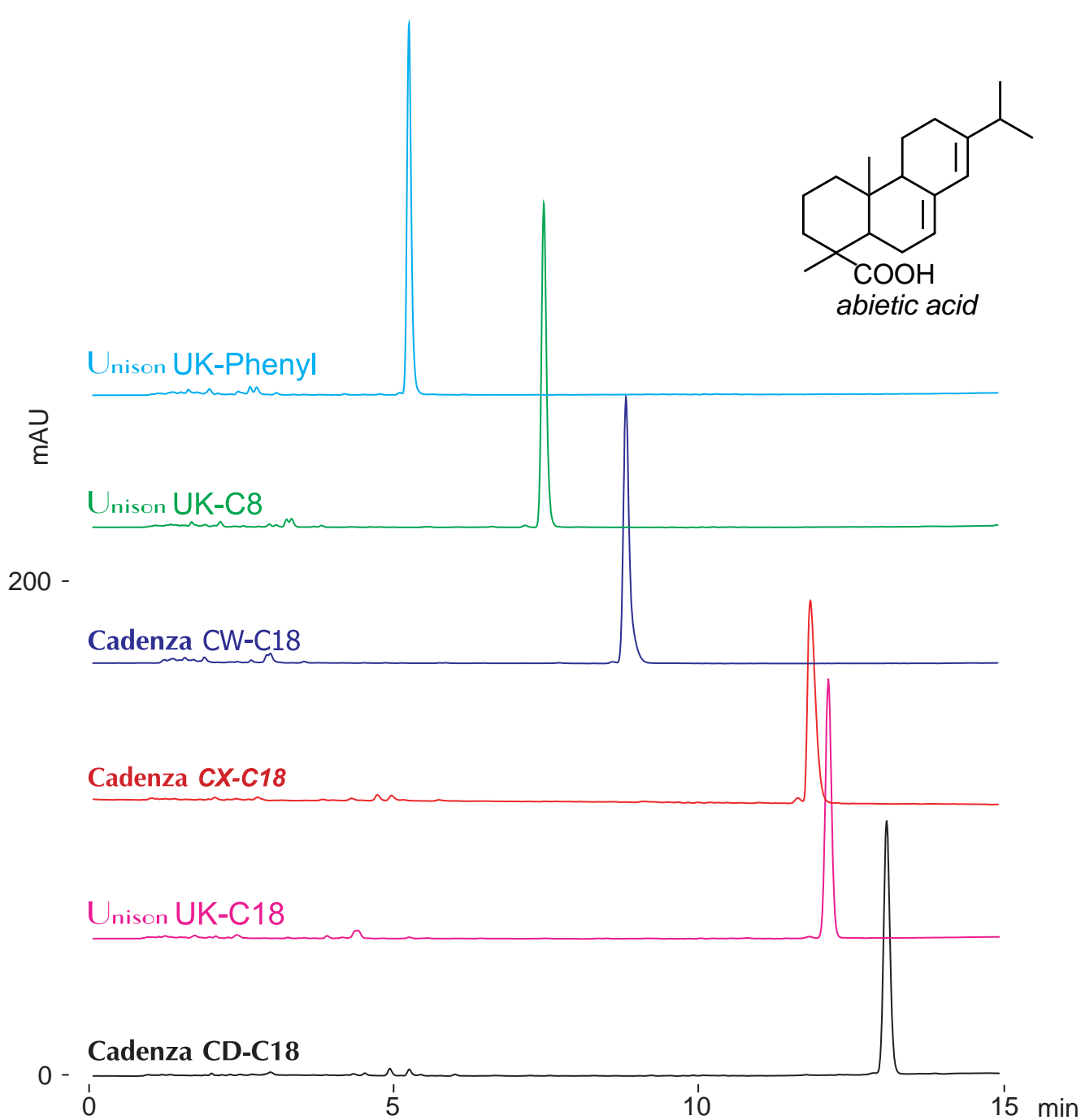


Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Abietic acid (Diterpenoids)
 アビエチン酸 (ジテルペノイド)

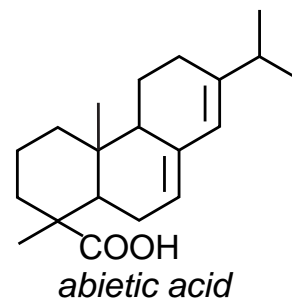

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 50 / 50

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 250nm

10µL (0.1mg/mL, 100mM ammonium formate / acetonitrile = 30 / 70)