


Unison UK-Phenyl Cadenza CL-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

2, 3, 4 - tert - Butylphenol (Butylphenol isomers)

2, 3, 4-tert-ブチルフェノール (ブチルフェノール構造異性体)

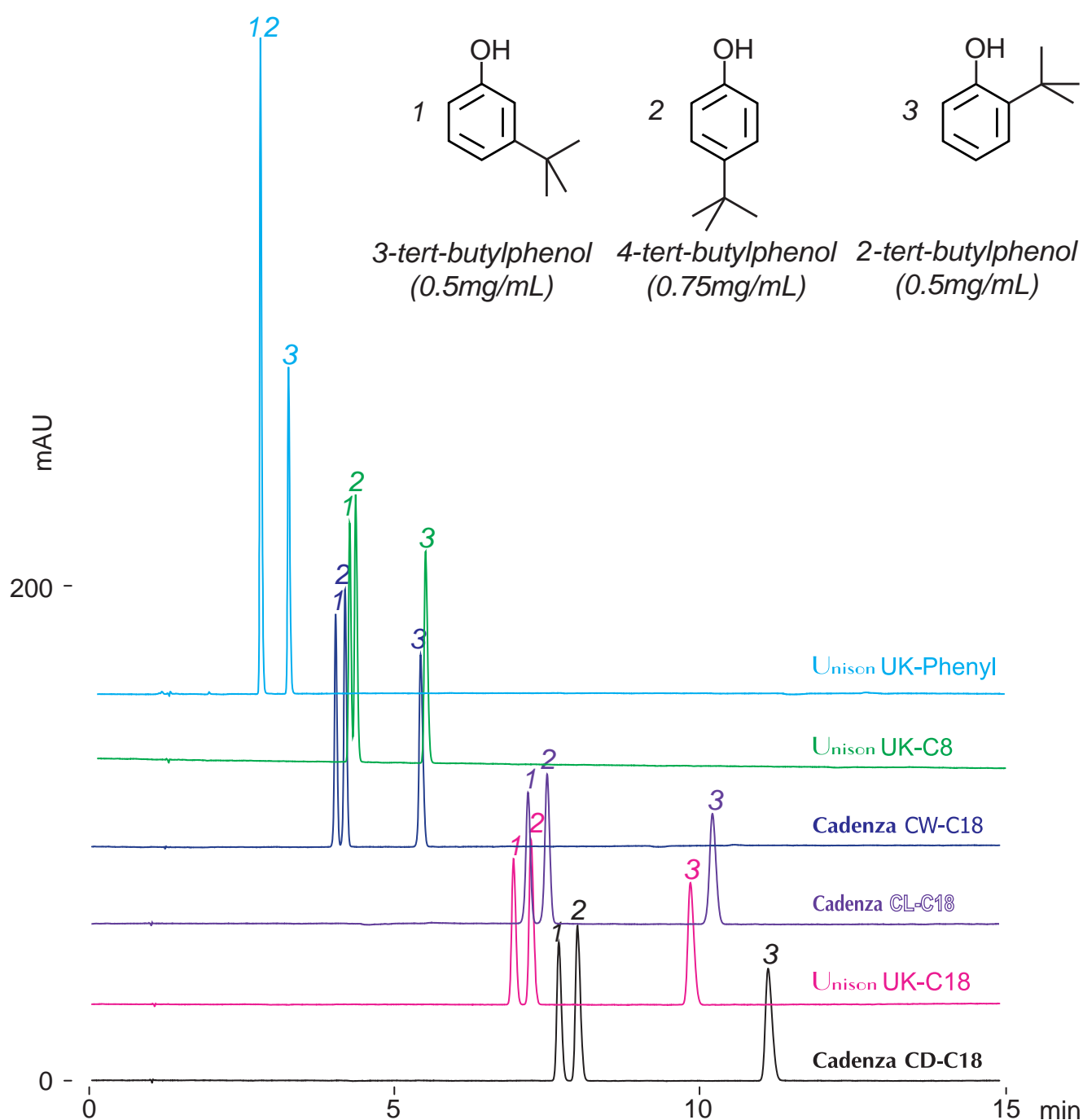

250 x 4.6 mm

A: 100mM ammonium formate / methanol = 90 / 10

B: 100mM ammonium formate / methanol = 10 / 90

57-62%B (0-40min)

1.0mL/min(20MPa), 50°C, 260nm

10µL (0.5-0.75mg/mL, 100mM ammonium formate / methanol = 70 / 30)