


Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Triclocarban, triclosan (Chlorophenyl derivatives)

トリクロカルバン, トリクロサン (クロロフェニル誘導体)

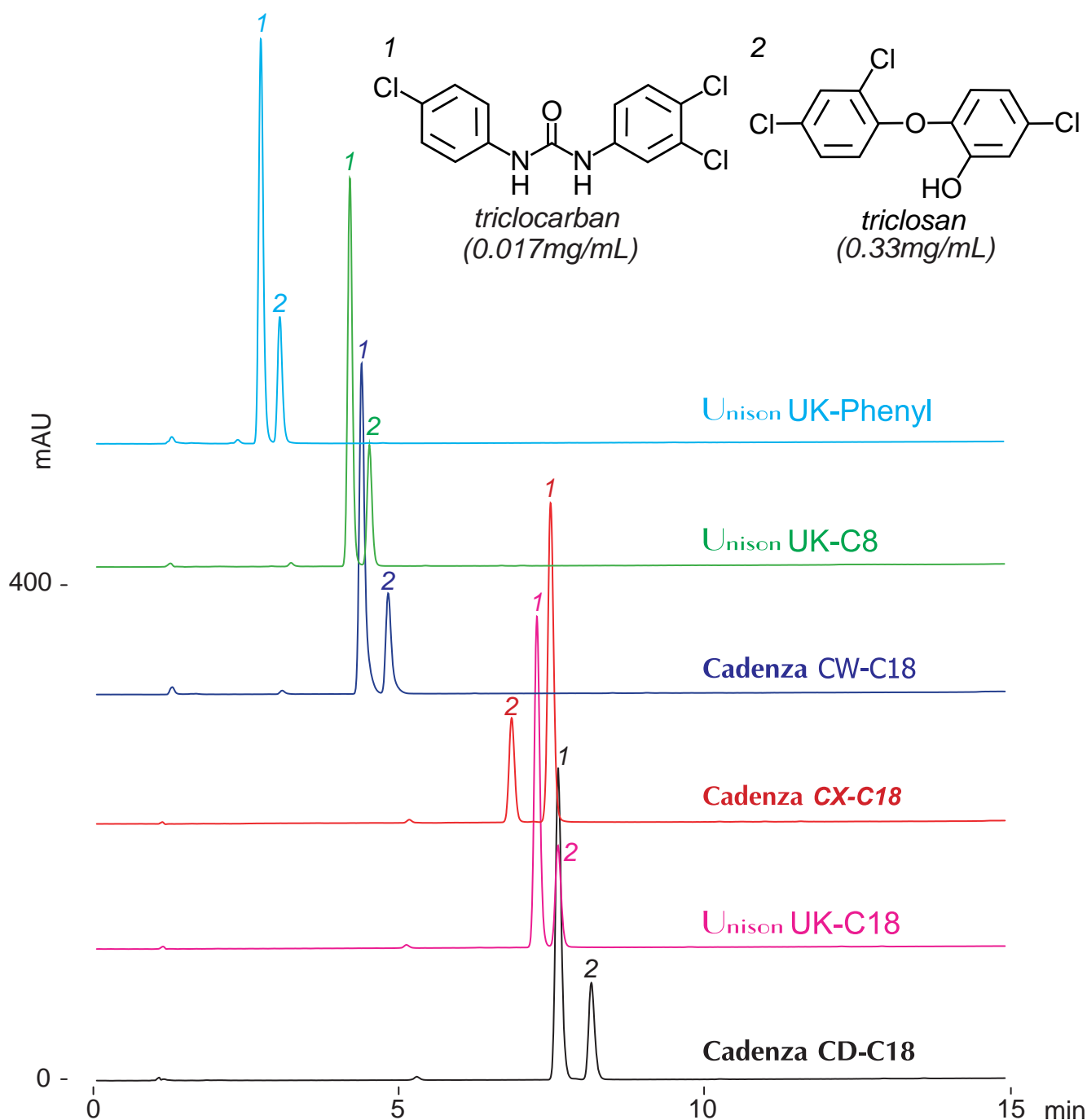

100 x 4.6 mm

A: water / methanol / formic acid = 50 / 50 / 0.1

B: water / methanol / formic acid = 10 / 90 / 0.1

50-100%B (0-15min)

1.0mL/min (15MPa), 37°C, 260nm

10µL (0.017-0.33mg/mL, water / methanol / formic acid = 10 / 90 / 1)