


Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Triclocarban, triclosan (Chlorophenyl derivatives)

トリクロカルバン, トリクロサン (クロロフェニル誘導体)

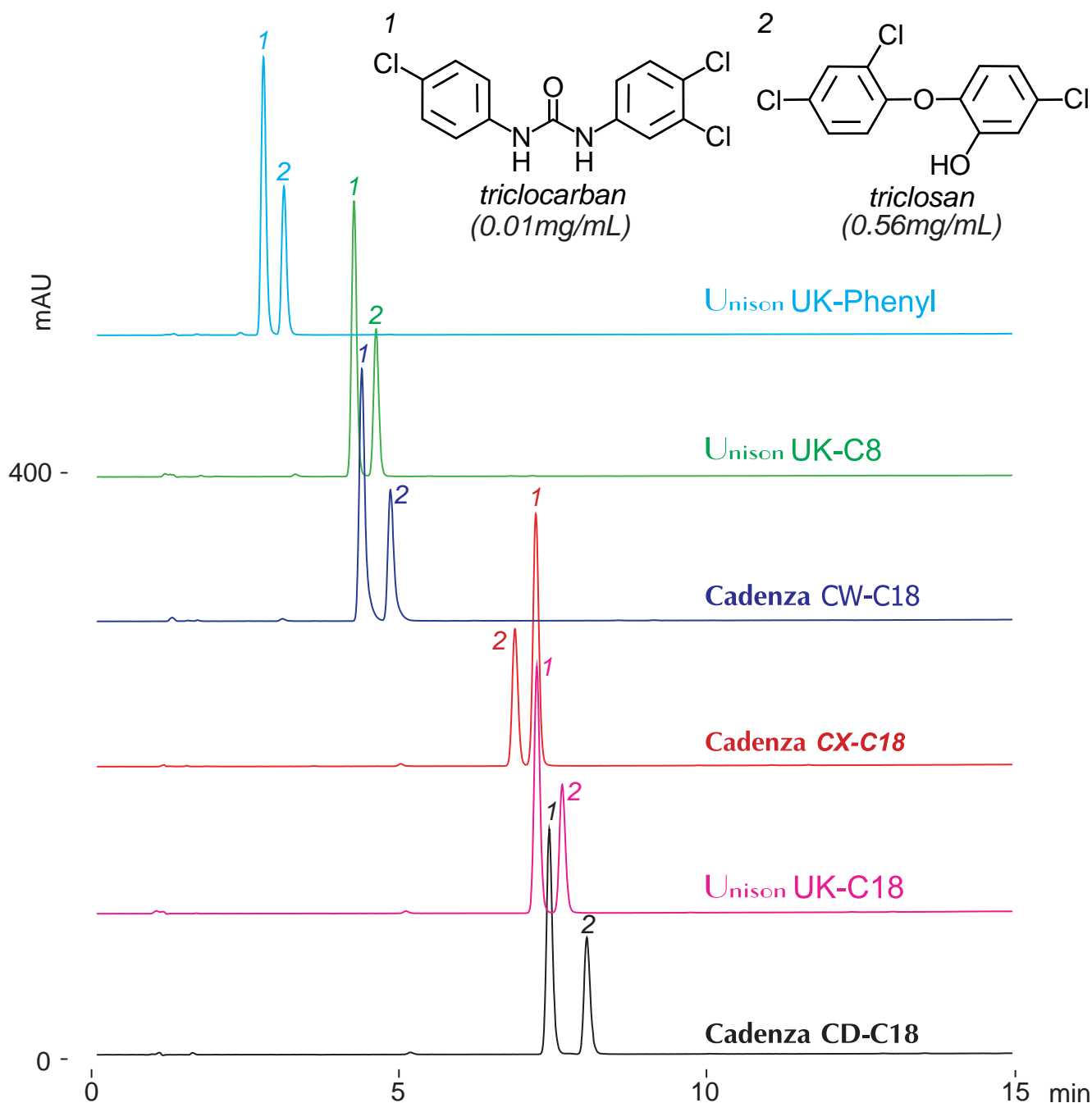

100 x 4.6 mm

A: 100mM ammonium formate / methanol = 50 / 50

B: 100mM ammonium formate / methanol = 10 / 90

50-100%B (0-15min)

1.0mL/min (15MPa), 37°C, 260nm

10µL (0.01-0.56mg/mL, 100mM ammonium formate / methanol = 10 / 90)