


Intrada Organic Acid

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: Acetic acid (Short-fatty acids)

酢酸 (短鎖脂肪酸類) の LC-MS 分析

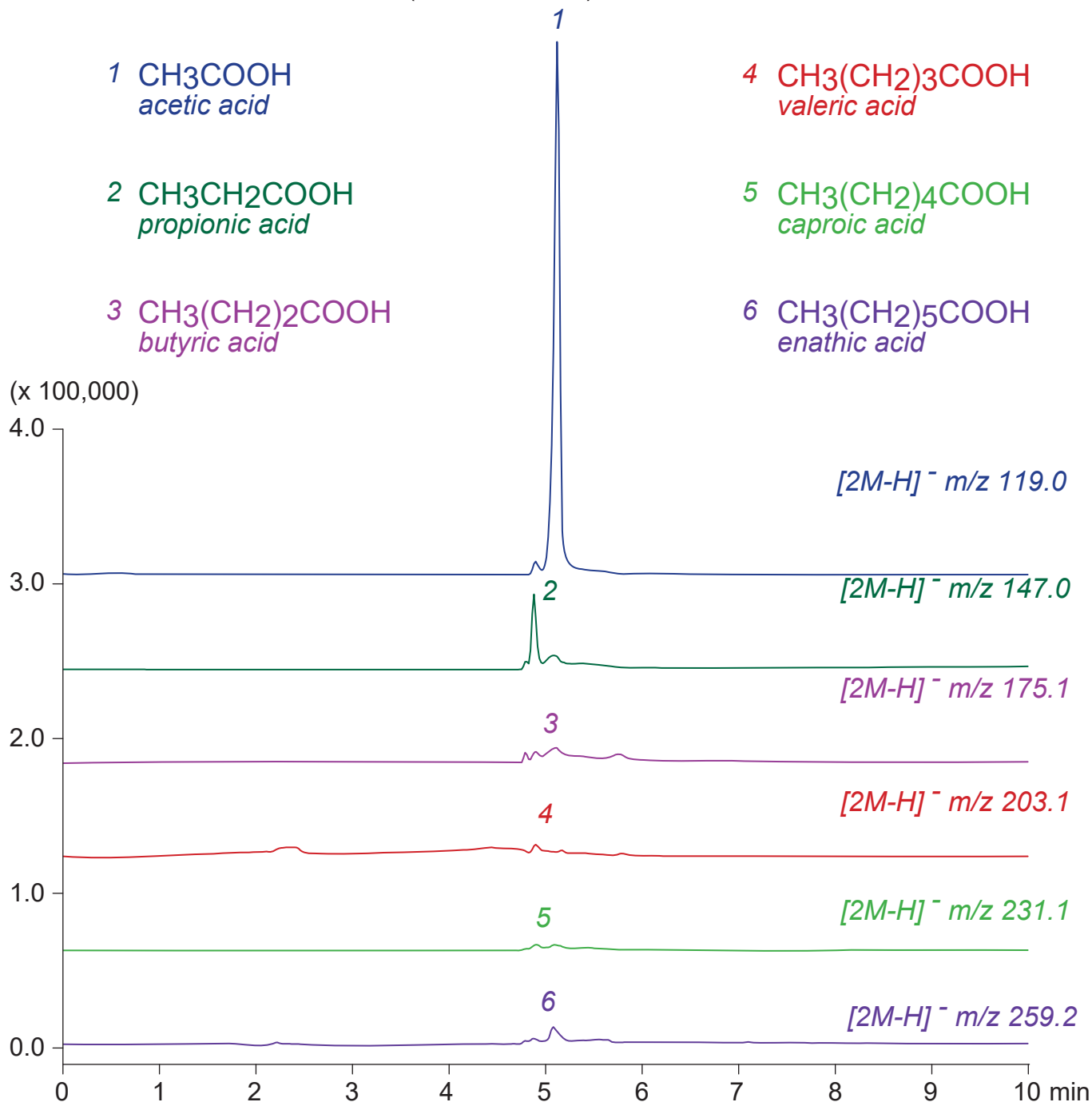

Intrada Organic Acid, 150 x 2 mm

A: acetonitrile

B: 10mM ammonium bicarbonate

0-80%B (0-7min), 80%B (7-10min)

0.2mL/min (5MPa), 37°C, Single Quad. MS (ESI)

1μL (0.4mg/mL, **ACN / 100mM HClO₄ = 50 / 50**)