


Scherzo SS-C18

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Short-fatty acids

短鎖脂肪酸類の LC-MS 分析

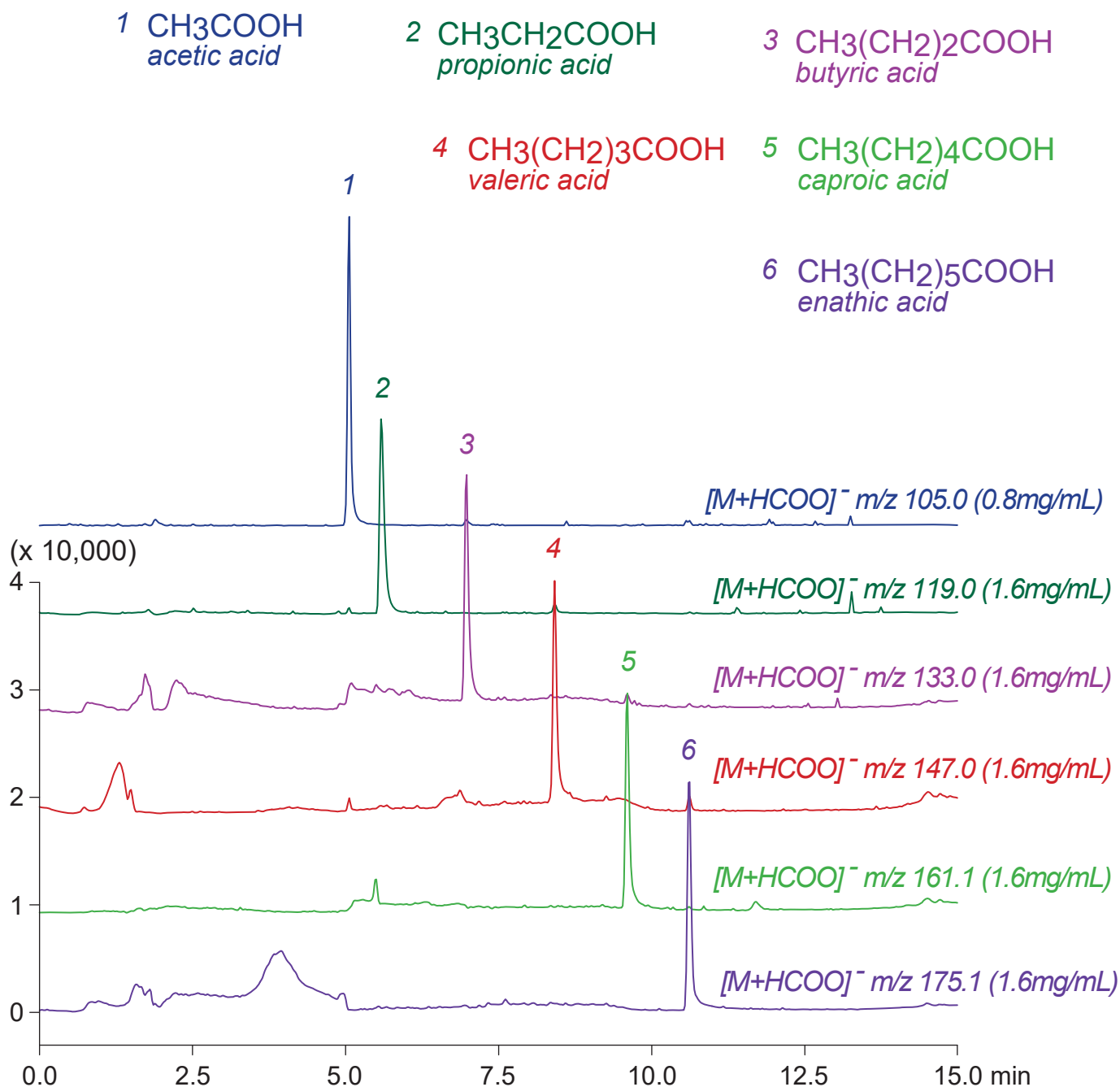

Scherzo SS-C18, 150 x 3 mm

A: 5mM ammonium formate

B: acetonitrile / water / formic acid =90 / 10 / 0.5

0%B (0-1min), 0-100%B (1-12min), 100%B (12-15min)

0.4mL/min (10MPa), 37°C, Single Quad. MS (ESI)

10µL (0.8-1.6mg/mL, 100mM ammonium formate)