


Unison UK-C18
Cadenza CD-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

Red yeast rice extract (Lovastatin, monacolin K)
紅麹抽出物 (ロバスタチン, モナコリンK)

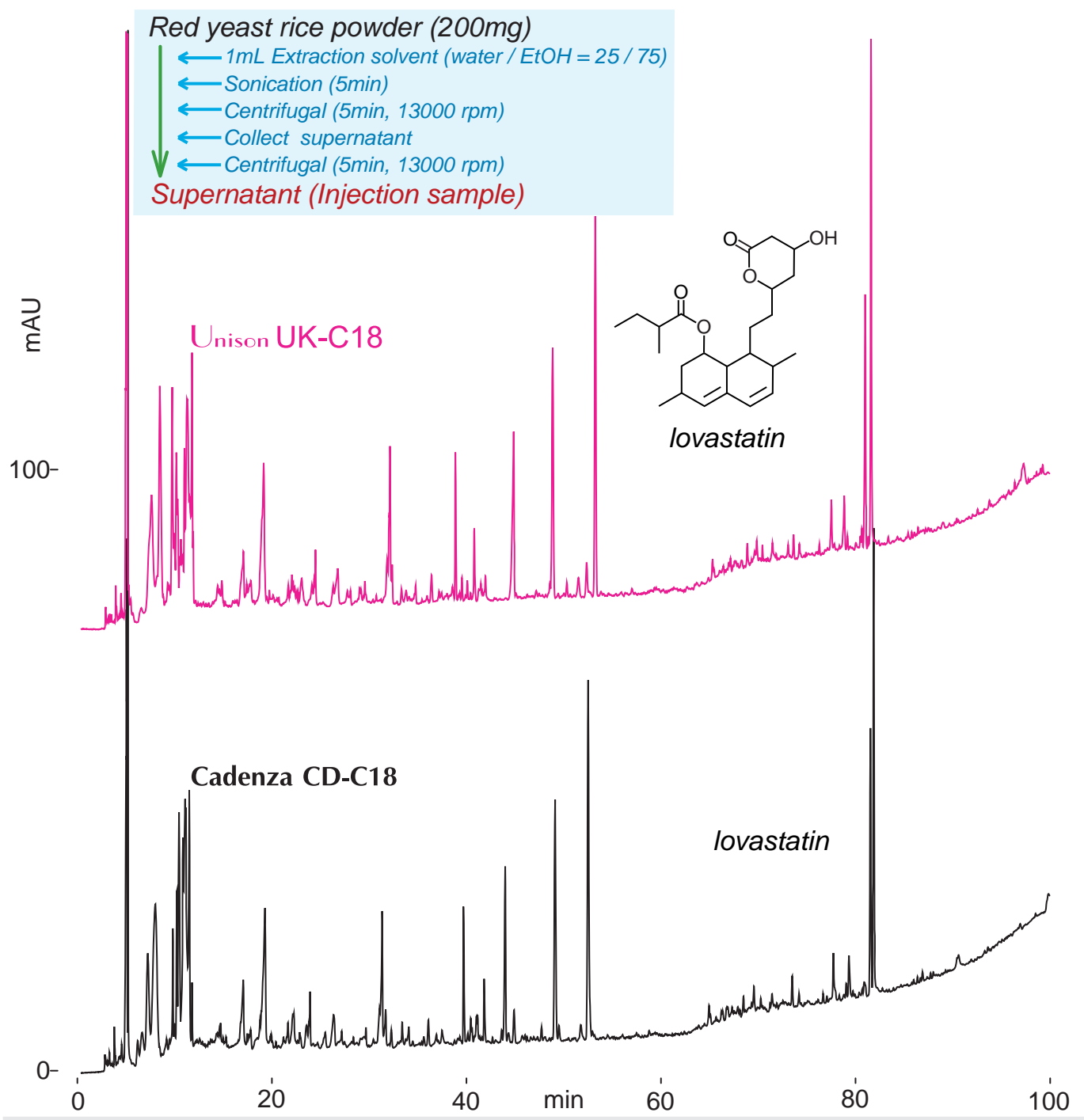

250 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0%B (0-1min), 0-40%B (1-60min), 40-100%B(60-95min), 100%B (95-100min)

1.0mL/min (16MPa), 250nm, 37°C

10µL (Red yeast rice powder 200mg/mL, water / EtOH = 25 / 75)