


Unison UK-Phenyl

100 x 4.6 mm


HPLC COLUMN APPLICATION

Organic Redox-Flow Battery positive electrolyte (TEMPO Free Radical)


有機レドックスフロー電池正極電解質 (TEMPOフリーラジカル)

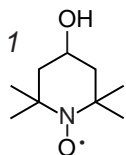


1
4-hydroxy-2,2,6,6-tetramethyl
piperidine 1-oxyl free radical

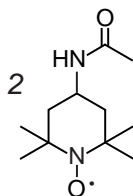


2
4-acetamido-2,2,6,6-tetramethyl
piperidine 1-oxyl free radical

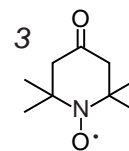


3
4-oxo-2,2,6,6-tetramethyl
piperidine 1-oxyl free radical

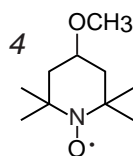

4
4-methoxy-2,2,6,6-tetramethyl
piperidine 1-oxyl free radical

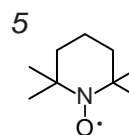

5
2,2,6,6-tetramethyl
piperidine 1-oxyl free radical

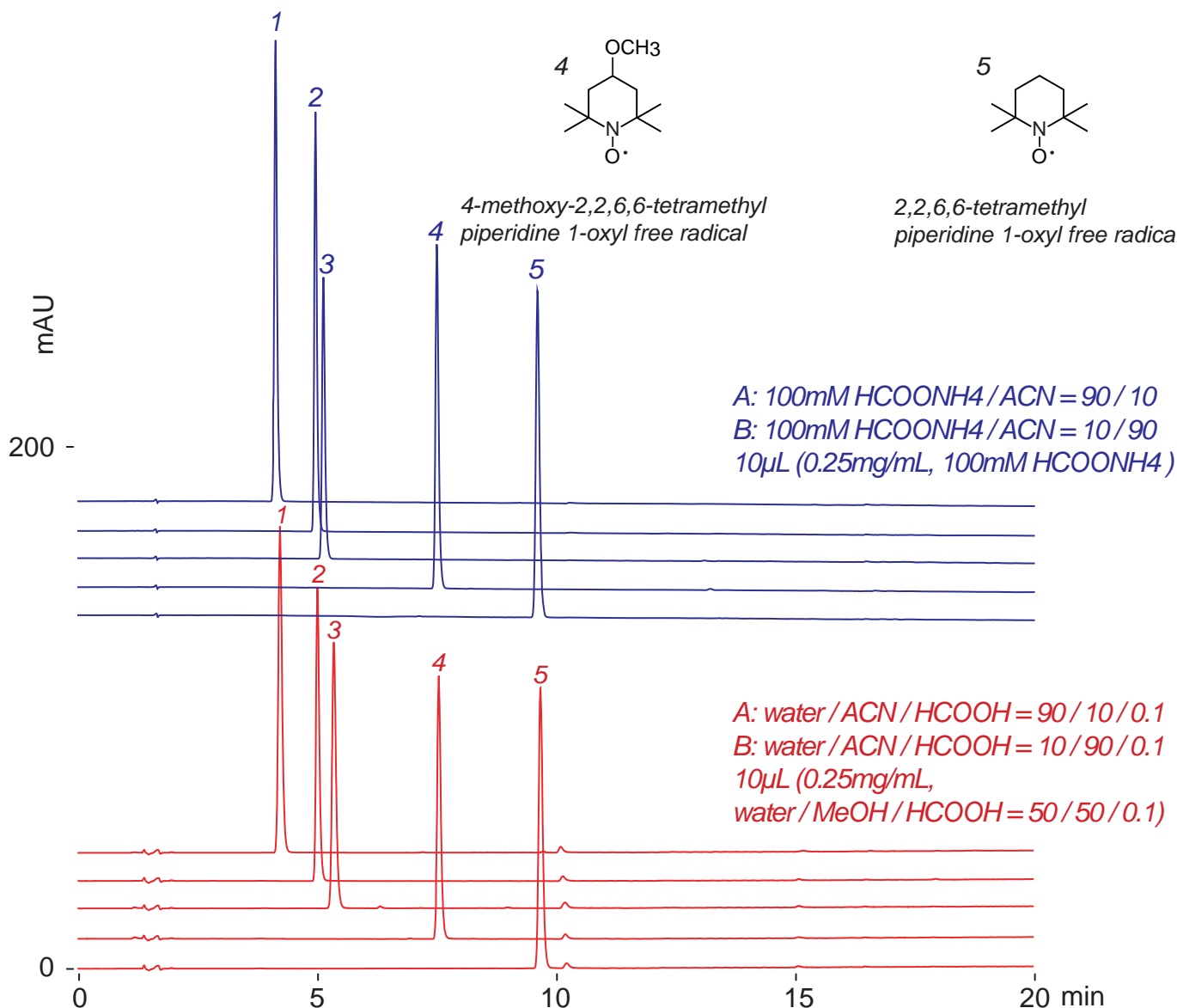

Unison UK-Phenyl, 100 x 4.6 mm
 0%B (0-1min), 0-80%B (1-20min)
 1.0mL/min (9MPa), 37°C, 260nm