


Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Organic Redox-Flow Battery positive electrolyte (10-Acetylphenothiazine)

有機レドックスフロー電池正極電解質 (10-アセチルフェノチアジン)

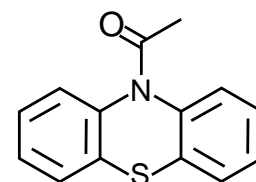

10-acetylphenothiazine

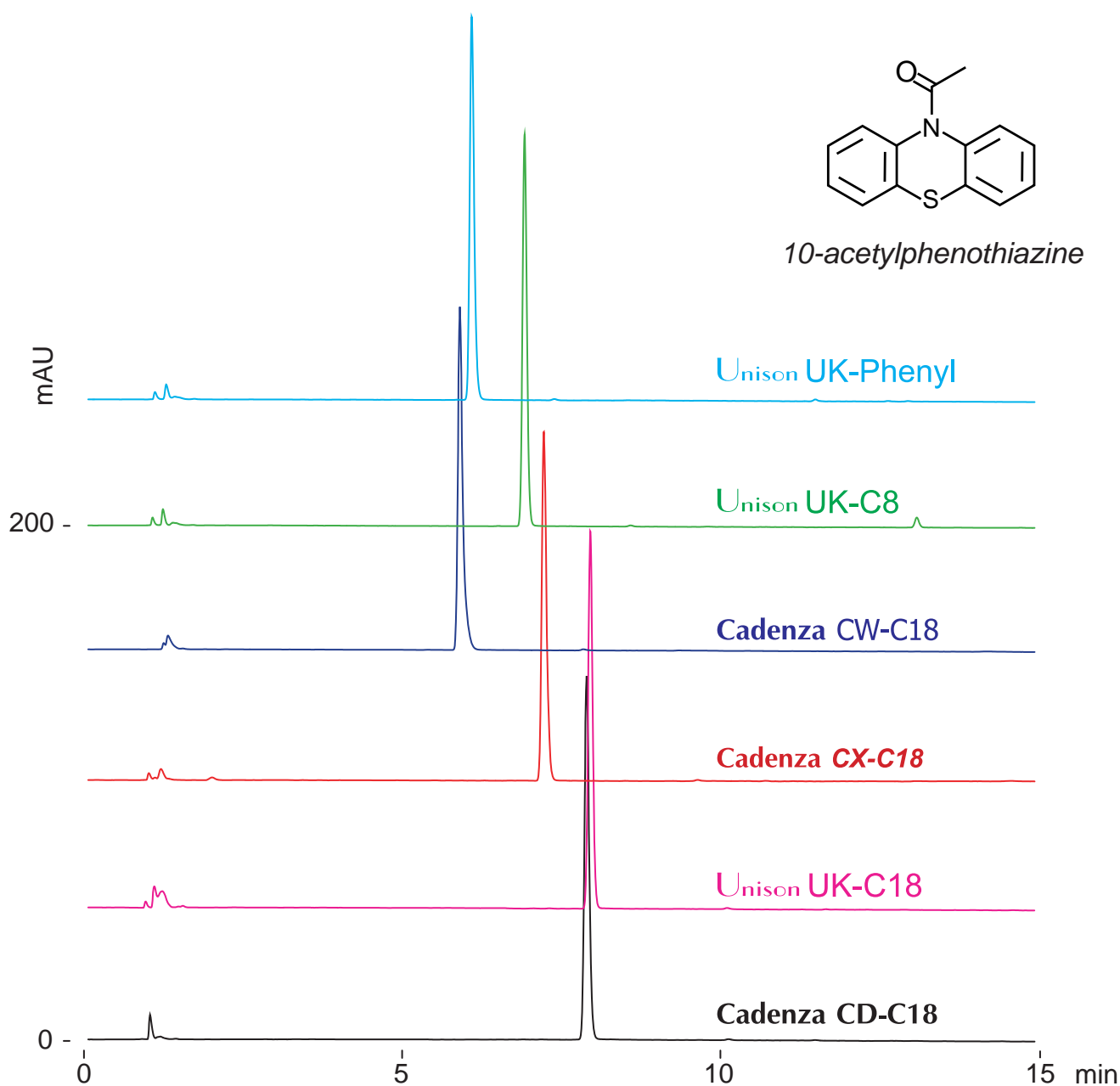

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

25-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.05mg/mL, 100mM ammonium formate / acetonitrile = 30 / 70)