



**Scherzo SW-C18MF** 150 x 2 mm **HPLC COLUMN APPLICATION**

LC-MS: Organic Redox-Flow Battery negative electrolyte (anthraquinonesulfonic acids)


有機レドックスフロー電池負極電解質 (アントラキノンスルホン酸類) の LC-MS 分析

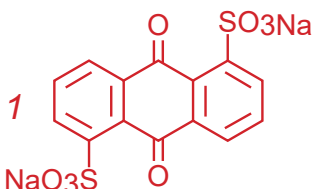


anthraquinone-1,5-disulfonic acid sodium salt (25µg/mL)

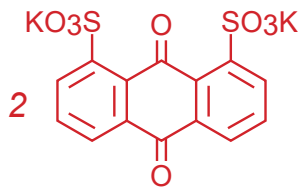

anthraquinone-1,8-disulfonic acid sodium salt (20µg/mL)

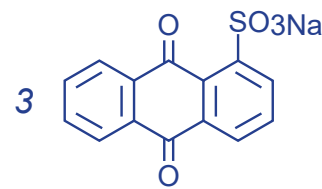

anthraquinone-1-sulfonic acid sodium salt (10µg/mL)

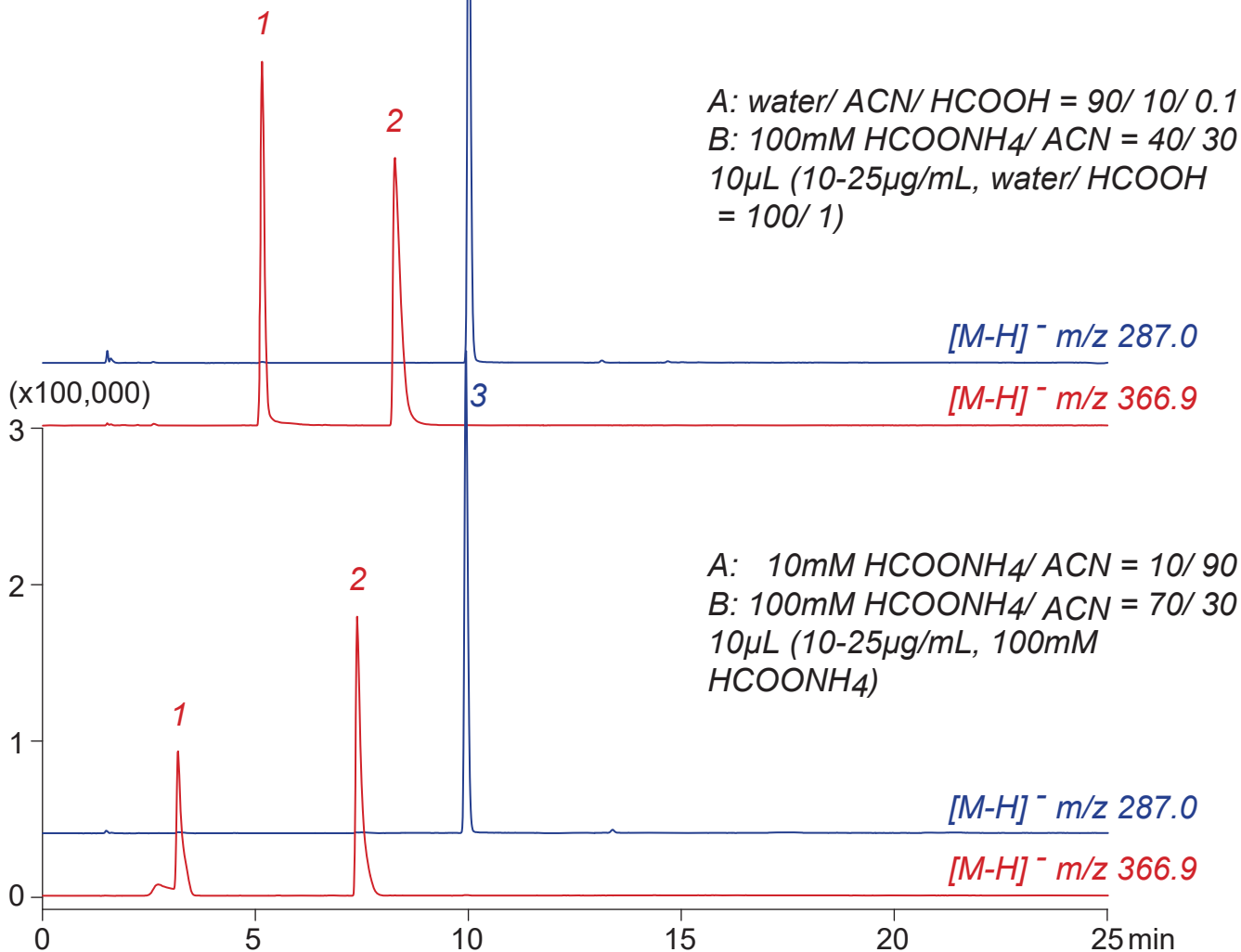

Scherzo SW-C18MF, 150 x 2 mm  
 0-100%B (0-10min), 100%B (10-25min)  
 0.2mL/min (11MPa), 37°C, Single Quad. MS (ESI)