

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18

Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Donor building block for organic semiconductor

(Dibromothiophene isomers)

有機半導体ドナービルディングブロック

(ジブロモチオフェン異性体)

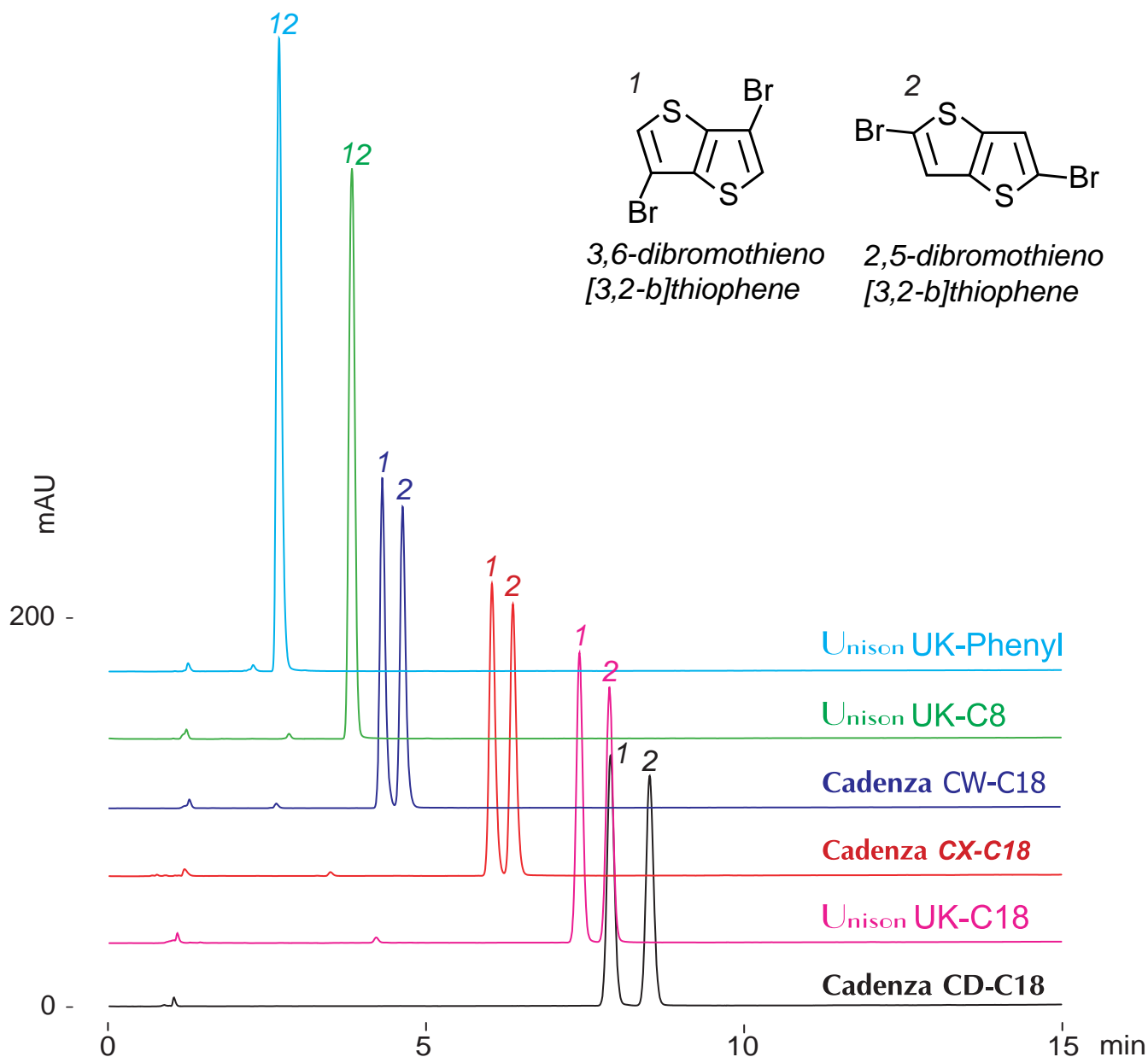

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

40-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (0.1mg/mL, water / acetonitrile / formic acid = 20 / 80 / 1)